




Pro for Pros

www.AEOLUSpet.com
www.KANI-STAR.com
www.TOEXpet.com

KANISTAR GLOBAL CATALOG

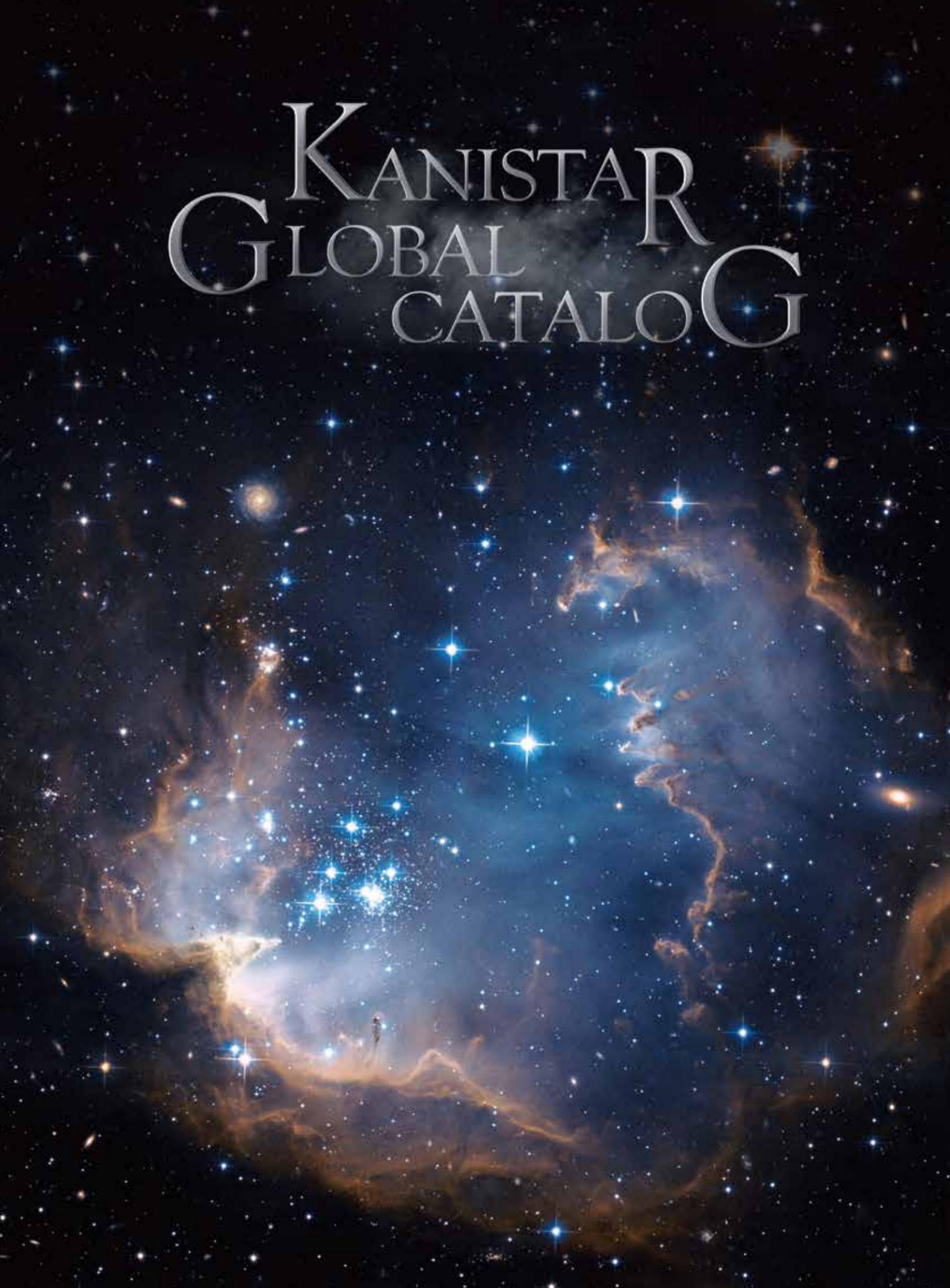



Professional Walk-in Kennel System



Professional Walk-in Kennel System

USKN Series is our Newest Walk-in Kennel System. All the components of this series can be customized to fit your needs and space. Various combination can be made of front door, side plates and accessories. We can make these of various materials. The kennels in your dreams are here. Various combinations of stainless steel, high density HDPE, tempered glass and/or powder coated steel leaves you with many different material choices. Swivel feeding door, one-click magnetic latch, easy-to-lift guillotine door, etc.; many user-friendly designs can be incorporated into this system. We welcome all inquiries, requirements and customized orders. Our engineers and designers will be happy to help lay out 2D and 3D renderings for your consideration.



Stainless Steel Feeding Bowl and Mount USKN-FDB-150 & USKN-FDSB-002



Stainless Steel Name Card Holder 5"× 3" or 8"× 5" USKN-NMHD-3/5 & USKN-NMHD-5/8



Guillotine Door USKN-GD71-1830

Custom Kennel Displays




Compared to normal cages, the kennel could bring larger space to bigger animals, reduce the stress on medium or small sized animals, provide opportunity to house family groups together, or satisfy the clients who are looking forward to giving luxury space to their pets.

Many details are taken into account during the boarding/breeding kennel facility planning, a good architect could help you a lot and reduce possibilities of many potential problems in the future. If you do not have such resource in hand or you are not fully ready for this yet, KaniStar may help you with some basic rendering layouts.

Our designing team has great flexibility and coordinating spirits about turning clients' dreams into reality.

Our engineering team is very experienced, not merely capable in making the kennel sturdy and long lasting, but also able to consider the actual conditions and get all the products adapted to their conditions.

Whether you're looking for a new kennel installation, a partial or complete renovation, special dog kennel runs, a new dog shelter set-up, a luxury pet hotel establishment, or a solution to challenging indoor/outdoor space, we can help!

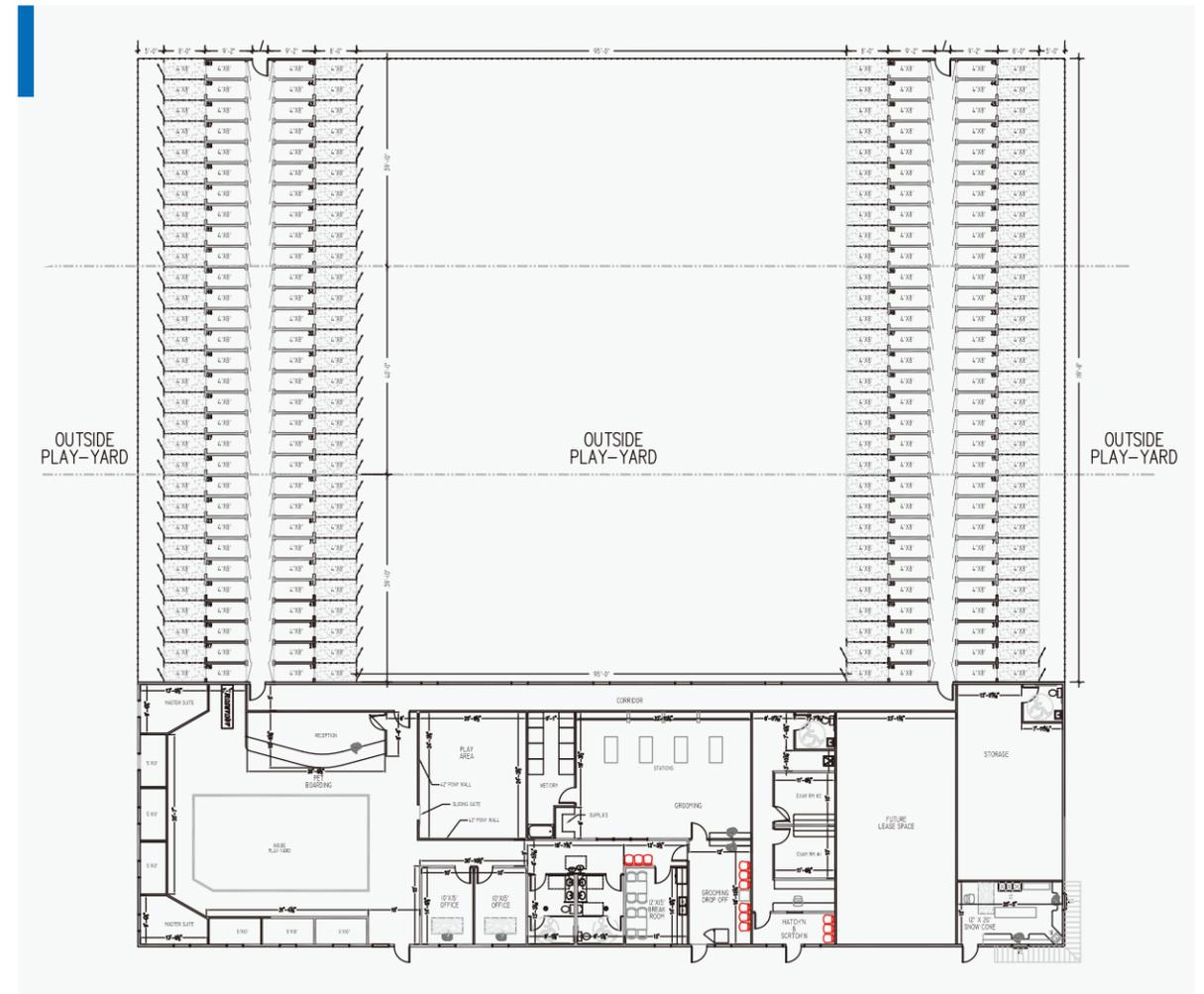

Custom Display





Professional Walk-in Kennel (USKN series)

Pet friendly: **Big place for pets to live and reduce the stress.**

Steel floor with Rubberizing coating: Environmental friendliness(no VOC), high tensile strength, elongation, flexibility, wear resistance, aging resistance and anticorrosion.

Drainage system: easy to clean the kennel and clear wastes.

Side wall: stainless steel frame and HDPE panel. Solid structure and easy maintenance.

Easy installation: Module assembly. Easy to assemble and connect cubicles.

Middle dividing guillotine gate can be controlled by a double-control pulley system to connect two rooms or separate the rooms. Easy to adjust room allocation according to pet population. Also make it an easy job to clean the kennel by easily moving the pets in one room to the other room.

Practical and user-friendly accessories: Swivel feeding door, one-click magnetic latch, etc.

Customizable for you to choose different material and displays that are most suitable to your place.



Comfortable space for animals.



Front view



Side view



Top grid



Double decker kennel



Drainage system available



Semi-finished kennel



Side wall: stainless steel frame and HDPE panel.



Steel floor with rubberizing coating: Easy to clean and maintain



Filter cover for protecting the drainage trough.



Middle guillotine door



Control knob of the double control pulley system with weight balancer.










Heavy duty hinge











One-click magnetic latch




Fully welded bars

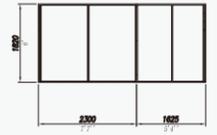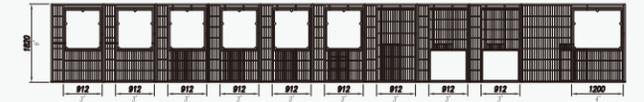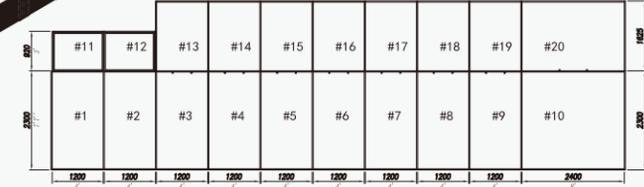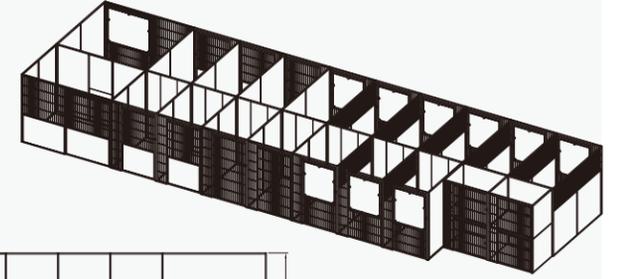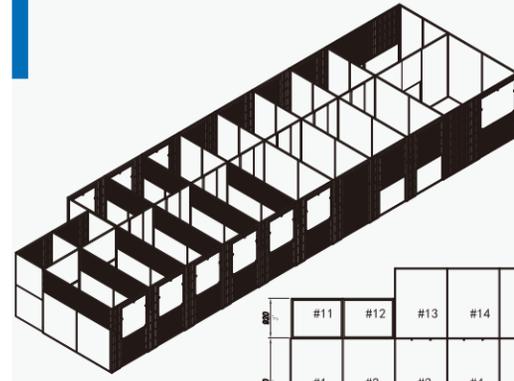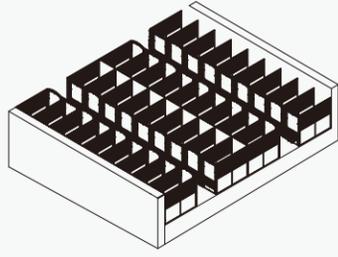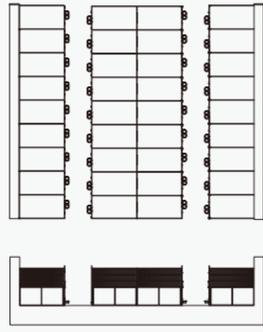



unit: mm/ft

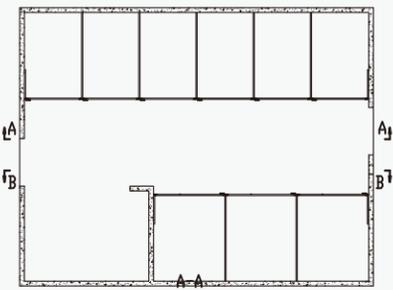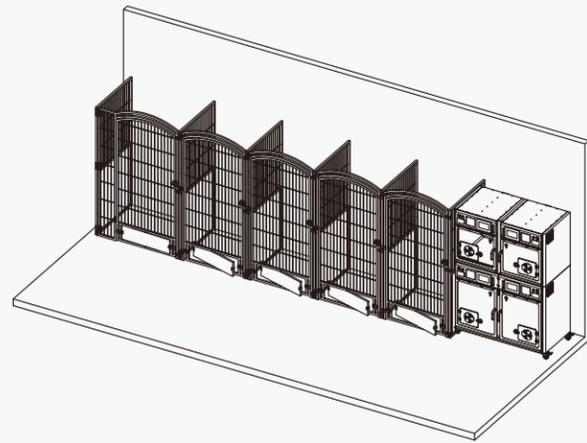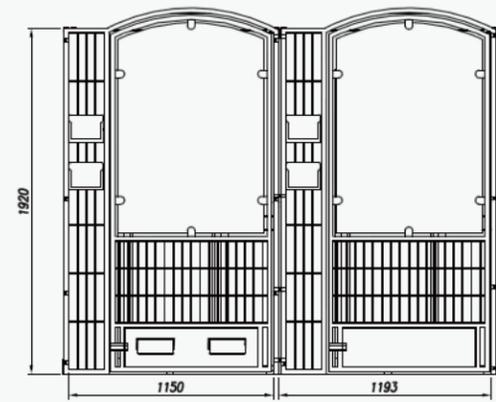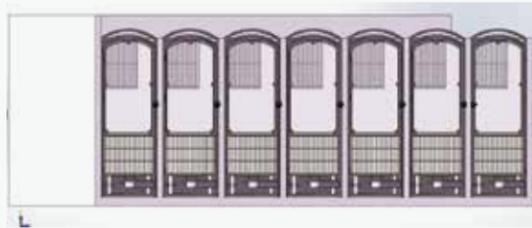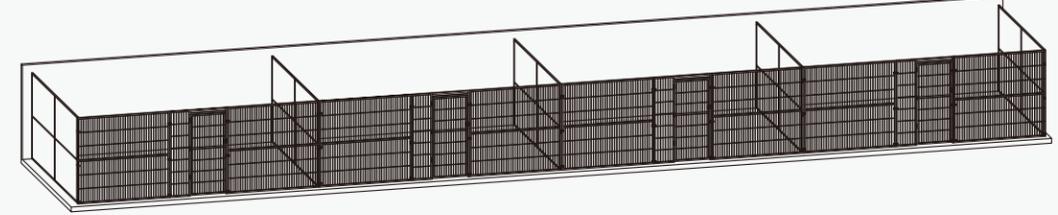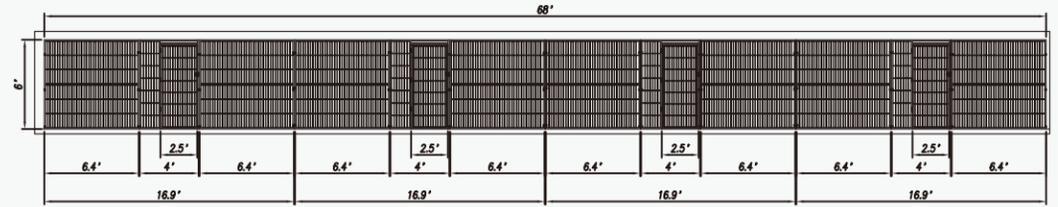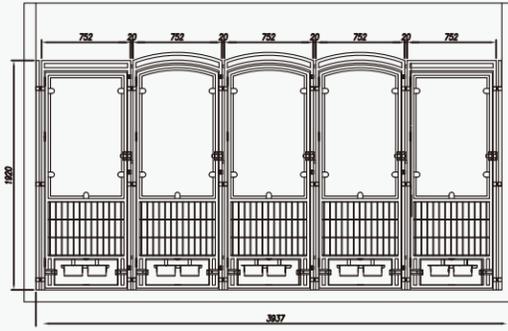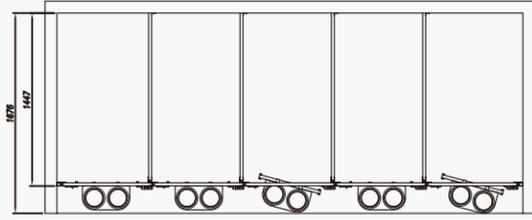


B-B

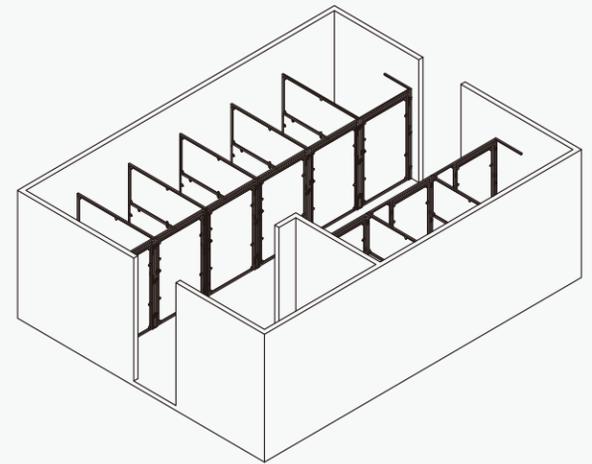

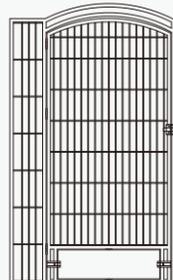Category	Item Number	Descriptions	Remarks	Dimension	Accessories	
Stainless Steel Grid Front Door	USKN-FD11-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Full Stainless Steel Grid, Frame made from 1"×1"(25mm×25mm) tube, 1/16" (1.5mm) thick.	Flat Top Frame	47.2"W×75.6"H (120cmW×192cmH)	All the front gate includes a swivel door underneath, exclusive of the feeding bowls and bowl mount, including two door-to-door connector clamps.	
	USKN-FD11-4872T	Front Gate, 47.2"W×75.6"H(120cm×192cm) Full Stainless Steel Grid, Frame made from 1"×1"(25mm×25mm) tube, 1/10"(2.5mm) thick.		47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD12-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Stainless Steel Grid + Arch Top, Frame made from 1"×1"(25mm×25mm) tube, 1/16"(1.5mm) thick	Arch Top Frame	47.2"W×79.8"H (120cmW×202cmH)		
	USKN-FD12-4878T	Front Gate, 47.2"W×79.8"H(120cm×202cm), Stainless Steel Grid + Arch Top, Frame made from 1"×1"(25mm×25mm) tube, 1/10"(2.5mm) thick		47.2"W×79.8"H (120cmW×202cmH)		
Stainless Steel Grid + Tempered Glass Front Door	USKN-FD31-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Stainless Steel Grid + 1/3"(8mm) Thick Tempered Glass + 1/10" (2.5mm) thick S/S Tube + Flat Top	Flat Top Frame	47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD51-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Stainless Steel Grid + 1/3"(8mm) Thick Tempered Glass + 1/16"(1.5mm) thick S/S Tube + Flat Top		47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD32-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Stainless Steel Grid + 1/3"(8mm) Tempered Glass + 1/10"(2.5mm) Thick Tube + Arch Top	Arch Top Frame	47.2"W×79.8"H (120cmW×202cmH)		
	USKN-FD52-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Stainless Steel Grid + 1/3"(8mm) Tempered Glass + 1/16"(1.5mm) Thick Tube + Arch Top		47.2"W×79.8"H (120cmW×202cmH)		
Powder Coated Grid Front Door	USKN-FD21-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Powder Coated Grid, Frame made from 1"×1" (25mm×25mm) tube, 1/16"(1.5mm) thick.	Flat Top Frame	47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD21-4872T	Front Gate, 47.2"W×75.6"H(120cm×192cm), Powder Coated Grid, Frame made from 1"×1" (25mm×25mm) tube, 1/10" (2.5mm)thick.		47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD22-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Powder Coated Grid + Arch Top, Frame made from 1"×1"(25mm×25mm) tube, 1/16"(1.5mm) thick	Arch Top Frame	47.2"W×79.8"H (120cmW×202cmH)		
	USKN-FD22-4878T	Front Gate, 47.2"W×79.8"H(120cm×202cm), Powder Coated Grid+ Arch Top, Frame made from 1"×1"(25mm×25mm) tube, 1/10"(2.5mm) thick		47.2"W×79.8"H (120cmW×202cmH)		
Powder Coated Grid + Tempered Glass Front Door	USKN-FD41-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Powder Coated Grid + 1/3"(8mm) Thick Tempered Glass + 1/10" (2.5mm)thick S/S Tube + Flat Top	Flat Top Frame	47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD61-4872	Front Gate, 47.2"W×75.6"H(120cm×192cm), Powder Coated Grid + 1/3"(8mm) Thick Tempered Glass + 1/16" (1.5mm) thick S/S Tube + Flat Top		47.2"W×75.6"H (120cmW×192cmH)		
	USKN-FD42-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Powder Coated Grid + 1/3"(8mm) Tempered Glass + 1/10"(2.5mm) Thick Tube + Arch Top	Arch Top Frame	47.2"W×79.8"H (120cmW×202cmH)		
	USKN-FD62-4878	Front Gate, 47.2"W×79.8"H(120cm×202cm), Powder Coated Grid + 1/3"(8mm) Tempered Glass + 1/10"(2.5mm) Thick Tube + Arch Top		47.2"W×79.8"H (120cmW×202cmH)		
Stainless Steel Side Panel	USKN-SW10-9672	Side Panel, 90.5"W×75.6"H(230cm×192cm), Stainless Steel Grid + 0.4" (10mm) thick HDPE Panel, Frame made from 1/16"(1.5mm) thick tubes	Stainless Steel Grid on top + Solid HDPE Panel at bottom	90.5"W×75.6"H (230cmW×192cmH)		
	USKN-SW10-9672T	Side Panel, 90.5"W×75.6"H(230cm×192cm), Stainless Steel Grid + 0.4" (10mm) thick HDPE Panel, Frame made from 1/10" (2.5mm)thick tubes		90.5"W×75.6"H (230cmW×192cmH)		
	USKN-SW11-9672	Side Panel, 90.5"W×75.6"H(230cm×192cm), Stainless Steel Top Grid + Solid 1/16"(1.5mm) Thick Stainless Steel Panel, Frame made from 1/16"(1.5mm) thick tubes		90.5"W×75.6"H (230cmW×192cmH)		

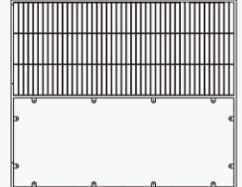Category	Item Number	Descriptions	Remarks	Dimension	Accessories	
Stainless Steel Side Panel	USKN-SW11-9672T	Side Panel, 90.5"W×75.6"H(230cm×192cm), Stainless Steel Top Grid + Solid 1/16"(1.5mm) Thick Stainless Steel Panel, Frame made from 1/10"(2.5mm) thick tubes	Stainless Steel Grid on top + Solid Stainless Steel Panel at bottom	90.5"W×75.6"H (230cmW×192cmH)	ground-panel connection and the anchor bolts for installation	
Powder Coated Side Panel	USKN-SW20-9672	Side Panel, 90.5"W×75.6"H(230cm×192cm), Powder Coated Grid + 0.4"(10mm) thick HDPE Panel, Frame made from 1/16"(1.5mm) thick tubes	Powder Coated Grid on top + Solid HDPE Panel at bottom	90.5"W×75.6"H (230cmW×192cmH)	All side panel includes 3 L shape connectors for side panel - wall connection, 4 side panel - front door connectors, 2 side panel - top grid connectors, one full length bottom beam for ground-panel connection and the anchor bolts for installation, also the screws for HDPE panel attached to the frame	
	USKN-SW20-9672T	Side Panel, 90.5"W×75.6"H(230cm×192cm), Powder Coated Grid + 0.4"(10mm) thick HDPE Panel, Frame made from 1/10"(2.5mm) thick tubes	90.5"W×75.6"H (230cmW×192cmH)			
	USKN-SW22-9672	Side Panel, 90.5"W×75.6"H(230cm×192cm), Powder Coated Grid + Solid 1/16"(1.5mm) Thick Powder Coated Panel, Frame made from 1/16"(1.5mm) thick tubes	Powder Coated Grid on top + Solid Powder Coated Panel at bottom	90.5"W×75.6"H (230cmW×192cmH)		
Guillotine Door	USKN-SW22-9672T	Side Panel, 90.5"W×75.6"H(230cm×192cm), Powder Coated Grid + Solid 1/16"(1.5mm) Thick Powder Coated Panel, Frame made from 1/10"(2.5mm) thick tubes	90.5"W×75.6"H (230cmW×192cmH)			
	USKN-GD71-1830	18"(45cm) wide, 30"(75cm) tall Guillotine door for walk-in kennels, with two durable stainless steel pulleys, counter weight, aluminum tracks and steel cable included.				
Accessories	USKN-DC-FLT-19/ USKN-DC-FLT-21	Straight Bottom Beam, for the connection between side panel and ground, will make the space completely isolated from neighboring kennels and prevent washing fluid leaking.	Suffix 19 means the tube is 19mm×19mm (3/4"×3/4"), which fits the 1/10" thick side panel; the suffix 25 means the tube is 21mm×21mm (13/16"×13/16"), fitting the 1/16" thick side panel	90.5"L×3/4"×3/4", or 90.5"L×13/16"×13/16"		
	USKN-DC-WDG-19/ USKN-DC-WDG-21	Wedge Bottom Beam, for the fluid draining kennel floor, the angle could be customized according to customers' requests.		90.5"L×3/4"×3/4", or 90.5"L×13/16"×13/16"		
	USKN-FDB-150	Stainless Steel Feeding Bowl, 6"(15cm) diameter and 3.5"(9cm)deep		Φ6"×3.5"		
	USKN-FDSB-002	Stainless Steel Feeding Bowl Mount, fitting two 6" feeding bowls and hooks attaching to the swivel door at the bottom of kennel front door			USKN-FDSB-002	
Spare Parts	USKN-STB-1504	Goods Storage Bucket, made from stainless steel wire, hook up to side panel for private belongings.		15"L×4"W		
	USKN-TPGP-45	Top Grid, prevents dogs escaping from kennels, made from stainless steel bars		45.6"L×22.4"W	USKN-STB-1504	
	USKN-DGR-SS-36	Stainless Steel Rest Bed, 36"		35.4"L×23.6"W×3.9"H		
	USKN-DGR-FB-32	Textilene Fabric Rest Bed,32"		32"L×26"W×3.9"H		
	USKN-NMHD-3/5	Name Card Holder, 5"×3"		3"×5"		
	USKN-NMHD-5/8	Name Card Holder, 8"×5"		5"×8"		
	USKN-DTCV-45	Draining Ditch Cover, made from stainless steel panel		45"×7.7"×5.3"		
USKN-FRD-LTC-SS	Double locking front door			USKN-FRD-LTC-SS		
USKN-FRD-LTC-MAG	Magnum style magnetic door latch			USKN-FRD-LTC-SS		
USKN-HDPE-CLP	Stainless Steel screw kit for HDPE panel installation			USKN-FRD-LTC-SS		
USKN-INSTL-CLP-S	Connector for side panel and front door or top grid			USKN-INSTL-CLP-S		
USKN-GLS-CLP-10	Clamp for tempered glass installation, good for 1/3"(8mm)			USKN-INSTL-CLP-S		
USKN-GLS-CLP-06	Clamp for tempered glass installation, good for 1/3"(8mm)			USKN-INSTL-CLP-S		
USKN-FDSB-LTC	Latch for the swivel door at the bottom of kennel front door			USKN-INSTL-CLP-S		

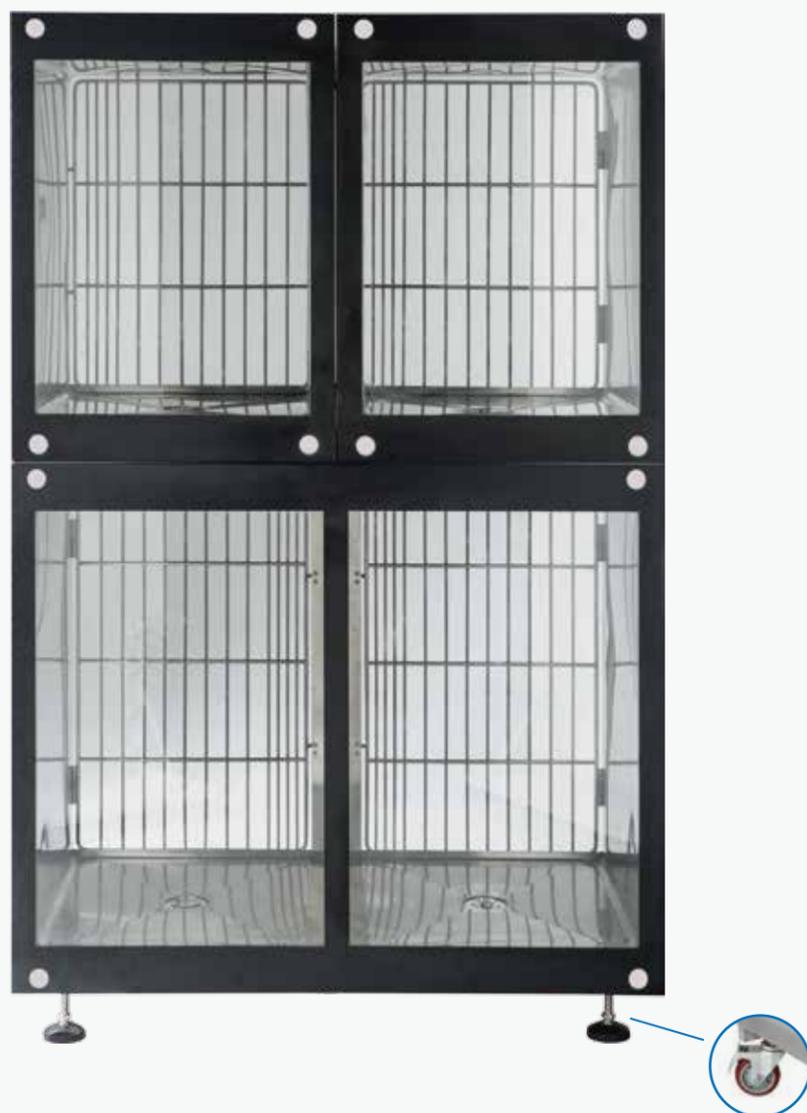

Optional casters available for better mobility

Stainless Steel Display Cage | KA-509D

Maximum visibility into the cages.

Thick, safe tempered glass, ensures incredibly long service life without any animal paw scratching.

Round cornered, fully welded cage; easy to clean and keep hygienic. Allows placement of cage units in front of existing windows. Excellent for young animal sales and adoption areas.


Custom configuration welcome. We can create the perfect single, double or triple stacked display and boarding options for your space. Optional parts available, including illumination LED light, central drainage on floor, water sprayer for cleaning, heated floor, ventilation window, etc.



The back door for staff to operate



Front view of a 2M+1L bank

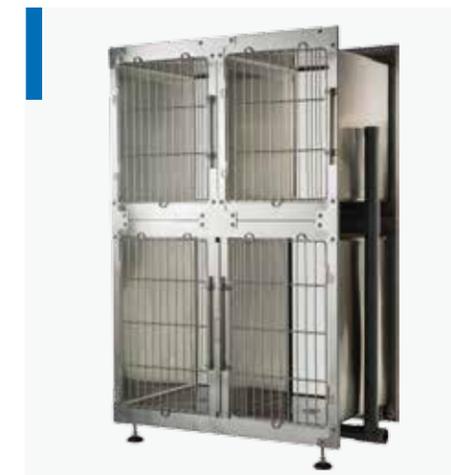

Back view of a 2M+1L bank



Optional LED light



Optional central drainage on floor



Secured double locking door latch system

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension		Volume (cbm)
			Stand ^(cm) (W×D×H)	Package	
KA-509D-M	22.0	26.0	59×66×80	65×85×71H	0.39
KA-509D-L	45.0	50.0	118×95×66	123×100×71H	0.87



Optional caster is available

Professional Stainless Steel Modular Cage with Rounded Corner | KA-509

The Best of Best Professional S/S Modular Cage

Rounded corners and solid walls eliminate areas where germs can grow, leading to a better and more hygienic environment; easy to maintain and keep the space sterilized and clean.

Top grade SUS304 stainless steel material provides a long service life. Single hand locking system makes it easy to handle a dog and lock the door at the same time.

Strong front door made of 8mm and 6mm stainless steel rods makes for a safe and secure area; door options are available (i.e. transparent doors for therapy purposes).



2L+4M+5S Bank




Magnetic Latch



Removable Door Hinge



T-joint Connectors

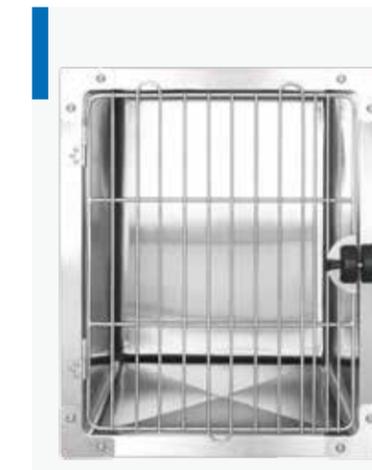


Professional Polishing and Finishing



Stack on Firmly and Securely



Optional: Drain Hole

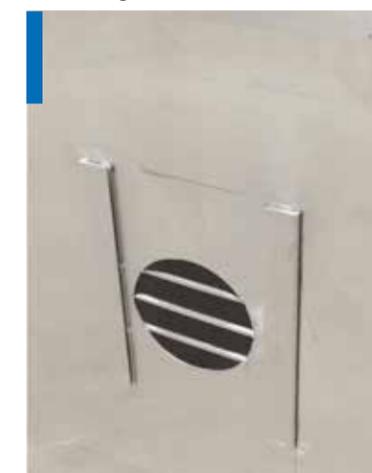

Optional: Air Vent



DRY BED is optional which is the best mate for this cage system.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
KA-509S	13.5	16.5	47.2×45×53	18-1/2×17-23/32×20-55/64	67×51×49	0.14
KA-509M	22.0	26.0	59×66×73	23-15/64×25-63/64×28-47/64	77×63×68	0.33
KA-509L	45.5	50.0	118×66×88	46-29/64×25-63/64×34-41/64	122×92×68	0.76
KA-509M-TH	22.0	27.0	59×66×73	23-15/64×25-63/64×38-47/64	77×63×68	0.33



KA-509-WET

Stainless steel waste tray for easy cleaning purpose
Fully welded floor grid, smooth and comfortable
Polished radius corners all around, easy to keep the whole compartment hygienic
Magnetic latch designed for noise dampening closing operation.




Stainless Steel Waste Tray



Stainless Steel Floor Grid

KA-509-DRN

Special drainage system removes your stress of cage cleaning in day to day maintenance.
65mm ultra big pipe ensures fluent waste fluid draining.
Based on classic KA-509 cage system with radius corners, the cage is super easy to be cleaned by water stream.
Sizes can be tailored according to actual space and practical needs.

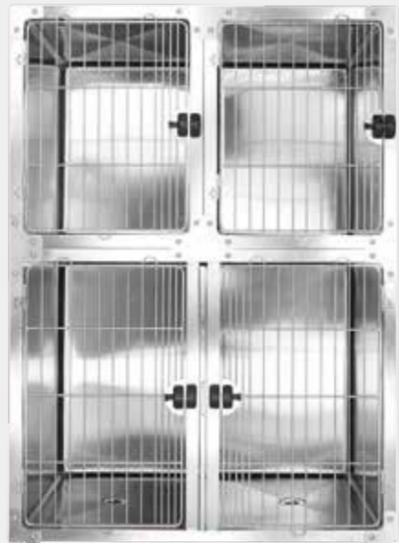

Actual Application





Optional caster is available

Stainless Steel Professional Modular Cage System


KA-508

SUS304 stainless steel material, bite proof and vibration proof, very sturdy material and construction.

Our smartly designed, strong, stainless steel joint connection makes for easy assembly.

A standard bank consists of four small cages, three medium cages and two large cages assembled into a stacked grouping.

Totally sealed structure eliminates droplet infection; the best choice for vet clinics and kennels.

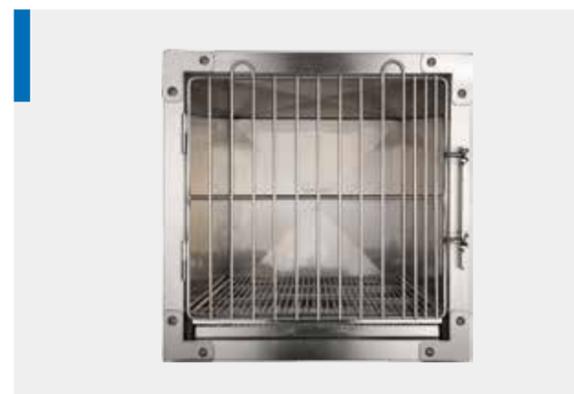

KA-508S



KA-508L



Pull-out Stainless Steel Waste Tray



Varied Connectors for different configurations




Fine Stainless Steel Floor Grid



The 750cm long custom medium size and 1500cm long custom oversize combination Display (KA-508-750 & KA-508-1500)

ITEM	Net Weight	Gross Weight	Dimension			Volume (cbm)
	(kg)	(kg)	Stand _(W×D×H) (cm)	Stand _(W×D×H) (inch)	Package _(L×W×H) (cm)	
KA-508S	24.5	30.0	60×72×60H	23-5/8×28-11/32×23-5/8	64×76×64	0.31
KA-508M	38.2	43.0	80×72×70H	31-1/2×28-11/32×27-9/16	83×73×76	0.47
KA-508L	43.6	48.0	120×72×80H	47-1/4×28-11/32×35-7/16	123×76×93	0.89
KA-508-1500	54.5	77.0	150×72×70H	59-1/16×28-11/32×35-7/16	154×76×74	1.22
KA-508-750	34.5	54.0	75×70×72H	29-17/32×28-11/32×27-9/16	79×76×74	0.50

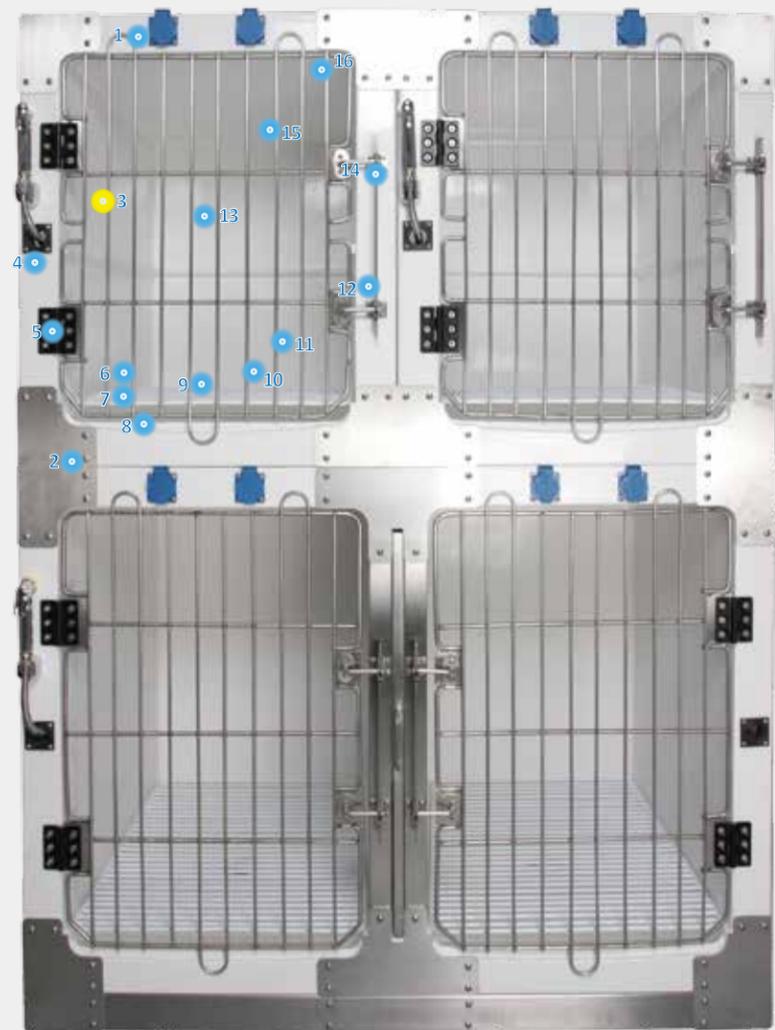

● Front ● Back

Fiberglass Modular Cage | KA-510

The newest addition to our cage line, KA-510 fiberglass cage, has been developed to look great, work functionally, and outlast the average modular cages in the market.

The fiberglass unit is constructed from smooth molds that provide added strength, resistant to scratching, easy to clean and provides a healthy living environment for the animals.

All our various sized fiberglass cages are fully assembled and ready for years of use. They are shipped as one piece.

Fiberglass is easy to clean and maintain. Less time and labor involved in maintenance.

5cm (2") large radius round corners and our fluid retaining lip make clean-up easy and allow for a much healthier kennel environment.

No waste trays to remove and clean. Spray water inside and let waste liquid run out through the drainage system. Cleaning has never been this easy before.

Sturdy; Long Lasting; Durable.

Constructed of 4-5mm thick, high strength fiberglass material.

The same thickness as fiberglass used in Yacht construction ensures extended cage life.

Reinforcing strips have been added on each side to ensure unbeatable strength.

Strong and Secure Door Locking System.


Comes standard with dual point self-latching doors; the narrow tolerance spacing between the vertical rod and latch opening makes these extremely quiet (even when dogs are pawing against the door).

8mm diameter frame, 6mm vertical rods and 8mm cross members; together create a safe bank of cages.

Very Functional, A Complete Working Station.

Pre-installed spray head and hose, drainage system with filter and power supply sockets, make it not merely a cage but a multi-functional working station.

Anesthesia machine, oxygen machine, heating panel, etc.; all electric appliances can be plugged into the power sockets and hung up to horizontal rods in the middle of the door.

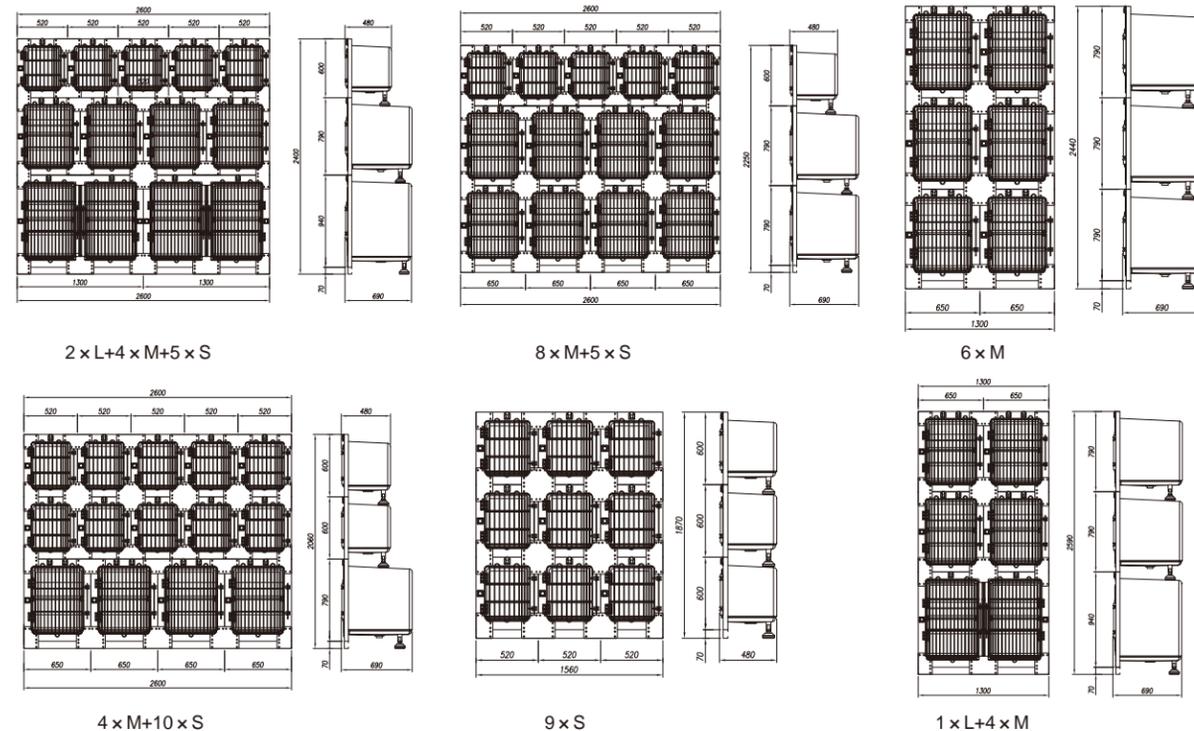

The Unique Features Which Make KA-510 Fiberglass Cage The Ideal Choice for Both Groomers and Veterinarians



1 Isolated wire pipes and water pipes, installed on the outside of the cage ensures a safe situation. Several power sockets pre-installed so all the electric appliances could be plugged directly into the cage banks.



2 Varied one-piece stainless steel connection board, in L, T or cross shape.



3 Reinforcing strips built in and an extremely smooth mold ensures rigidity on each wall.




4 Pre-drilled hole for water sprayer head and hose.



5 Removable front door affixed by strong zinc die cast hinges could be replaced by a transparent acrylic door and transform the cage into an oxygen therapy station.



6 65mm diameter waste pipe and 100mm main pipe will not get blocked by hair or dirt.

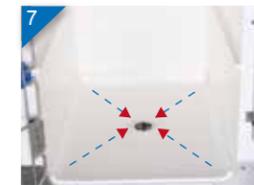

7 Basin shaped base, 10mm lower in center than each corner, ensures all waste liquid drain fluidly.



8 Fluid retaining lip, ensures no leakage from the front lip.



9 Built-in chrome plated bronze drainage system with fine gauge filter and stainless steel capped rubber plug.



10 Optional fine grid on the bottom, to ensure liquid drains away from the animal.



11 5 cm radius round corners makes cleaning easier than ever.



12 Dual point self-latching door system.



13 Step lips on 3 sides could divide the cage into two compartments with a grid to contain more dogs, or change it to a cat cage with a 1/4" grid rest shelf.



14 Narrow tolerance between latch opening and the rod on door, make the door really quiet even dogs are pawing at it.



15 Pre-threaded on 8mm thick ultra sturdy fiberglass, will not cross thread or strip out.



16 Nicely welded and polished round cornered SUS304 stainless steel front door.



Fiberglass Medium Size Oxygen Therapy Cage

KA-510M-TH

An optional acrylic door easily transforms the unit into a therapy cage. Thermometer and hygrometer preinstalled on the door, for accurate temperature and humidity monitoring.



Ventilation Holes Predrilled on the Acrylic Door



Classic Double Locking Latch




Pre-drilled Hole For Water Sprayer Head and Hose




Precision Thermo-Hygrometer




KA-510-2M-HRZ

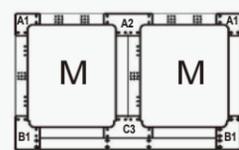


KA-510-4M-BK

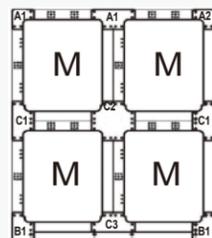


KA-510-2M-VER

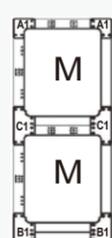

KA-510-1L2M-BK

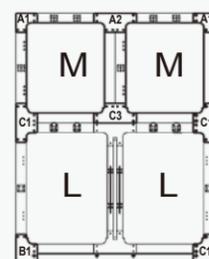

KA-510-2L-BK

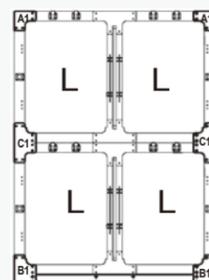

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume (cbm)
	(kg)	(kg)	Stand(W×D×H) ^(cm)	Package	
KA-510S	18.0	21.0	52×48×60H	56×77×64H	0.20
KA-510M	25.0	29.0	65×69×79H	69×77×83H	0.44
KA-510L	60.0	65.0	130×69×94H	134×77×98H	1.01
KA-510M-TH	22.0	27.0	65×69×79H	69×77×83H	0.44



Stainless Steel Medium Size Oxygen Therapy Cage

KA-509M-TH

Strong, durable 8mm thick acrylic front door fixed on 8mm diameter stainless steel frame. Very important facility for animals with cardiovascular or respiratory disease. This proven design safely, reliably, and conveniently keep pets under the oxygen condition you prescribed. Adjust the flow meter to allow oxygen into the unit. Ventilation holes predrilled on both front door and back board eliminates worry about pet inside. Thermometer and hygrometer installed on front door, display the humidity and temperature in the cage accurately. Stainless steel grid and waste tray are available.



Classic Single Hand Lock



High-precision Thermometer and Hygrometer



Adjustable Air Vent (Front and Rear)





Door unit (door and frame) is made of SUS 304 Stainless Steel.
Other parts are made of SUS 201 Stainless Steel

Stainless Steel(Hybrid) Cage | KA-505-SS

Cost-effective Stainless Steel Modular Cage Bank System. Economical and Functional.

Specific Functions of High-end Modular Cage Bank: Classic single-hand double-point lock; Dividable Large-size Cage Unit (divided into 2 medium-size cage units); finer and narrower spaced wires at latch area (preventing pets from opening the latch by extending their paws).

Optional mobile platforms (with a bottom frame and casters).

Medium-size platform is used for 1L+2M; Large-size platform is used for 2L.


Multiple combinations to choose: 5S+4M+2L, 4M+2L, 2M+1L, 3S+3S+3S.

Flat packed, easy to reduce transportation costs.

Perfect design, quick installation and stable structure.



Magnetic Latch

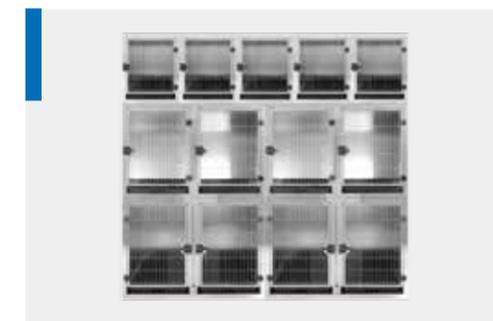

2L+4M+5S Set



Fine plastic dip-in treated floor grid



Increased density wire at latch parts



Dividable large cage



Double layered plastic waste tray



Coin-driven screws



Large-size platform

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension		Volume (cbm)
			Stand ^(cm) (W×D×H)	Package	
KA-505S-SS	24.5	26.0	44×50×56H	58×52×17H	0.06
KA-505M-SS	38.2	40.5	55×66×75H	67×76×16H	0.11
KA-505L-SS	43.6	46.0	110×66×80H	144×82×18H	0.19



Optional mobile platforms (with a bottom frame and casters)

Flat Packing Professional Modular Cage with Solid Walls

KA-505

This professional modular cage system includes everything you need (pans, grates and cage connectors) to build a kennel area unique to your space and needs. Specially engineered and designed for safety, strength, durability and longevity. Unique auto-closing latch system for ease and convenience of operation. Plastic dip coated, narrow tolerance heavy duty floor grid; provides a safe and long lasting surface. Add mobility and versatility with a frame and large cage dividers. Cages are made of heavy-duty 20 gauge zinc plated steel with electrostatically powder coated exterior and interior. Off-white color matches any decor and always look attractive. Removable doors, dividers, floor grates and pans allow for easy cleaning. Mobile platforms are the perfect solution for when kennel banks need to be moved or rearranged. They are equipped with 4" heavy duty zinc plated casters which rotate 360° for easy maneuverability. The front wheels lock, allowing the assembled kennel to remain stationary. Plastic liquid guide plates have been added to prevent liquid waste leakage along the side walls. The tray is formed from plastic blow mold. Nice polished, no sharp edge at all.



Magnetic Latch



Rubber seal, plastic treated grid and plastic tray



Coin-driven screws



Dividable large cage



2L+4M+5S Set



Optional mobile platforms (with a bottom frame and casters)

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension		Volume (cbm)
			Stand ^(cm) _(W×D×H)	Package	
KA-505-S	13.6	15.6	44×56×50H	60×55×17H	0.056
KA-505-M	19.5	21.5	55×66×75H	78×69×18H	0.097
KA-505-L	28.3	31.5	110×66×80H	119×87×18H	0.186




New top LED lighting system (optional).

Stainless Steel Display Cage | KA-505D

Stainless Steel Modular Display Cage with Solid Walls.
 New vision, wonderful vision!
 Maximized visibility into this display cage system.
 Thick, noise dampening and safe tempered glass, ensures incredibly long service life without any animal paw scratching concern.
 A very good solution for puppies and kittens display at front windows.
 Excellent for young animal sales and adoption areas.
 Cost-effective modular cage, easy to assemble and maintain.
 Add mobility and versatility with a frame and large cage dividers.
 Removable doors, dividers, floor grates and pans allow for easy cleaning.
 Polished and finished with no sharp edges.



Back View

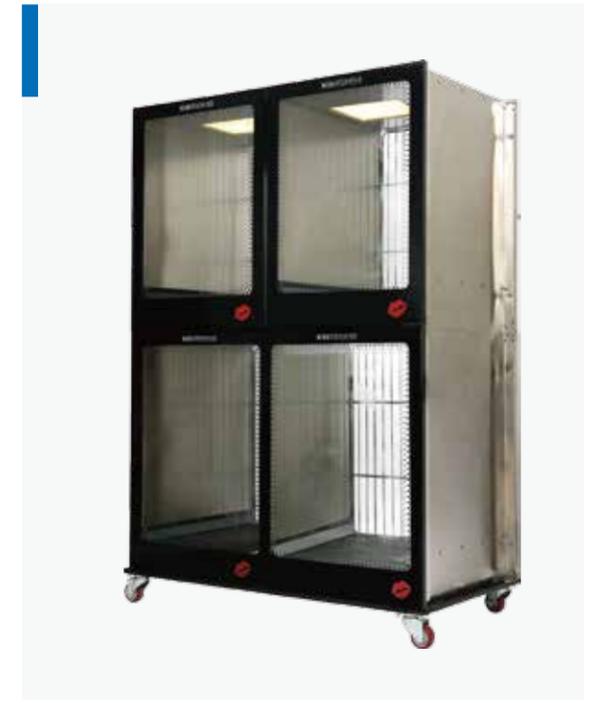

Front View



Double layered plastic pan



Plastic dip-in treated floor grid



Magnetic Latch



Optional wheeled platform



KA-503
Power-coated steel



KA-503SS
Stainless steel

Professional Modular Cage (Wire) | KA-503

These cages are designed to ship Flat Packed, resulting in compact master carton ensuring economical transportation cost.

Double layer ABS waste pan, durable and long lasting.

Unique connection system; easy to install in both vertical and horizontal directions.

Solid divider can be utilized in the middle of both the large cages and between small/medium cages on adjacent small/medium cages.

Stainless Steel double locking latch provides security.

Plastic dip-in treated grid with narrow spacing for pets' comfort and safety.

Wire construction provides good ventilation and unobstructed view.

The design discourages dirt and is easy to keep clean.

Industrial strength coated finish in hygienic hammer stone white pattern.

Removable doors for easy cleaning, with finer, narrower spaced wires at latch area, prevent pets from opening the latch by extending their paws.

Multitudes of modular combinations available; perfect fit for any and all your possible spatial arrangements.

Optional base frame with casters allows for ease of mobility for this fabulous cage. Use 4-caster half frame for one large cage or 6-caster full frame for two large cages.



Horizontal Affixing System



Quality Powder Coating Treatment



Assembly Tools Come Standard In Package



Vertical Connection Between Up-down Decks



Plastic Dip Coated Floor Grid



Classic Stainless Steel Latch System

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			Stand _(W×D×H) ^(cm)	Stand _(W×D×H) ^(inch)	Package _(L×W×H) ^(cm)	
KA-503S	17	19	55×58×53H	21-21/32×22-53/64×20-55/64	58×22×58	0.08
KA-503M	22	24	55×71×80H	21-21/32×27-61/64×31-1/2	84×67×23	0.13
KA-503L	41	44	112×71×80H	44-3/32×27-61/64×31-1/2	116×80×24	0.22
KA-503S-SS	17	19	55×58×53H	21-21/32×22-53/64×20-55/64	58×22×58	0.08
KA-503M-SS	22	24	55×71×80H	21-21/32×27-61/64×31-1/2	84×67×23	0.13
KA-503L-SS	41	44	112×71×80H	44-3/32×27-61/64×31-1/2	116×80×24	0.22



Brand New Cat Condo Available Now

CC-1201

More suitable conditions to cats. The left half of the condo is a living area equipped with resting panel; right half is an area for rest and excretion (cat litter boxes available).

Rounded corners; Smooth and easy to clean and keep hygienic.

Clear Stainless Steel mesh and transparent acrylic door.

Easy for users to observe and free of depression for cats.

Open/close the horizontal tunnel to excite and interest cats.

Size customized for Group Customer.

Perfect Performance!





Cat Restraint Cage

|KA-509-RES

Specifically designed for the easy handling of cats during injections, examinations, or anytime absolute control is a must.

Squeeze panel hung on two rolling tracks ensures very smooth back-forth movement. Easy to grip strong holding bars at front of the squeeze panels make the operation fairly easy.

Magnetic latch is the key point for gentle closing/opening; it automatically locks to eliminate the risk of cat escapes.

Built-in little box and waste tray helps reduce the stress of daily cleaning.

Hang-on feeding trays for water and food included.

Two-stage stopping lifting gate ensures safe and comfortable operation on different sizes cats.



Specialized Isolation Cage

|KA-509-ISO

Specially designed reverse sloped floor contains all the infectious liquid in cage; the draining pipe on the back guides the fluid away fluently.

Semi-sealed acrylic front door cuts the infectious droplet spreading.

Raised stainless steel floor keeps the patient elevated, minimizing the risk of contact with their waste.

Sitting on a raised rack for ergonomic observation and operation.

Standard size at 900mmWx660mmDx750mmH.



S/S Cat Cage

| CC-1301

Space saving solution for cat holding

Slim size maximizes the space usage, versatile combinations realize varied applications.

Horizontal crawling holes pre-drilled, with S/S doors equipped for animal activity management.

Thick rubber edge protection added to make sure there's no sharp edge in the whole system.

Durable front door made from S/S and transparent acrylic materials, ensures good ventilation and reliable strength.

HDPE middle resting platforms included too.





Cat Condo Series

| KC-801

Wonderful Fish shaped decorative patterns realized on this fantastic cat condo. Luxury vertical living space on top, with scale shaped ventilation holes on front door and playful mezzanine platforms installed in the space. Tall scratching post encourages cats to discover more in vertical direction and have more fun in the compartment.

Specially designed sunken toilet and sleeping room are the favorite space for many cats; they love to have a little privacy there.

Functional storage cabins at the bottom of whole unit for keeping the personal belongings handily.





New Cat Condo Bank

| KC-702

UVC sterilization light and illumination light installed.

Newly invented natural ventilation roof equipped, it keeps air refreshed, powered by the heat generated by animals inside.

Classic central playground in the middle of the whole unit, offering great social and playing space for cats, encouraging them to discover more between compartments.

Noise dampening safe magnetic latches on the front doors ensure quiet and secured door operation.





Optional mobile platforms (with a bottom frame and casters)

Brand new model in white is coming

Cat Condo

KC-701 Series

A wonderful, specifically designed feline living habitat with playful and functional space arrangement. Spacious living area with separate litter and sleeping rooms attached. Each space could be managed by the doors from either side. Equipped with terrific hardware, including double locking door latch, stainless steel high polished handle, magnetic crawling tunnel door lock, etc.


Constructed out of lightweight, durable MDF board with scratch-resistant melamine facing. Optional cabinet underneath provides substantial storage. Optional caster base available for easier cleaning operation. Customization available – i.e., tall cat condos for narrow areas – in colors to complement your decor.



Versatile combinations



Customization welcome

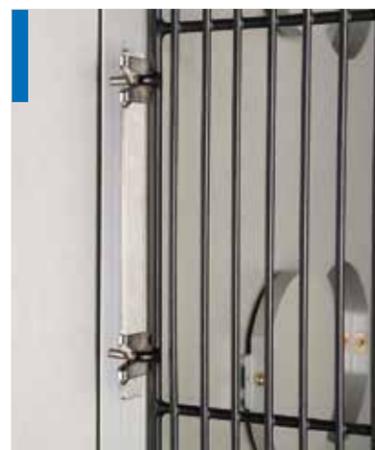

Two types of door latches available S/S double lock or magnetic latch



Horizontal and vertical portals open and close to promote cat-pleasing exploration activities




Separate kitty litter section keeps the living area clean and hygienic



Easy managing portal door system with magnet positioner



Useful storage space for personal belongings

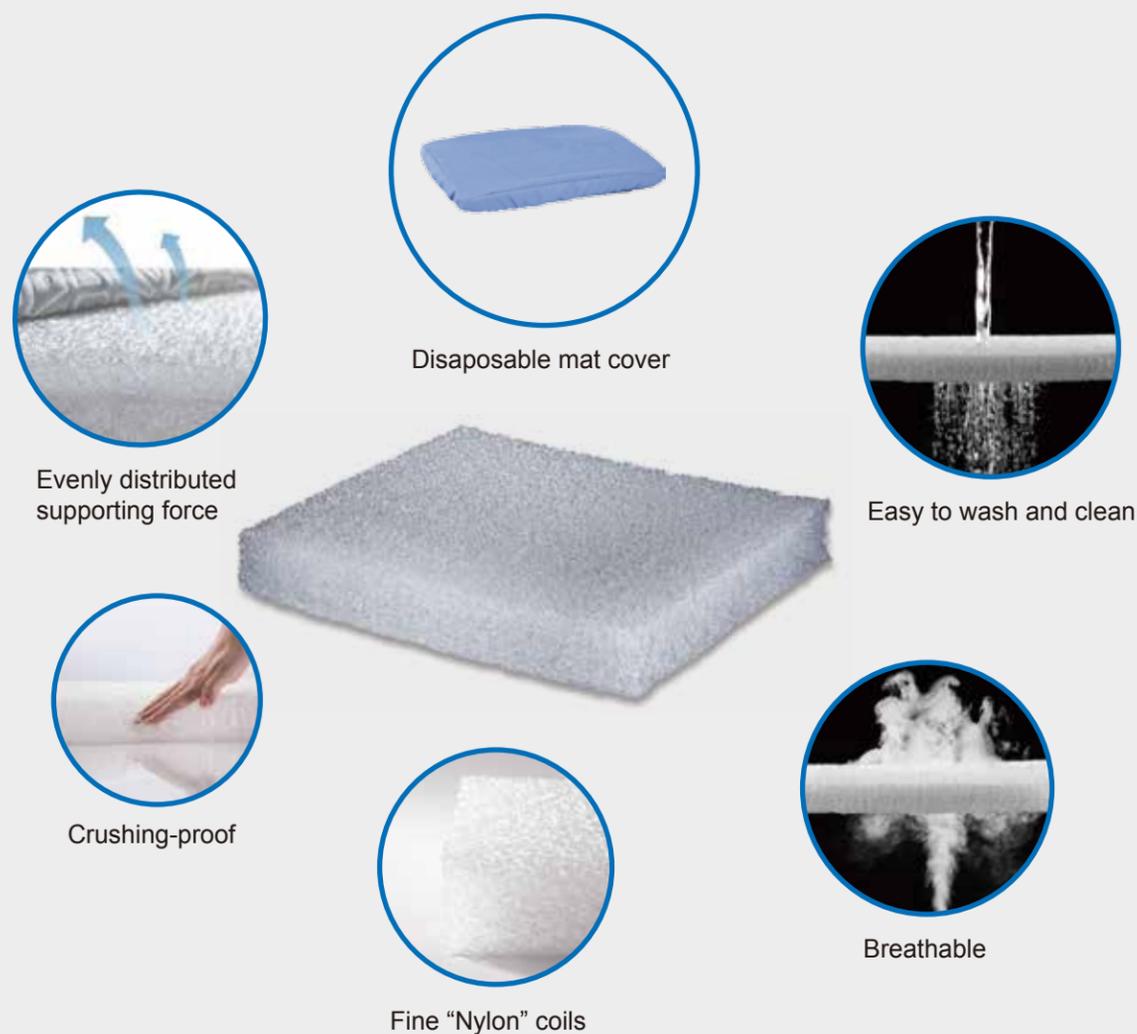

Nylon Mattress

Made of polyolefin fibers.
 Safe, non-toxic, and recyclable.
 Breathable, washable, and easy drying.
 Thousands coil structures provide evenly distributed supporting force to the whole animal body, super comfortable.
 Very elastic, would not be crushed even after over 60,000 times fully pressing or squeezing!
 Warm in winter and cool in summer, a natural atmosphere adjuster for pet users.
 Providing very good ventilation to patient body after surgery, preventing bed-sores and skin softening.
 Certificated zero formaldehyde during manufacturing.
 Little residual moisture thus certificated bacteriostatic efficacy and anti-mite result.
 Perfect elasticity of fiber makes pets have comfortable and healthy sleep.
 The patented fiber has passed BS5852 flameproof test and met European REACH rules.
 Size and shape at your discretion. Sold at an entire roll or fixed length.
 The best mate for Stainless Steel cages



Dry Bed

Super-high quality warm bed, super-density fiber, heat preserving and dampproof.
 Made of unidirectional conducting fiber, which can conduct most of moisture to the bottom of the dry bed and keep the upper surface dry for quite a long time.
 Washable, reusable, durable.
 Size and shape at your discretion. Sold at an entire roll or fixed length.
 The best mate for Stainless Steel cages.

Accessories and customized products

Floor grid (separating the wet from the dry)

Fully welded stainless steel mesh, have no hygienic dead corner, easy to clean and maintain.
 Made of stainless steel rod, comfortable foot feel, protect the health of pets' feet.

Waste Tray

Used to collect the waste liquid of pets' excretion and washing.
 Enclosed weld Stainless Steel, preventing any leakage, easy to take out and clean.

Stainless Steel Tray-cleaning Basin

Designed for cleaning of floor grid and waste tray in the cage.
 Fully welded Stainless Steel Round-cornered, easy to clean and maintain.
 Size customized according to your need;
 Structure design (like water guard, drain, connection between wall and ground, etc) can be made according to actual conditions.

Other optional accessories

Stainless IV pole, name card, feeding bowl, water bowl, etc.

Casters and base frame

Casters work with KA-509 and KA-508 series.
 Base frames with casters work with KA-505, KA-510, KC-701 series.
 Casters can be installed in the pre-drilled threaded holes on the bottom.
 These cages without pre-drilled holes must work with the base frame with casters.





Car Cage

KA-601

A breakthrough design in the portable car crate realm. Great looking, light weight and sturdy construction. Made of 100% aluminum material; light weight for ease of transport. Smart flat packing design is only 7cm thick when folded up. No bolts or nuts required, just easy spread and snap on assembly. Constructed with slim, sturdy aluminum tubing, yields an attractive, strong rigid cage. Fully sealed bottom base prevents waste fluid running off and contains waste products.

ITEM	Dimension		Volume (cbm)
	W×D×H(cm)	Package	
KA-601S	52×52×52H	56×52×10	0.03
KA-601M	77×43×51H	78×52×10	0.04
KA-601L	88×51×58H	90×60×10	0.06
KA-601XL	100×60×65H	100×68×10	0.07



Clamp holding top



Double secured front door latches



Fully sealed base



Easy to fold



Travel Pet Kennel

VK-CD Series

Portable pet carrier for air travel, car travel or visit to the vet for small dogs, cats, rabbits, hedgehogs, chinchillas, piggies, etc. Easy to stack up with pre-molded matching feet on both upper and lower sides of different sized compartments. Assembling with clips, no screws or tools needed, very convenient design. Super strong structure, anti-crushing even with heavy goods loaded on top. Ventilation holes on 4 sides with quick-latching front door. Stable structure, safe latching, ergonomic handle, stylish shape, excellent workmanship, and abrasion resistance. Pets have a comfortable and relaxing trip in the kennel.

ITEM	Dimension		Volume (cbm)
	W×D×H(cm)	Package	
VK-CD-03	67.5×51×47H	68×51×47H	0.053
VK-CD-04	80×56×58.5H	80×56×59H	0.087
VK-CD-05	100×67×75H	100×67×65H	0.150



Crate with door open



Safe latching















Side view



Top handle



Detailed Dimension Drawing of Cage Bank

Unit : mm

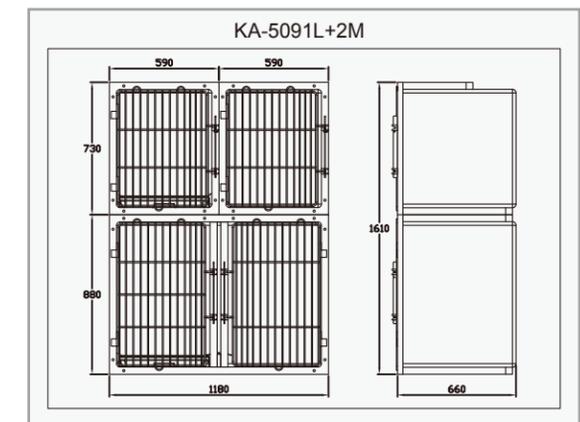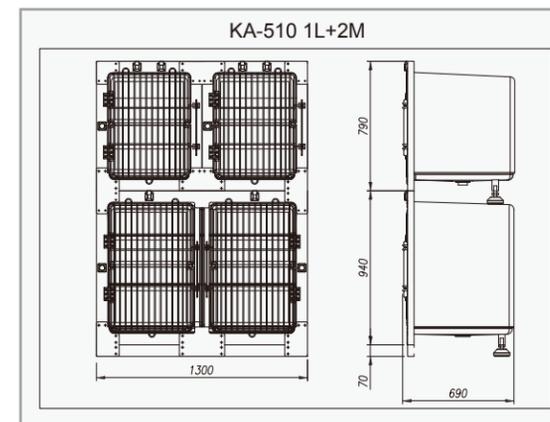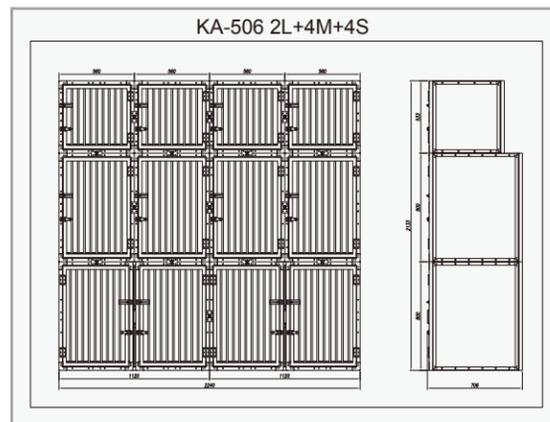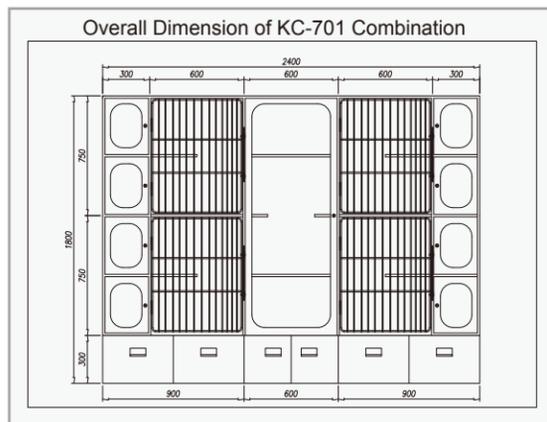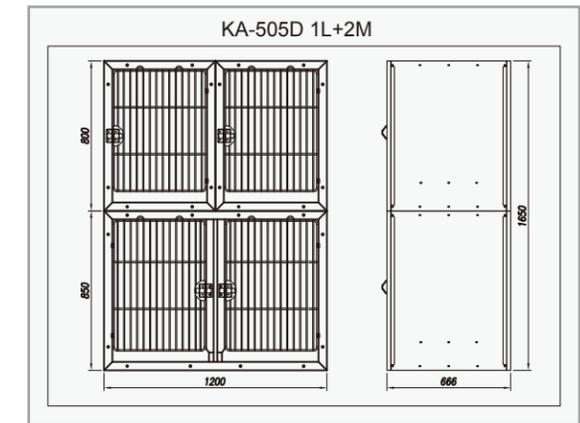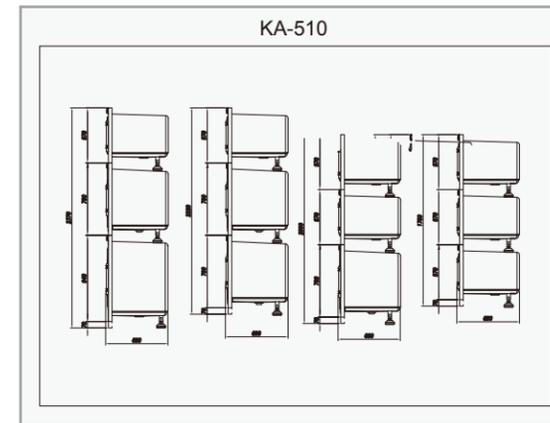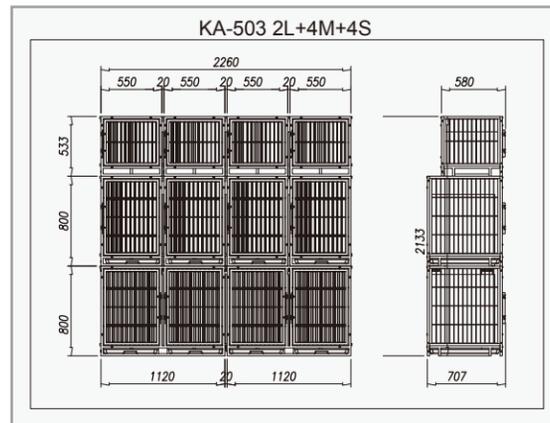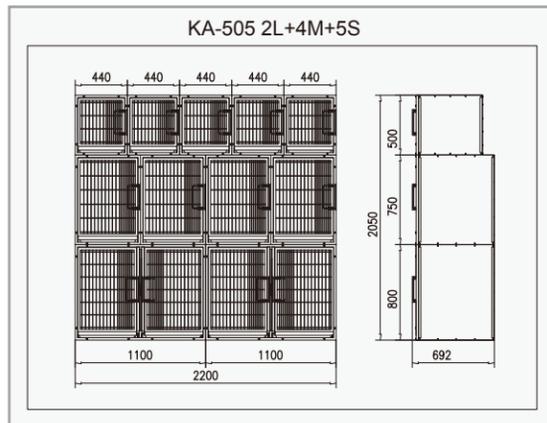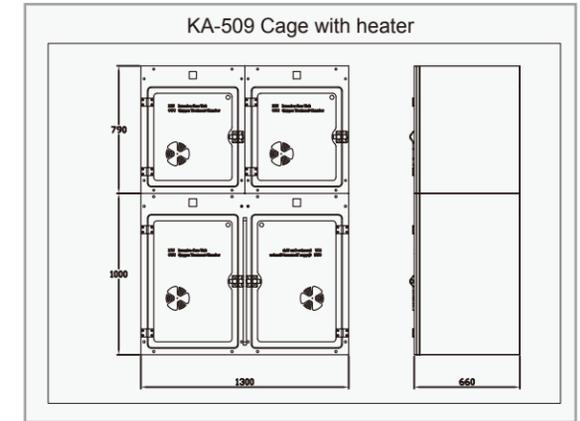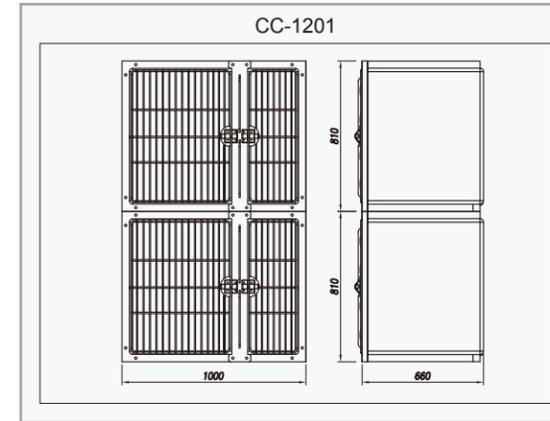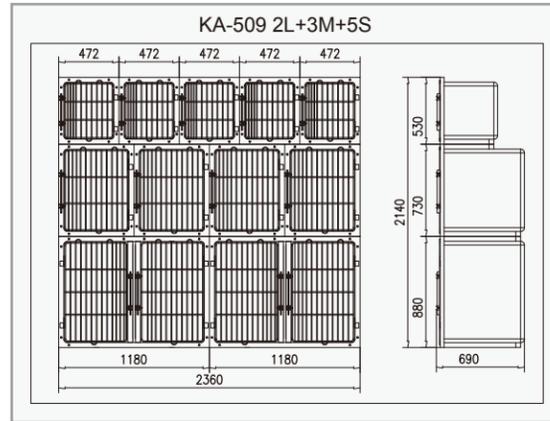
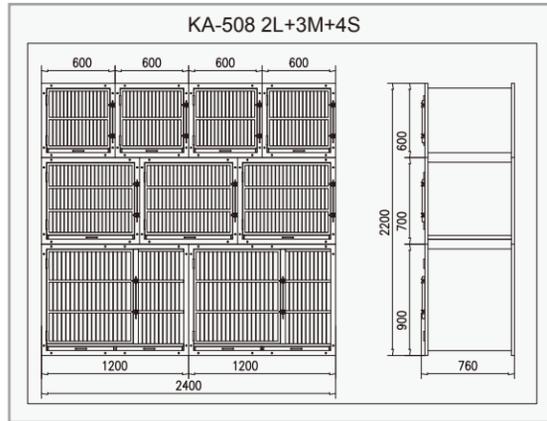


Unit : mm

Detailed Dimension Drawing of Cage Bank

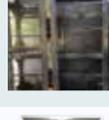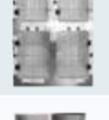Feature	Material		Customers' demand							Customization		
Model	Stainless steel	Other material	Pet store	Pet hospital	Pet show	Pet adoption	Kennel and breeder	Special for cats	Transit pet	Size	Advertising window	Optional drainage system
 Kennel	☆					☆	☆	☆		☆		
 KA-509D	☆		☆			☆				☆	★	☆
 KA-509	☆		☆	☆		☆	☆			☆		☆
 KA-508	☆		☆			☆	☆			☆	☆	☆
 KA-508C	☆		☆	☆		☆			☆	☆	☆	☆
 CC-1201	☆		☆	☆		☆	☆		☆	☆	☆	☆
 KA-510		Nautical grade fiberglass		☆			☆					★
 KA-506		Light weight Aluminum composite board	☆		☆		☆					
 KA-505-SS	☆		☆	☆	☆		☆					
 KA-505		Quality power-coated steel	☆		☆		☆					
 KA-503		Quality power-coated steel	☆			☆	☆		☆			
 KA-701		Environmental friendly multilayer HDF	☆	☆			☆		☆	☆		

Optional accessories: ☆ shows that the accessories need extra payment. ★ shows that the accessories are standard ones without extra payment.

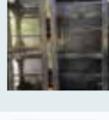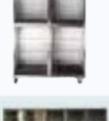Feature	Optional accessories							Package and Transportation		
Model	YangCang cage	Floor grid and tray	Medical record card and name card	IV pole	LED pole	casters	Frame and casters	Carton	Crate	Platform truck
 Kennel			☆						☆	☆
 KA-509D		☆	☆		☆		☆	☆		
 KA-509	☆	☆	☆	☆	☆		☆	☆		
 KA-508		★	☆	☆	☆		☆	☆		
 KA-508C			☆	☆	☆		☆	☆		
 CC-1201			☆	☆	☆		☆	☆		
 KA-510	☆	★	☆	☆		☆		☆		
 KA-506		★				☆		☆		☆
 KA-505-SS		★		☆		☆		☆		☆
 KA-505		★				☆		☆		☆
 KA-503		★				☆		☆		☆
 KA-701						☆			☆	

Optional accessories: ☆ shows that the accessories need extra payment. ★ shows that the accessories are standard ones without extra payment.



ACE Deluxe Grooming Table | FT-889

The new flagship of our grooming table fleet.

A real all-in-one design, with great looks and fabulous functions.

LED light top provides soft, effective illumination. The safety glass top is designed to be antislippery and has a soothing matte finish.

The entire table top revolves 360 degrees for easy and comfortable positioning of animals.

The heavy duty grooming posts adjust up to 48" high and rotates 180 degrees around both ends of the table. Position and lock the arm exactly where you need it. Unlock it with a gentle pulling down motion.

Two swing arms could be connected and formed into an overhead grooming arm system when needed.

Stable and sturdy electric lift column, brings smooth movement of up to 40cm /16 inches height adjustment range.

Practical and helpful power supply receptacles installed on the operation panel, accommodates up to 15A electric appliances.

Transparent sliding drawer keeps your small tools and accessories handy.


Two USB ports included, for convenient charging of your electronic devices (like cell phones).

Hanging hooks on the other side of tabletop keeps all your equipment (nail grinder, clipper and dryers) close and convenient.

ACE ♠



Rotate 360
Degrees Desktop

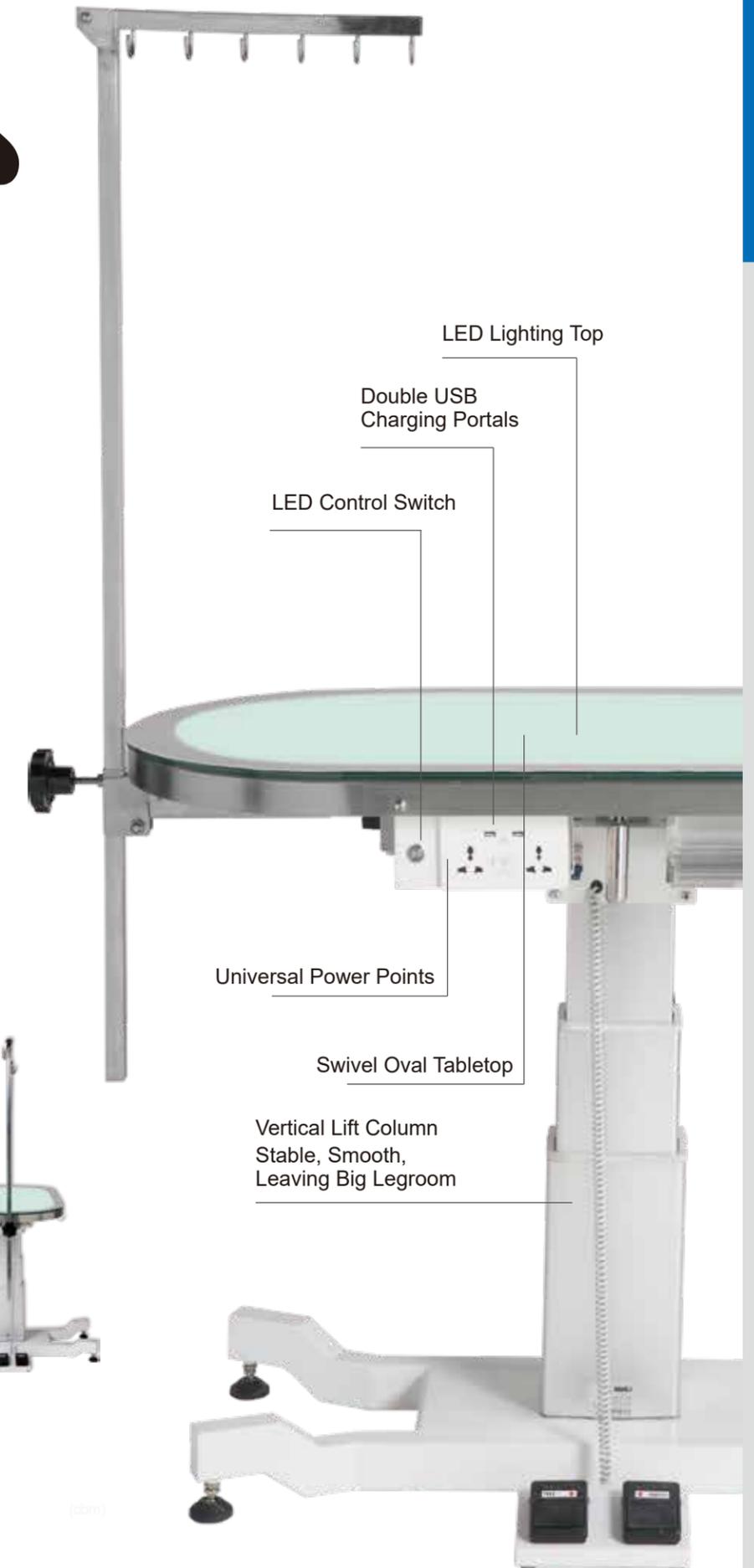

ACE♠



Double swing grooming arms




Reliable rotating top lock



Useful Power Points for up to 15A Power Supply



Two USB charging Ports with tool-organizing drawer

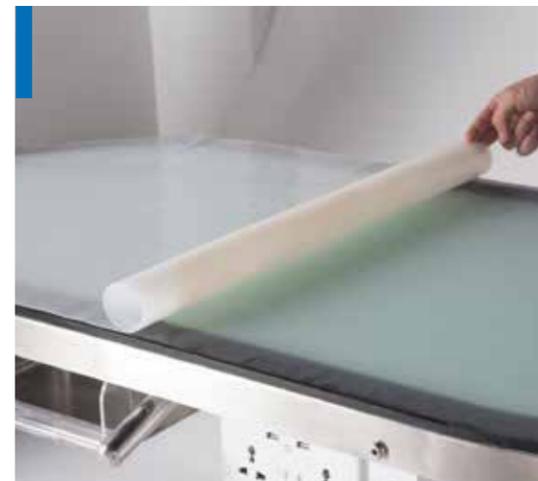

Optional Silicon Mat for Anti-slippery



Flexible Foot Pedal Control

ITEM	Net Weight	Gross Weight	Dimension(cm)		Volume (cbm)
	(kg)	(kg)	Stand(W×D×H)	Package	
FT-889	77	117	120×60×(60-102H)	125×65×74	0.60



Rectangle Electric Pedestal Grooming Table

FT-891

A Real Back Saver!
Great knee space for long time close grooming work.
Super heavy duty grooming arm made from over 1.6" tubular.
Overhead grooming arm system covers over 48" high range, adapting any sized dogs.
Specially designed hooks on the super big grooming arm for easy loop clipping.
Water repellent foot pedal ensures very long lasting service life.
Tool organizer bar, power supply receptacles and handy drawers all included.



Rectangle Pedestal Table with Stainless Steel Top

FT-892

Unbreakable stainless steel table top, good for very long termed service life.
Equipped with designated designed tabletop mat for pet grooming, it does not catch loose hair, has wonderful water repellent and anti-slippery features. The printed patterns are customizable too.
Lifted by Linak high quality column ensures very good stability and smooth running, with a 200kg+ loading capacity too.
Grooming arm and restraint loops included.
Tool organizer bar, power supply receptacles and handy drawers all included.



Oval Top Electric Pedestal Grooming Table

FT-891-OV

Stable Lifting column finds your most ergonomic height without any stress.
Oval shaped top offers great working convenience and save precious surrounding space.
All Pro accessories included, like tool organizer bar, power supply receptacle and handy drawer.
World top class pedestal column by industry leader Linak.
Water repellent foot pedal set gives great durability.
Removable super heavy duty grooming arm with over 48" height adjustment, adapting any sized dogs.



LED Lighting Pedestal Lifting Table

FT-893

Fancy table is the best mate for your fashionable grooming salon!
LED table top not merely provides great appearance, but also makes all details more visible especially on animals' abdominal portions.
Newly invented special lifting base equipped with heavy duty castors realizes great mobility.
All Pro functional parts included in this wonderful table, with power supply receptacle, heavy duty grooming arm, handy drawer and tool organizing bar.



FT-838

FT-838REC


Air Lift Grooming Table | FT-838 Series

60cm (23-1/2") diameter table top swivels 360° smoothly on central shaft, making it possible to view the dog from any angle. Easy-to-Control air lift base offers a wide adjustable operation range from 80cm (31-1/2") -120cm (47-1/4"). You will be able to find a comfortable height no matter how tall either you or the dog is.

Smartly designed metal base with scratch resistant white veined powder coating offers secure stability of operation and durable product wear.

Comes equipped with a set of aluminum clamp and stainless steel dual purpose grooming arms. 4 colors available for the table top rubber mat.

FT-838 has round tabletops with two diameters, either 60 or 70cm; and FT-838REC are referring to the versions with rectangular table, at either 60×40cm or 75×45cm.

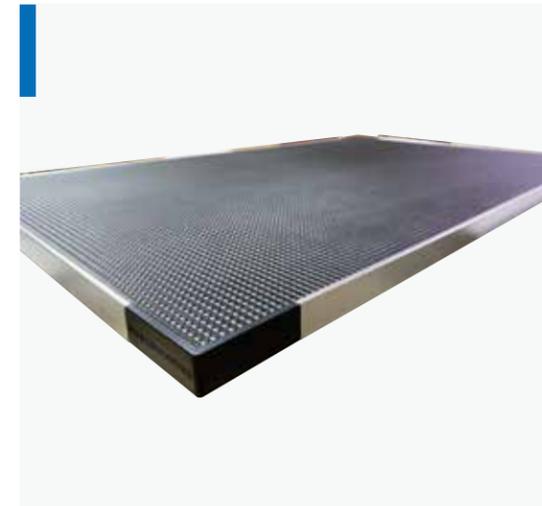

Aluminum Clamp



Rubber Tabletop W/ Metal Frame



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)		Volume (cbm)
			Stand ^(cm) (W×D×H)	Package	
FT-838	21.5	23.5	Φ60×(68-98)	69×63×26	0.11
FT-838-REC	17.5	19.5	60×40×(68-98)	69×45×26	0.08
FT-838L	23	25	Φ70×(68-98)	71×71×26	0.13
FT-838L-REC	19	21	75×45×(68-98)	47×77×26	0.09



FT-808Pro



FT-808

Super stable electric lifting base




Stainless steel new edition launched



Ultra-low initial height



Earthing foot and breaker with leakage detector

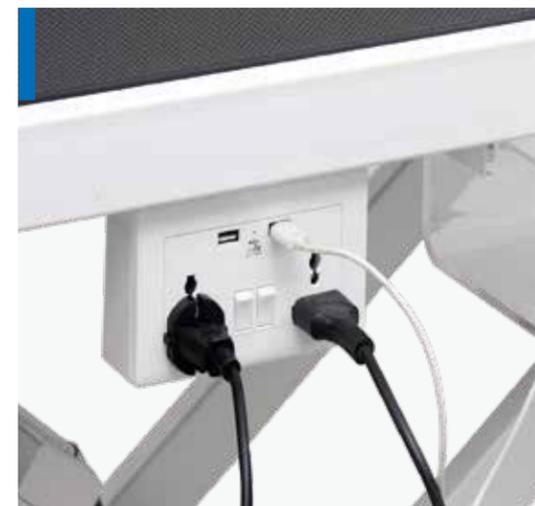

Power receptacle with USB port



Storage drawer

Low-Low Electric Lifting Table | FT-808/808Pro

Keeps the legend going!
After 8 years successful selling in pet grooming industry, our low-low table has won great honor from the global market. Now it welcomes its moment to get the whole system re-engineered, to launch a more powerful, functional and reliable table, and to keep writing the legend!

Practical and helpful power supply receptacles installed on the operation panel, accommodates up to 15A (220V) electric appliances.
Sliding drawer keeps your small tools and accessories handy.
Two USB ports included, for convenient charging of your electronic devices (like cell phones).
Hanging hooks on the other side of tabletop keeps all your equipment (nail grinder, clipper and dryers) close and convenient.
With adjustable leveler feet added.

The whole table is earthed by a metal foot leveler.
The power supply system is protected by a top class breaker with leakage detector included too.
Also the power cable itself is earthed too, by the earth prong.
Three layered protection, ensures 100% safe operation on this wonderful piece.

The lowest point of this table is only 22cm (8.5"), providing easy step-on access for pet; including mature, large and infirm animals.
World top class linear actuator provides outstanding loading capacity of 150kg (300lbs), powered by Linak.
Innovative X-style frame is stable, powerful and quiet.
Removable tabletop, easy to clean, enables you the ability to change the color, better suiting your personal style and preference.

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume	
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	(cbm)	
FT-808	64	67	125×65×25-95	49-7/32×25-19/32×9-27/32	37-13/32	129×69×26H	0.23
FT-808Pro	70	74	125×65×25-95	49-7/32×25-19/32×9-27/32	37-13/32	129×69×26H	0.23



FT-801

X-style Electric lifting table

FT-801

Fully covered solid base, leaves hygienic and easy to clean working environment.

Ultra strong single driving system, providing up to 150kg loading capacity and very comfortable legroom under table top.

Super big height adjustable range, from 22cm to 96cm, allows convenient walk-on for most large breeds.

The world top class linear actuating system, Linak, brings a super velvety smooth movement for table lifting up and dropping down, giving the nicest comfort to the animals on top.

Removable tabletop, easy to clean, enables you the ability to change the color, better suiting your personal style and preference.

Practical and helpful power supply receptacles installed on the operation panel accommodates up to 15A (220V) electric appliances.

Transparent sliding drawer keeps your small tools and accessories handy.



Convenient power supply and tool hooks.



Power supply receptacles



Ultra-low initial height



Powerful linear actuator



Earthing foot leveler



Power box with breaker and leakage detector



Electric Mobile Grooming Table | FT-833E

The hottest new table offering is the ideal choice for grooming shows. Wide range height adjustment, from 22cm to 100cm. Unique electric driving system. Four heavy duty casters ensure great mobility and easy moving. Stainless steel matte finish hand control panel for accurate height setting. Large tabletop, accommodating most large breeds.



Ultra big height adjustable range from 22cm to 100cm.



Easy-to-drive pulling bar



Great mobility



Stainless steel manual control



Aluminum clamp and stainless steel grooming arm included



Secure Caster Brake



ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-833E	34	36	120×60×31-103	47-1/4×32-5/8×12-13/64~40-35/64	120×60×30H	0.22



Big Z Electric Lift Table | FT-898 Series

Wide height adjustable range, from 10" to 48". It has a huge loading capacity of 350lbs. Built-in electric receptacles (installed on one end), make for a very convenient power source solution for your nail grinder, dryer and clipper. Four versions available, FT-898 comes with standard flat top, FT-898-R is equipped with a round rotational turntable board in the middle of the rectangular top, FT-898R-SW is the one with fantastic swing arm on one end and also a round rotation plate. All of them could be used with strong overhead grooming arm.




Super Big Lift Range, from 10" to 48"



Simply operating swing arm



Rotary grooming-arm locking system

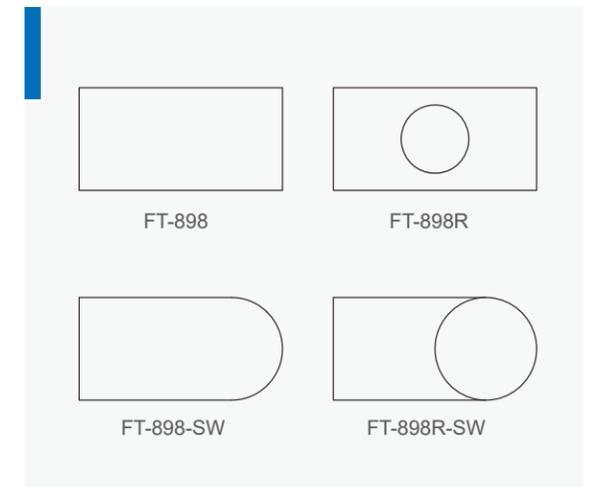

Multiple tabletop configurations

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-898	62	65	120×60×36-108	47-1/4×23- 5/8×14-11/64~42-33/64	123×63×40H	0.31



Accordion Table

FT-899

Extremely stable and sturdy, this table has an exceptionally broad lifting range, from a low of 10" to high of 48". It has a huge load capacity of 280lbs.
 Removable table top, could be cleaned and refreshed easily.
 Welded grooming arm clamps fixed on both ends of top frame ensures the arms are firmly planted in place.
 Overhead grooming arms included.



Ultra-low initial height



Smooth lifting motion



Stable at the maximum height



Water-resistant up-down Foot Pedal



Overhead grooming arm and adjustable leveler feet included



MC-5025

SpongePop Magic Mat Available

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand	W×D×H (inch) Package	
FT-899	64	67	125×65×27-123	49-7/32×25-19/32×10-5/8~48-27/64	129×68×32H	0.28



Single lifting-Column Grooming Table | FT-721

Foldable manual switch control.

Talented invented mechanism allows ultra-high lifting speed of up to 11" (29cm) per minute.

Compact design of the column leaves maxi comfortable legroom for groomers.

Non-slip rubber matted tabletop, scratch-resistant and easy to clean.

Up to 265lbs (120kg) loading capacity.

Grooming arm included.



Lowest height-- highest height



Big loading capacity



Foldable manual control

Aluminum Plastic Compound Table Top

The base board of this new table top is made from aluminum plastic compound. Lighted weight, sturdy, durable, corrosion resistant, easy to keep and maintain. A lighter weighted top always helps lower the center of gravity of the whole table, improving the table stability dramatically too. Aluminum plastic compound table tops are now used on FT-721 and FT-838 series products.





Single Column Electric Lifting Table | FT-721Pro

Slim, stable and smart.

An ultra slim single column mechanism system provides unbeatably comfortable legroom.

Very fast lifting speed with great loading capacity up to 80kg.

Super stable metal base, made from 10mm thick solid metal panel.

All functional Pro series accessories included in the complete table, such as power supply receptacle, tools organizing bar and handy clear drawer.

Fantastic rubberizing coating top, anti-slippery, super durable and very long lasting.



Handy Clear Drawer



Power receptacle



Tools Organizing Bar



Fully Wrapped Rubberizing Coating



Double Column Electric Lifting Table | FT-722

A wonderful double column electric lifting grooming table.

Soft fully wrapping rubberizing coating on top provides nice anti-slippery and comfortable experience for pets.

Ultra high speed linear actuator system allows quick adjustment of height, achieving the most ergonomic point in a second.

Heavy and super stable base made from 10mm metal board.


Optional work station accessories available too, including power receptacle, handy drawer and tools organizing bar.



Super Functional Rubberizing Coating



Foldable Manual Control

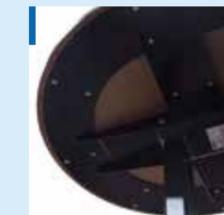

Strong Reinforcement Under Tabletop



Ultra Big Height Adjustment Range



FT-802Pro

Classic Z Lifting Table | FT-802 Series

An economical solution for height adjustable electric grooming table. Slim base leaves comfortable legroom, yet still stable without any wobbling. Two sizes available, 36"×24" for medium, 48"×24" for large; height adjustable ranges 24" - 44". Practical and helpful power supply receptacles installed on the operation panel, accommodates up to 15A (220V) electric appliances. Transparent sliding drawer keeps your small tools and accessories handy. Two USB ports included, for convenient charging of your electronic devices (like cell phones). Hanging hooks on the other side of tabletop keeps all your equipment (nail grinder, clipper and dryers) close and convenient. With adjustable leveler feet added.



FT-802LE

Electric Version



Strong linear actuator



FT-802E Pink/802LE Blue

Color customization welcome



FT-802

Hydraulic version



FT-8101 grooming arm included



Water-resistant up-down foot pedal

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume
	(kg)	(kg)	Stand ^(cm) (W×D×H)	Package	
FT-802E	30	34	60×90×(52-100)	94×64×54H	0.30
FT-802LE	33	39	60×120×(52-100)	124×64×54H	0.40



Japanese Style Grooming Table | FT-875

Japanese Style grooming table

Japanese stylish hydraulic grooming table.

Always comes with a handy tools organizing cart as a combo.

Not merely good for the practical functions, also becomes a decorative part of your grooming room.

Height adjusts from 75cm to 90cm, with lockable rotational table top, comes with grooming arm, clamp and loop.





FT-806

Super Stable Electric Lifting Table | FT-806 Series

Two table options, FT-806 and FT-806S are metal framed with wooden table top; FT-806P comes with a plastic table top. FT-806 has large size frame at 126cm by 65cm; FT-806S is a smaller size at 110cm by 55cm. The best feature of this table is its STABILITY!

Four vertical bars support the table top. When lifted to maximum height, the table is still stable. The center of gravity remains stable and centered, making this the safest electric lifting table in the market.

The minimum height is 36cm (14") when lowered down to bottom.

The height could be adjusted from 36cm to 104cm (14" to 40"), easily satisfying most daily grooming requirements.

Control Unit: Two Pedal Plates (optional manual switch available).



FT-806P



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand	W×D×H (inch)	
FT-806	53.0	57.5	126×65×38-103	49-39/64×25-19/32×14-61/64×40-35/64	129×69×36H	0.32
FT-806S	40.0	43.3	110×55×38.5-100	43-5/16 ×21-21/32×15-5/32×39-3/8	114×59×37H	0.25
FT-806P	42.5	47.5	126×65×39-103	49-39/64×25-19/32×15-23/64×40-35/64	129×69×38H	0.34



FT-804

Adjustable Grooming Table | FT-804 Series

Rock-solid at any adjusted height, from 20 1/4"-40".

Attractive silver toned hammer stone powder coating for durability.

Classic "Z" lifting design eliminates wobbling and tipping.

Chilean made MDF board tabletop.

Pliable PVC edging, resists scratches and mars.

Complete metal-edged frame, durable and long lasting.

Removable top that is easy to clean.

FT-804/804L are hydraulic. FT-804E/804LE are electric. For FT-804E/804LE, we can upgrade them to be Pro version with all the accessories too.



FT-804L-SW



FT-804E



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand	W×D×H (inch)	
FT-804	46.0	49.7	96×65×51.5-101	37-51/64×25-19/32 ×20-9/32-39-49/64	99×68×52H	0.35
FT-804L	48.2	52.7	126×65×51.5-101	49-39/64×25-19/32 ×20-9/32-39-49/64	129×69×49H	0.43
FT-804E	46.0	49.7	96×65×65.5-106	37-51/64×25-19/32 ×25-25/32-41-47/64	99×68×65H	0.43
FT-804LE	49.2	53.7	126×65×65.5-106	49-39/64×25-19/32 ×25-25/32-41-47/64	129×69×64H	0.57



Grooming Table with LED Light | FT-829

This table will be the featured performer, the showcased table in your grooming salon. It is a work of art; the choice of many world class Master Groomers.

The first LED lighting table in the pet industry provides mild and gentle light for your grooming needs.

World class LED chips, made in Germany, ensures long lifespan for lighting panel.

Extremely stable electric lifting base gives you a sense of comfort and ease of operation.

Safety tested matte surface glass table top, silicon sealed edges, clear, skid-proof mat (does not affect or diffuse the light below), makes the surface safe and soothing for pets to stand in comfort.



ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-829	75.50	96.00	120×60×25.5-95	47-1/4×23-5/8×10-3/64×37-13/32	132×71×44H	0.412



Cost Effective Round Hydraulic Table | FT-807


Very affordable hydraulic grooming table.

Lockable swivel table top made from high quality MDF board with black anti-slippery rubber mat attached.

Height adjusted in over 14cm, easily gives the most ergonomic height groomers desire.

With good loading capacity at over 150kg on this table.

Flat packing, economic in transportation and easy to install.

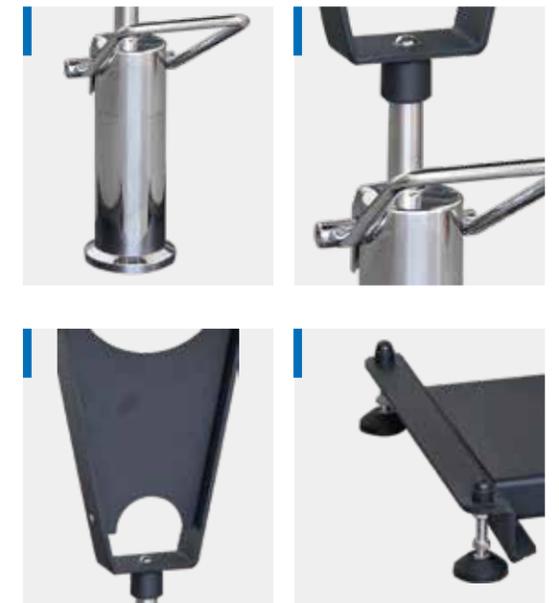



Ultra-Light Weight Competition Table

FT-821 Series

The oxidizing treatment of our custom colors remain vibrant. Numerous color options available, purple, pink, black, and silver.

Skid-proof, comfortable surface texture is soothing to both dogs and handlers.

Ultra-light only 3.7kg (8.2lbs), allowing junior handlers to transport their own ringside tables.

Aluminum, long lasting; the perfect product for dog shows.

Big enough to hold most medium breed dogs while they're being primed for the ring.

Simple two step folding process, pull it apart and it folds down for easy carrying.



Multiple color options



FT-822H



Telescopic version FT-821 is on the market already, which is especially suitable for mobile grooming service.



KB-8107

KB-Clamp

Optional competition grooming arm



FT-821HL
90cm×60cm

Large size available too.



Portable Competition Table

FT-820Series

Wooden top with static proof rubber mat; anti-slip surface. Simple folding method, light weight, easy to pack. Several tabletop colors to choose, black, blue and purple, etc.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)				Volume (cbm)
			W×D×H _(cm)	Stand W×D×H _(inch)	Fold	Package	
FT-820	6.4	7	60×45×71	23-5/8 ×17-23/32×27-61/64	82×45×7H	84×47×9H	0.036
FT-820H	7.2	8	60×45×81	23-5/8 ×17-23/32×31-57/64	89×45×7H	94×47×9H	0.038
FT-820M	8.4	9	75×50×71	29-17/32×19-11/16×27-61/64	89×50×7H	91×52×9H	0.040
FT-820HM	9.3	10	75×50×76	29-17/32×19-11/16×29-59/64	97×50×7H	94×52×9H	0.044
FT-821H	3.7	4.2	60×45×81	23-5/8 ×17-23/32×31-57/64	89×45×7H	92×50×8H	0.036
FT-821HM	4.8	5.3	75×50×76	29-17/32×19-11/16×29-59/64	97×50×7H	92×53×8H	0.039
FT-821HL	6.0	7.3	90×60×78	35-7/16×23-5/8 ×30-45/64	97×50×7H	92×93×7H	0.060



Dog Show Table with Casters

FT-812W

Perfect for shows, travel and overflow grooming. Reinforced fixed brackets under top board ensure extraordinary strength and stability. Four 4" sturdy build-in casters ensure reliable, easy rolling transportation.



Space Saving Folding Table

FT-841

The best solution for generating flexible working space. Each inch counts in the working area. Easy operation with reliable stability.



Useful storage cabinet included.



Round Hydraulic Table with Cabinet

FT-805C

Practical design with storage compartment to help groomers stay organized. Two drawers and one storage cabinet are useful for keeping the tools of the trade (handheld dryers, combs, shears, clippers and any other tools) easily accessible. Hydraulically controlled, adjusts from 78.5cm to 93cm (31" -36-1/2"). Grey cabinet is standard (available in variety of colors).



Plastic Top Folding Grooming Table

FT-712

King of cost effective, high performance grooming tables. Much more durable and stronger than a comparative wooden topped table. Very light weighted, portable, with carrying handle.



Reinforced support system under top board ensures strength and stability.



Round Hydraulic Table

FT-805

This table rotates 360 degrees for easy access to the entire pet, and brakes when placed in your desired position. Blue, Pink, Purple, Green, Grey and Black tops available upon request. Foot pump operated table adjusts from 71.5cm -86cm (28-1/8" -33-7/8") for maximum comfort and grooming efficiency.



Height Adjustable Folding Table

FT-818/819

Exclusive height adjusting system, easy to handle and offers perfect stability. 4 stages of height setting, finding the most ergonomic positions for most groomers. Two sizes available, 120cm×60cm and 90cm×60cm.



Patented locking system, easy to operated and very tight positioning.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)	
			W×D×H(cm)	Stand	W×D×H (inch)		
FT-812W	17	18.7	90×60×76	35-7/16×23-5/8	29-59/64	94×63×16H	0.095
FT-805C	35	38.5	ø70×78.5-93H	27-9/16×30-29/32	36-39/64	71×71×10H	0.12
FT-805	17	18.7	ø70×78.8-93H	27-9/16×28-5/32	33-55/64	71×71×10H 41×40×40H	0.09

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)			
			W×D×H(cm)	Stand	W×D×H (inch)				
FT-712	14.3	16.5	126×65×76H	49-39/64	25-19/32	29-59/64	128×67×14H	0.12	
FT-818	15.80	17.9	90×60×81-93	35-7/16	23-5/8	31-57/64	36-39/64	94×63×16H	0.09
FT-819	20.10	22.1	120×60×81-93	47-1/4	23-5/8	31-57/64	36-39/64	124×63×16H	0.12
FT-841	19.0	21.0	110×60×80	43-5/16	23-5/8	31-1/2	112×62×24H	0.167	



FT-811

Folding Grooming Table with S/S Legs

FT-811/812/813

Classical folding table with non-slip rubber surface.
Stable and durable stainless steel legs.
Dual functional grooming arm.



FT-812



FT-813

Several Tabletop Colors Available



FT-814

Folding Table with Powder Plated Legs

FT-814/815/816

Economical option for folding grooming table.
Square steel tubes bring great loading capacity to this series.
Non-slip rubber mat surface.

Reliable Support Bracket

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)				Volume (cbm)
			W×D×H(cm)	Stand	W×D×H (inch)	Fold	
FT-811	10.7	11.7	75×47×82	29-17/32×18-1/2×32-9/32	83×47×10.5	87×54×13H	0.06
FT-812	13.3	17	90×60×76	35-7/16×23-5/8×29-59/64	90×60×10.5	94×64×13H	0.08
FT-813	14.8	19.4	120×60×68	47-1/4×23-5/8×26-49/64	120×60×10.5	124×64×13H	0.10
FT-814	12	14	75×47×80	29-17/32×18-1/2×31-1/2	83×47×12	79×50×13H	0.05
FT-815	14.4	16	90×60×74	35-7/16×23-5/8×29-9/64	90×60×12	94×63×13H	0.08
FT-816	16.8	19	120×60×65	47-1/4×23-5/8×25-19/32	120×60×12	124×63×13H	0.10

Customized table top printing can be tailor made now!





Table Top Mat

TM-4824-003

The best mate for any grooming tables, giving the best comfort ever to pets. Anti-slippery, anti-fatigue, water repellent and does not catch loose hairs from pets.

Very easy to clean, just simple brushing and dryer blowing will do the job.

High duty rubber mattress at the bottom ensures very good durability.

Versatile top printings available, customization welcome too.



TM-4824



001



002



003



101



TM-3624



001



002



003



101



TM-2418

- 102 
- 103 
- 105 



TM-60RD

- 104 
- 101 



Grooming Hammock

GH-2216


A great helper for nail trimming.

Measures 22"x16" (55cmx40.5cm), good for most medium and small breeds.

Two strong length adjustable bands included for hanging height adaption.

Width adjustable at 4" (10cm) by zipping/unzipping the middle strap to fit different shaped breeds.

Just place dog's legs in holes, wrap around and use the s-hooks to hang. Will hold up to 30 pounds.

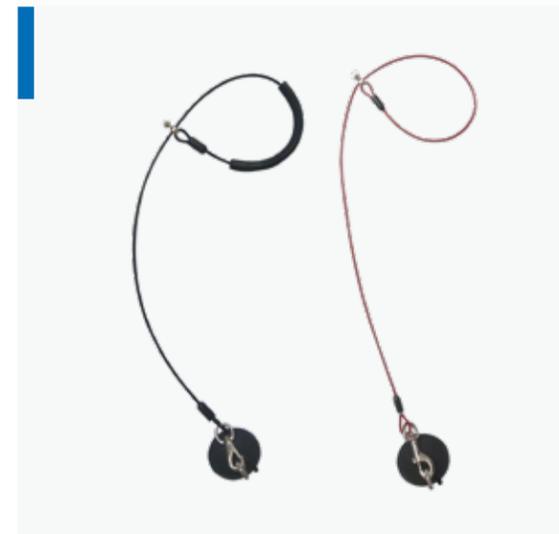

Restraint Cable with Rubber Sucker

SR-10

Industrial grade rubber suction cup holds the cable firmly on bathtub.

Specially designed composite rigging sleeve ensures comfortable animal contacts.

Measures 90cm for big adjustable control range.

PVC coated 4mm sturdy cable, holding up to 500lbs force.

5/8" swivel eye hook build in rubber sucker for added flexibility.

Thumb-turning bolt snap for easy fixing.

SR series wire padlock



SR-01

SR-02

SR-03

SR-04

SR-05

SR-06



Colorful Nylon Loop

NL-LG20





European Style Grooming Arm

FT-8106



Round Rotating Tabletop

FT-830



The Magic Carpet

MC-5025

Brand new product.
 Designed for daily maintenance of large electric lifting grooming table.
 Scratch-resistant material, even the pets' nails cannot cut the surface.
 Composite material inside, absorbent enough to keep extra water during drying from corroding table top.
 Reasonable design: an easy-to-clean hole left for easy usage of grooming arm.



“Birch” Solid Grooming Arm

FT-8105M/L



Dual Functional Grooming Arm

FT-8101



Aluminum Portable Grooming Arm

FT-8107A



Aluminum Clamp

FT-8102



Grooming Harness

FT-8104



















Classical Overhead Grooming Arm
















FT-80402



The “Max Miracle” Super Heavy Duty Grooming Arm

FT-80403

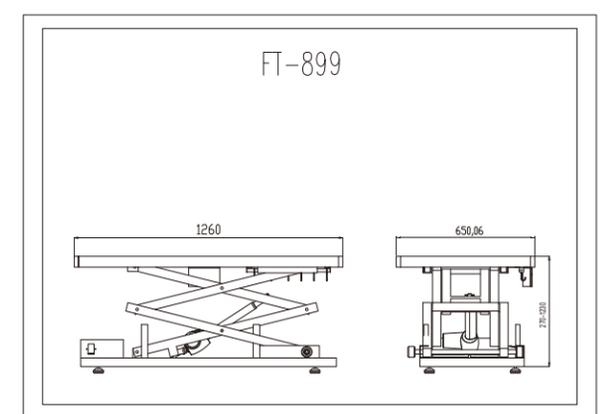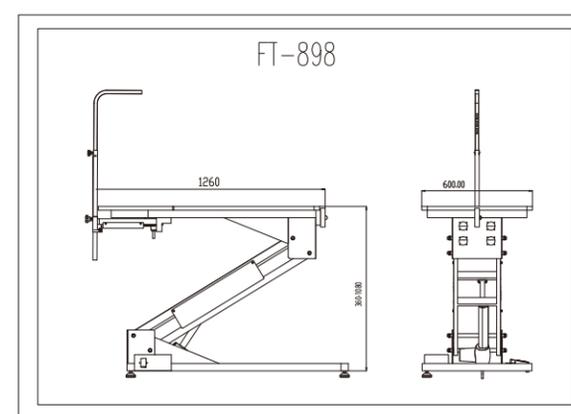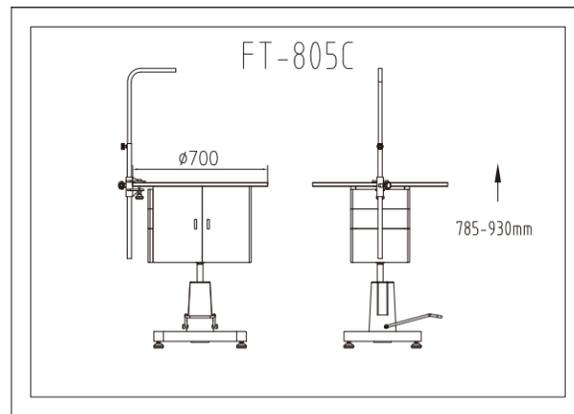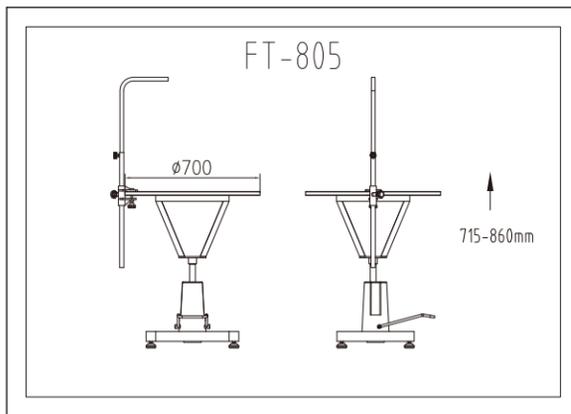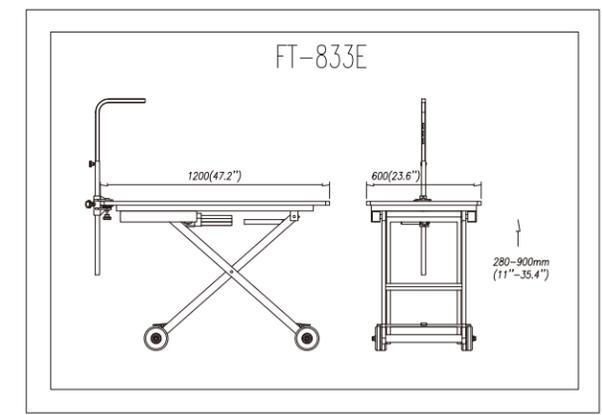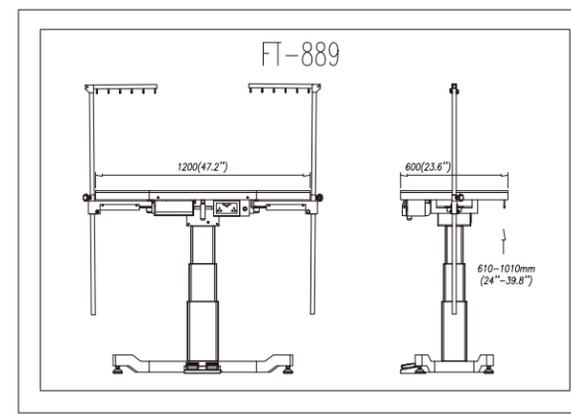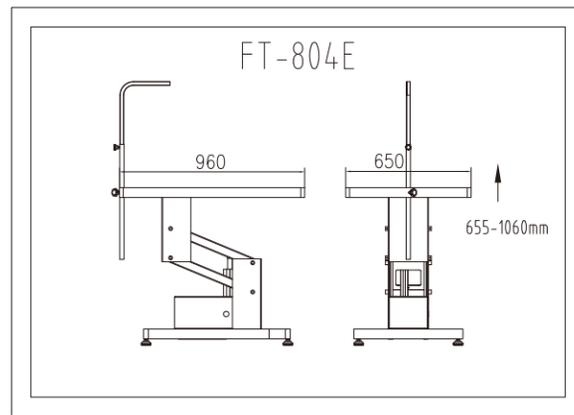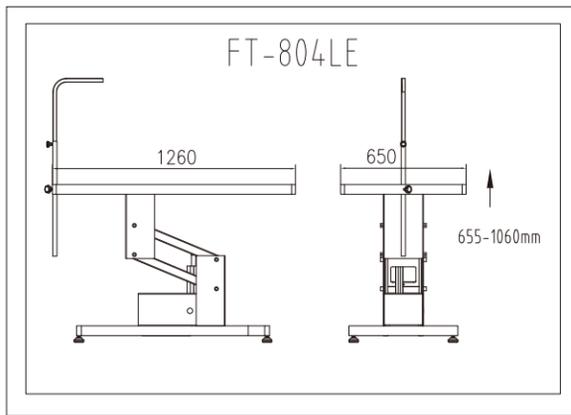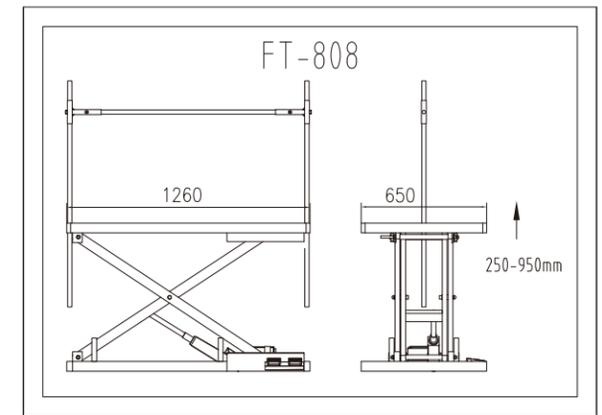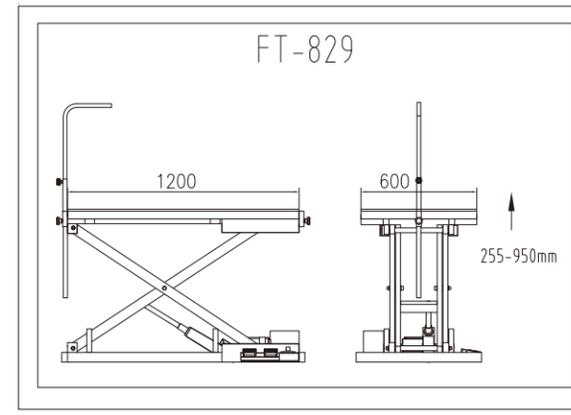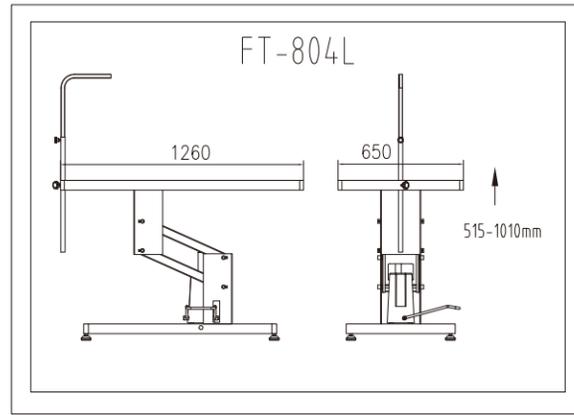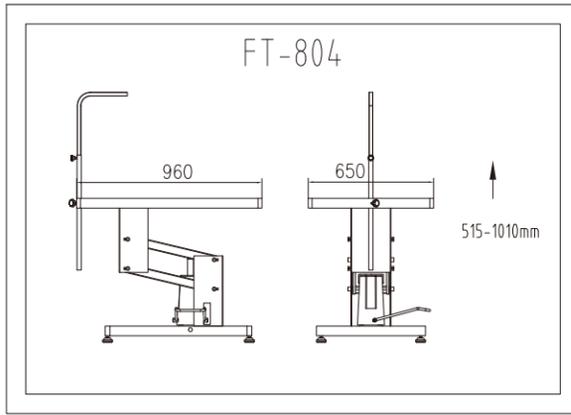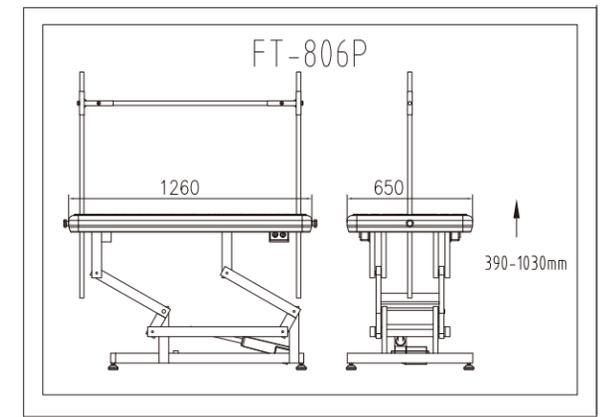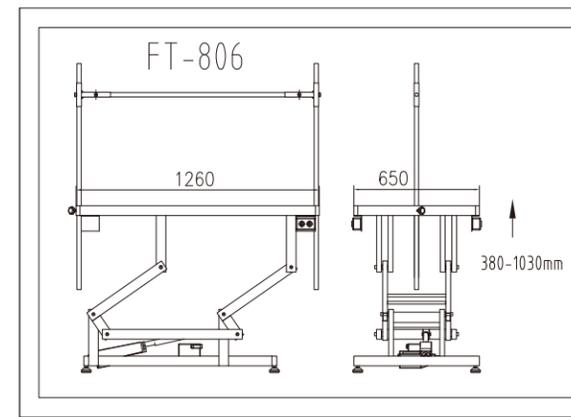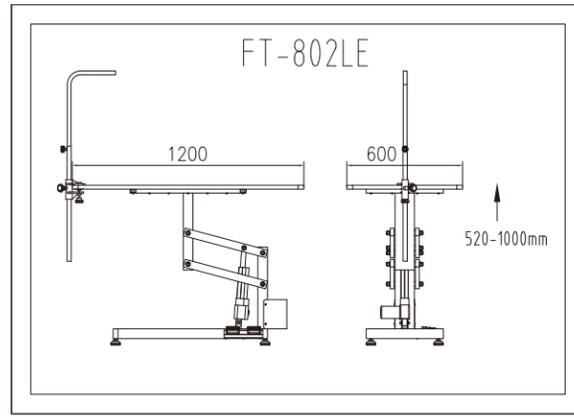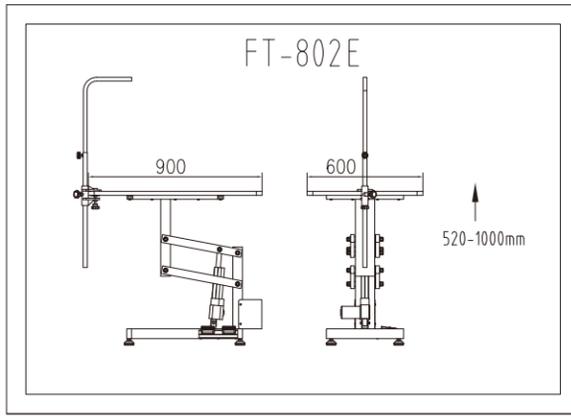

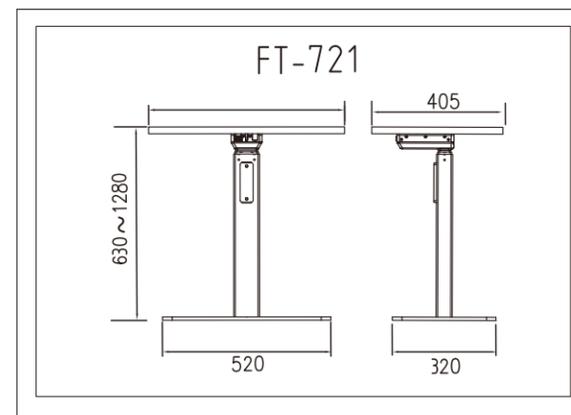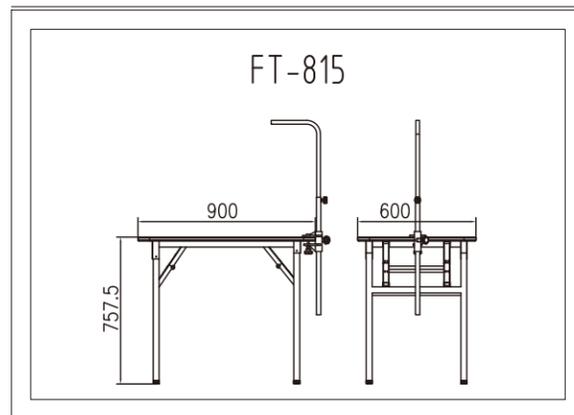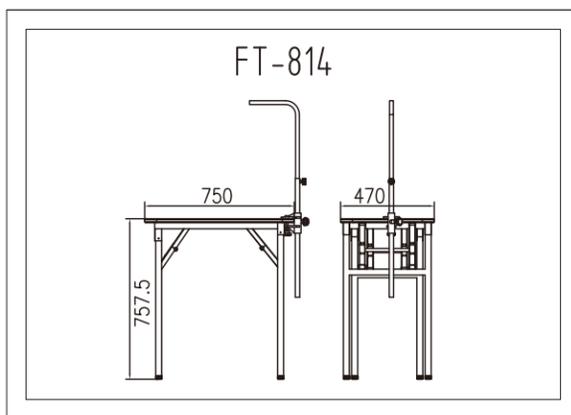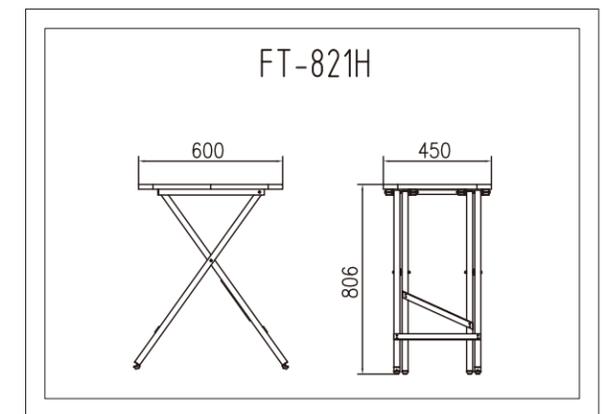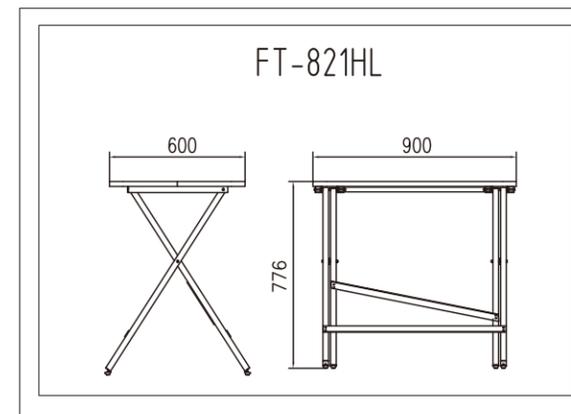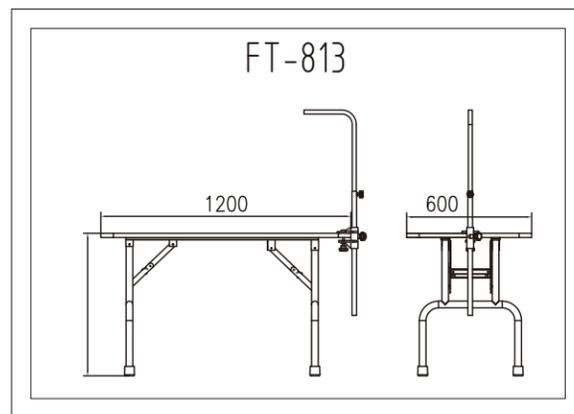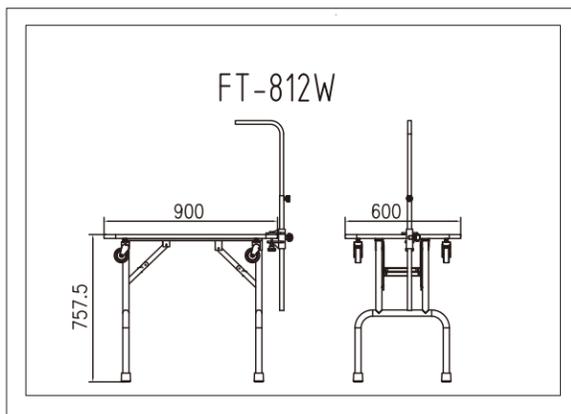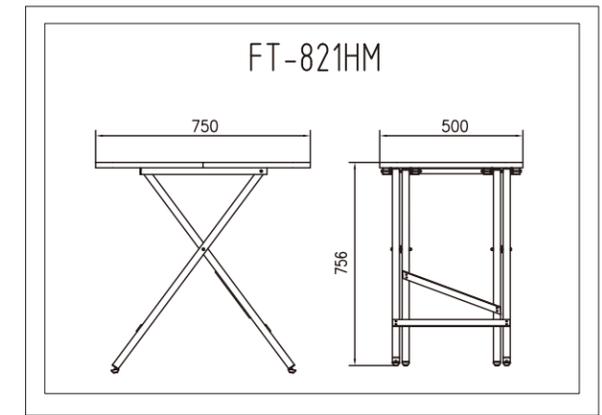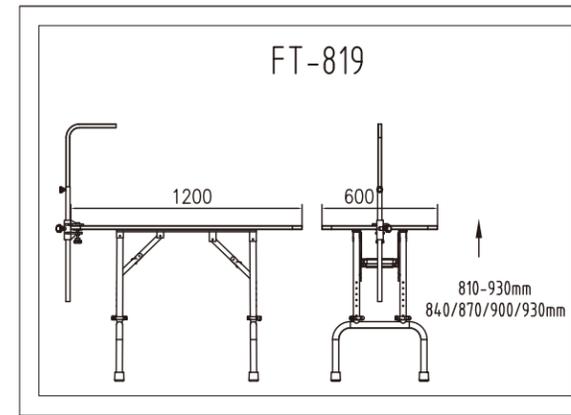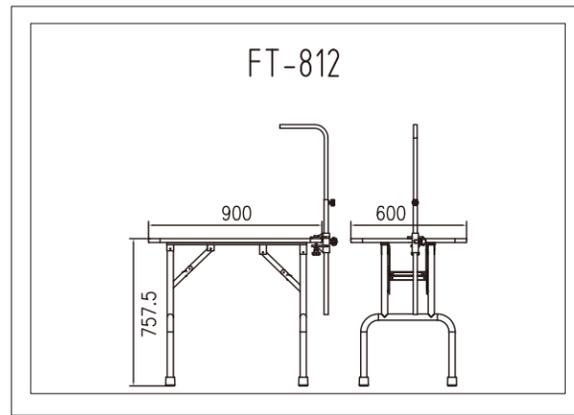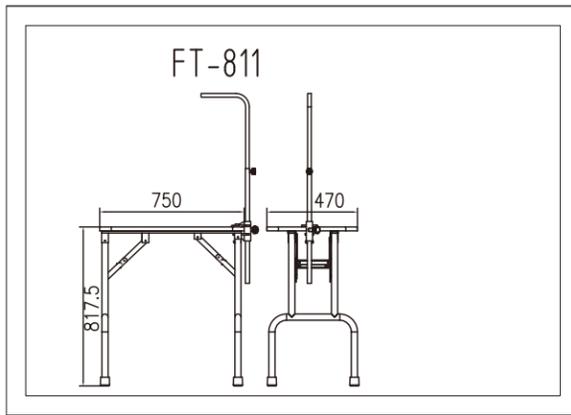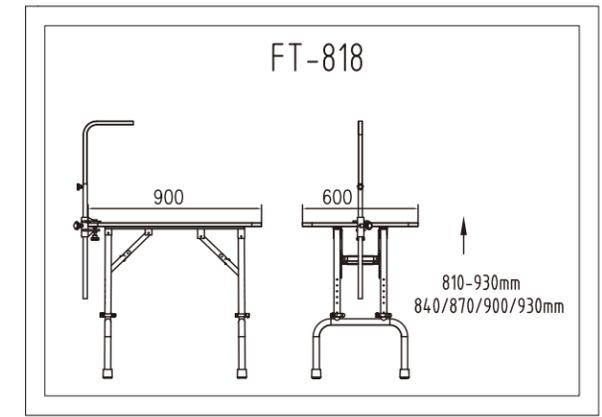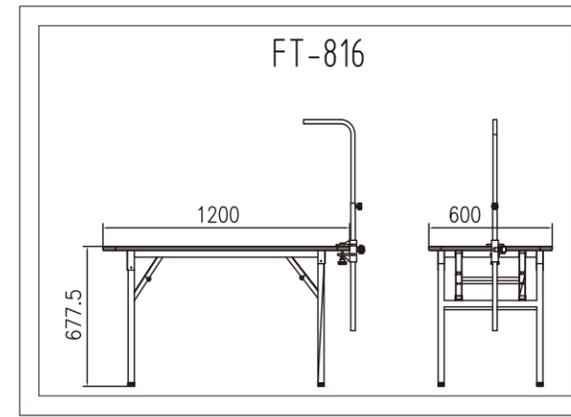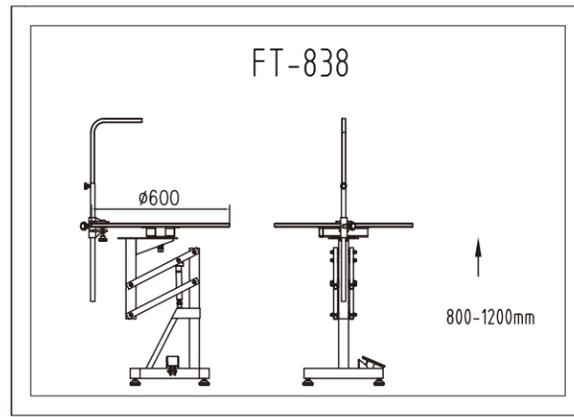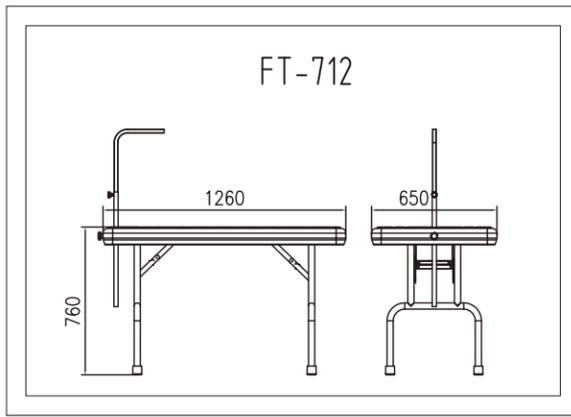

Pro for Pros

KISS Grooming TD-907 Compartment Dryer

There have been a few incidents in the recent past where neglect has led to pets being endangered and even to being fatally injured.

We are all pet people and we all have the same concerns about these issues.

We have taken a scientific approach to refining our dryer and making it safe for use. It is a very valuable tool, but like all tools the user must be educated and experienced.

Is there a risk of asphyxia?

KISS Grooming TD-907 Compartment Dryer has eight louvered air holes on the backboard of the standard model.

Why did we choose louvered (not circular, square or irregular shape)?

Why there are eight holes (not seven or nine holes)?

After intensive testing, we found that the louvered air hole offers the best ventilation and air flow. It helps to maintain the desired interior temperature and ventilation. Again after extensive testing, we determined that eight air holes balance the supply of air and the rate it moves (with the aid of two side draught fans) to realize enough fresh air under the prescribed heat setting.

When the temperature is considered unsuitable for safe operation of the dryer?

KISS Grooming TD-907 Compartment Dryer has adopted an internationally acclaimed, advanced performance digital temperature.

It accurately controls and governs interior temperature of the dryer with a minimal varying within 1 degree of the set temperature. PTC ceramic material heater blocks help resist unwarranted increases of temperature due to surface heat absorption. The heater will automatically power off to guarantee pets safety when the surface temperature is higher than the set threshold value. This eliminates the risk and prevents the machine and pets from overheating.

Will pets be anxious in the dryer?

KISS Grooming technology Research & Development team has investigated the reason of pets' performance anxiety when placed in the compartment dryer.

Is it because of fear of enclosed spaces; is it the fear or discomfort caused by the difference of temperature and humidity?

Our long term research yielded an interesting result---the protection grid/fence. The earlier edition compartment dryer had this stainless steel protection grid/fence behind the glass as a safety measure. This protection fence was the cause of anxiety in some animals.

The second edition of TD-907 adopts thicker fiberglass and eliminates the protection fence. The results are a much safer and less stressful and inviting environment.

Does noise makes pets feel uncomfortable?

The answer in many cases is Yes! KISS Grooming Global's Research & Development team worked intensely to corral, control and abate the amount of noise generated by the compartment dryer. The unique design of our air flow system, strong stable fan and air outlet devices combine for one of the quietest compartment dryers made.

When you open the front door of TD-907, you will hear soft low-frequency sounds of wind only without any disagreeable discord. When the door is closed, you will be pleased with the lack of noise.

Is it safe for pets to be subjected to a full body blow dry?

Yes. It is safe. Exactly calculated inner measurements of heat and air flow are controlled and directed by a swing fan. It provides a highly efficient circulation air system. The bottom space serves a dual purpose. Allowing for easy waste management as well as directing air flow to the lower portion and underside the pets carriage. Even when pets are lying prone, their abdomen hair is being




dried. Extra moisture will be carried off by the ventilation system, always maintaining a low humidity, aiding and enhancing the effects of drying.

How long will it take to dry a large breed?

An often asked question without a fixed answer.

All dogs have different thickness and types of coats (single coat or double coat, etc.). There are variances even within specific breeds.

Generally, we suggest using a forced air dryer, operating on a somewhat high wind speed and heat setting for about 3-5 minutes. After the majority of the excessive moisture is removed, then put pets into the compartment dryer and set 1-50 minutes drying time according to hair type and size/breed of dog. Check the

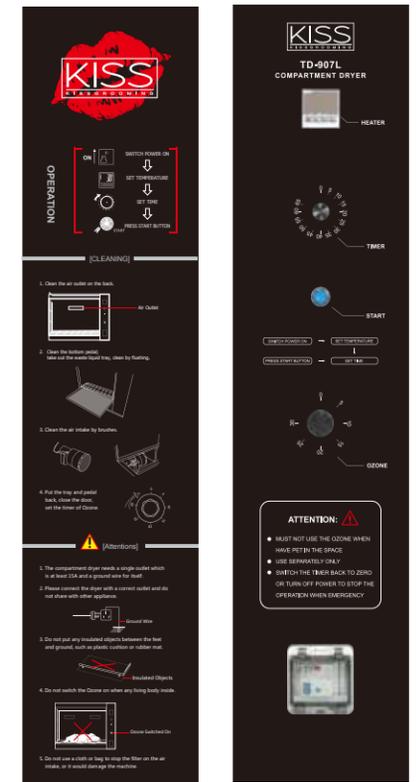

concealed parts such as ear, inner thigh. For most golden retrievers, we use a forced air dryer for about 3 minutes; then use the compartment dryer 20-25 minutes; and finish (if necessary) with a handheld or forced air dryer on the difficult to dry parts; about 2-3 minutes.

Is the labor savings justified?

Experienced groomers and shop owners know that it takes relatively short time to bath a dog. It can take an hour or more to dry some dogs (especially double coat and thick coat dogs such as the Samoyed and Alaskan Husky).

This compartment dryer replaces a lot of labor. In many situations the saving of labor and time can result in dramatically increasing your grooming output and the amount of dogs groomed per day.

Compartment dryers are highly appreciated, approved, and coveted in the international market.

They come close to being a substitute staff member in most shops, proving added labor and increasing revenue as a result.

Is daily maintenance and repair complicated?

After over 10 years research and development by our team, we've upgraded several component parts and greatly improved performance. It is easy and simple to maintain and repair this machine.

Minimal daily maintenance is required. Clean the waste liquid tray; remove the hair from filter on the fan mouth; clean the dryer interior regularly after operation.

The compartment dryer has a switch to prevent voltage/current excess surges. It will power off if short circuit or electric leakage occurs. If a problem occurs, open the control panel; everything is clearly marked on each electric wire. Follow the manual trouble shooting instructions to fix most common problems easily and quickly.

How do I clean inside of the compartment dryer?

The entire compartment dryer is made of high quality SUS304 stainless steel, discouraging rust buildup and making cleaning easy. Depending upon frequency of use, we suggest a thorough cleaning once a month. The two interior side grids are removable. You should also remove the screws periodically; also deep clean the air outlet and fan to extend the life of the machine. Use neutral detergent to clean the dryer interior. Use ozone generator (which located in the dryer) to eliminate the odors.

Are there any other advantages (besides saving time and labor)?


The compartment dryer's value is greatly enhanced when working on difficult animals (especially cats). Most groomers have been bitten by cats during the drying process. They are more sensitive and less tolerant during the drying process than dogs. Drying them in the compartment dryer is the answer. Most cat hair is easily dried in the compartment dryer. Aggressive, ill and older animals that may not like the repeated contact needed during the drying process also will benefit from this wonderful compartment dryer.





Compartment Dryer | TD-907/907L

Automatic shut-off timer, could be set from 5 to 30 minutes.
 Suitable for all breeds of dogs and cats, especially older and nervous animals.
 When on, this dryer is amazingly quiet, virtually vibration-free. It appeases the most sensitive animals.
 High air volume helps to remove and dry water from the coat.
 Specialized PTC heating element keeps the temperature safe and steady from 35°C (95F) to 50°C (122F).
 Two powerful motors generate a huge volume (2,700 cubic feet per minute) of air. Enough to dry a thick coated golden retriever in about 25 minutes.
 Efficient and reliable air exchange system makes sure the humidity is low and enough fresh air circulates for dogs' comfort.
 Stainless Steel faceplate and grill are durable and easy to clean, helping to maintain a sterile environment.
 Ozone function removes odor from the environment and refreshes air inside the compartment.
 (Warning: Ozone should never be used when there is a dog inside the dryer).
 Ease of operation makes this an incredibly valuable and time saving machine.

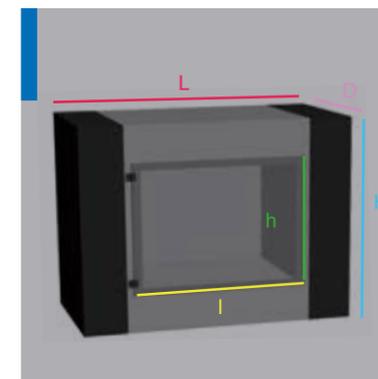

TD-907: inner size L81×D65×H69cm
 TD-907L: inner size L106×D65×H105cm




Over flow & leakage breaker




OMRON intelligent temperature controller



Removable middle divider

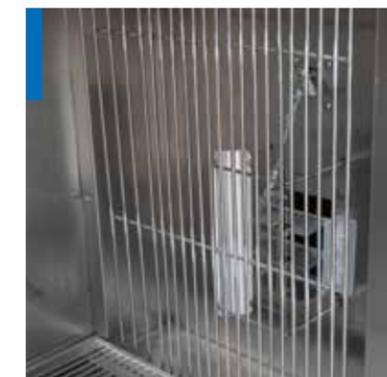

Stainless steel interior

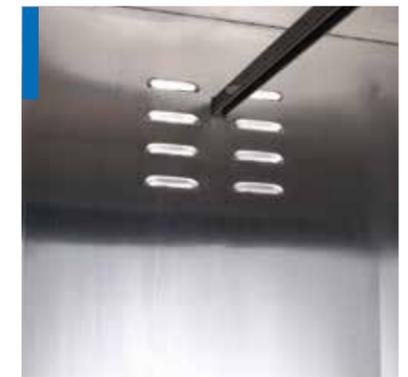

Fresh air intake

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand	Package	
TD-907	110	165	141×66×115	49-39/64×25-63/64×42-33/64	132×74.5×120H	1.19
TD-907L	120	175	176×66×148	69-11/16×25-63/64×55-29/32	183×73×157H	2.09



Black for optional color.

Pet Drying Cabinet

| TD-907T

This drying cabinet is a revolutionary product!

Groomers just need put pets in and press some buttons to start drying. When this pet is being dried, groomers can turn for another drying.

Cabinet drying brings groomers a totally new drying style. No tables, no blowers, no noise, no hot air, no close touch and scratch, no pets out of control, no floc and hair everywhere.

What you need to do is just to open, put, close, press, open and take out.

To clear and clean the filter and internal surface now and then is also a piece of cake.

Less labor and expense, more customers to serve.

Especially good for sensitive species like cat, avoiding scratching accidents significantly.

To change or not to change, that is not a question.



Quiet Running

TD-907T is up to 90% quieter than conventional blowers. Even cats, which are more sensitive than dogs, will feel more comfortable, be relaxed and quiet with high probability.



Super Efficiency

TD-907T saves labor, expense and time. The enclosed circulation of heated air makes a fast and energy-saving drying. When one pet is being dried in the cabinet, groomers can start drying another. Much greater efficiency will thus be achieved.



Safe Grooming

Many cat/dog scratch accidents happen during drying process. These potential dangers are eliminated by TD-907T. The vent holes at the bottom of the door ensure pets' refreshed breath. Anti-overheating system will shut down the machine when the temperature is detected over 40 C.




Near Infrared Treatment

The Infrared Therapy stimulates animal body recirculation, helps raise tissue cell vitality, reduces granulation edema, improve partial tissue nutrition and promotes recovery greatly.



Negative ion Therapy

TD-907T will dispense hot air with highly-concentrated negative ions to make the coat smoother, more moisturized and easier to comb and brush. Negative ions can also purify the air in the cabinet.



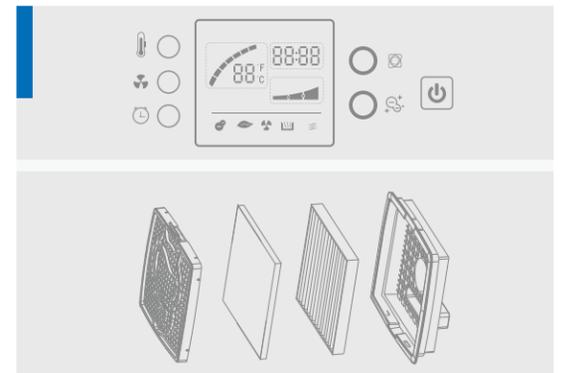


Simple Operation and Maintenance

What you need to do is just to choose a suitable working mode, set temperature, wind speed and the running time. The inside of the cabin and air inlet filters are easy to clean and be replaced.



TD-907TW-XL

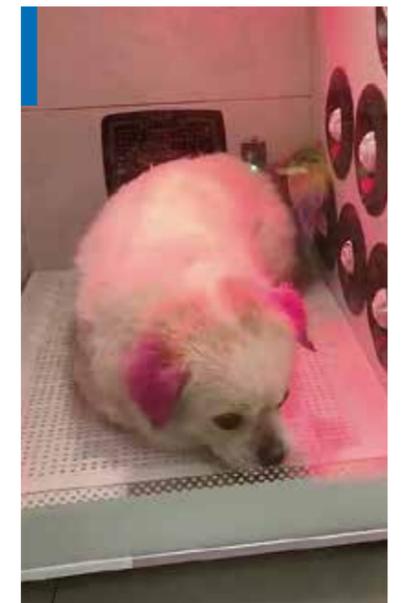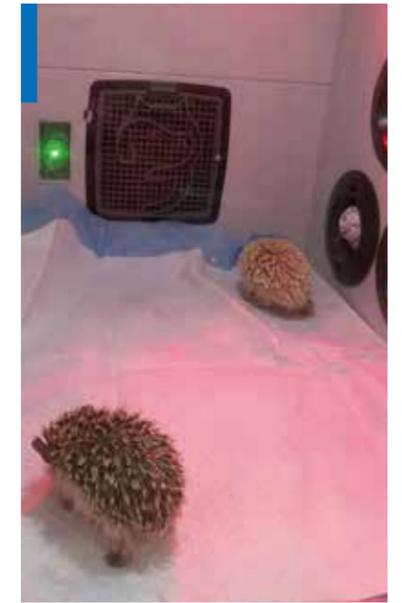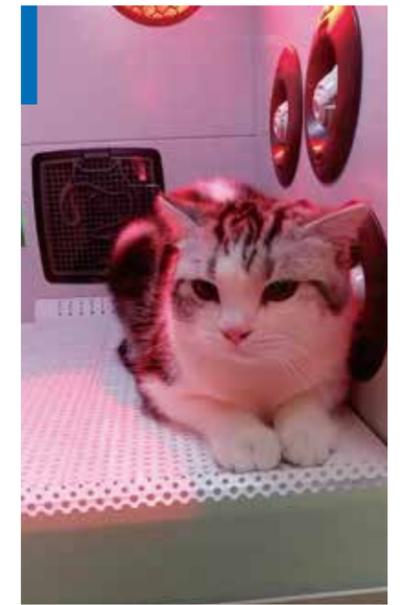

Pet Drying Cabinet

| TD-907TW

- Big room enough for drying a large-size dog, or drying two middle-size pets when the dividing panel is inserted;
- Quiet and cat-friendly;
- Two modes, three-stage wind speed, multi-stage temperature;
- Three-layer filtering system collects hair and floc, making cleaning an easy job;
- Super cost-efficient and labor-saving;
- Safe and comfortable both for groomers and pets;
- Negative ion, aroma and oxygen therapy make a pet-friendly room;
- The infrared therapy serves pets with added value;
- Easy operation, simple maintenance;
- Fashionable and beautiful appearance;
- Technology transforms pet-drying and makes a simpler grooming.






Ionic Wall-hanging Dryer | TD-905

Variable wind speed and temperature control.
 The maximum wind speed is 38m/s.
 Japanese made heavy duty 1HP fan, air-cooled brushless motor, reliable heating element and air adjustment.
 Temperature range is from 25°C (77F) to 75°C (167F).
 Maintenance-free brushless motor, does not require oiling or greasing of internal parts.
 Energy efficient, economic to run.
 Extremely quiet (lower than 40 db).
 Fashionably designed by Italian Master Designer.
 ABS body compartment, resistant to harsh chemicals and scratches.
 Double layer housing on the heater prevents over-heating.
 The outlet tube is made from heat resistant PPS material (resisting over 260°C (500F) temperature).
 Durable hanging arm could be lifted and lowered thousands of times.
 Advanced tourmaline brushes generate negative ions that aid in containing internal cuticle moisture and taming unruly hair and straightening curly coats.
 Negative ions also help in purifying the air and creating a healthier environment.



The 3rd Generation of Classical Ionic Dryer

Highly Recommended by Jennifer Lee and Olga Zabelinskaya, ICMG, Member of Hall of Fame in Groom Team USA.

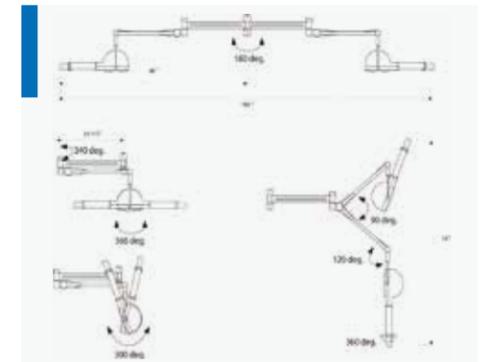

Specialized Pendant Available For Ceiling Mounting

2010

World is not Black and White.
Vivid and decorative colors were painted.



2014

A worldwide reorganization.
This dryer series has covered over 100 countries in this world.



2014

A "geek" dryer.
Ultra strong ionic generator, anti-bacterial plastic housing, brushless motor, etc, who said a "traditional product" like dryer cannot be "high-tech"?



2017

TD-906
And, Now in **RED**



2008

The beginning of the whole story.
Bringing maybe the first ionic styling dryer into this market.



2011

A Must in grooming salon.
The concepts in groomers' mind were turned from "could have", to "better have", now to "must have".



2016

New generation with new art touches.
Lighted touching buttons, beautiful new mechanism, shining flamingo pink and pearl white new colors.
Refreshing art elements added.



2019

We are Back,
We are Black.





Ionic Stand Dryer

| TD-906

Experience hands-free fluff drying.
 This gentle, powerful finishing dryer produces wind speed of 38m/s.
 Temperature range is from 25°C (77F) to 75°C (167F).
 Japanese made heavy duty 1 HP fan, air-cooled brushless motor, reliable heating element and air adjustment.
 Variable air speed and temperature controls, suitable for all kinds of breeds.
 Sturdy powder coated stand with durable locking casters.
 Balanced weight design provides stability.
 Available in two lovely colors, **PINK FLAMINGO** and **PEARLIZED WHITE SWAN**.
 As a result of our advanced technology, the noise level is greatly reduced.
 The washable foam filter is easy to remove and clean.
 Nozzle rotates 360°, easily positioned at any angle.
 Advanced tourmaline brushes generate negative ions that aid in containing internal cuticle moisture and taming unruly hair and straightening curly coats.
 Negative ions also help in purifying the air and creating a healthier environment.



Separate Controlled Negative Ion Generator and Heating Element.



Tension Adjustable Ear Mount



Anti-rusting Plastic Wheels



Easy Adjustment for Ergonomic Height Setting



The 3rd Generation of Classical Ionic Dryer

Noise Terminator Central Dryer System

| TD-908

The medical world has long known that excessive noise exposure triggers hearing loss and contributes to other illnesses (i.e. heart disease, anxiety, and stress related illnesses).

As a means of helping to curtail the catalyst of these physical ailments and to protect and improve groomer's health and working environment, we have spent years researching the noise situation faced by groomers.

A breakthrough has been made!

The result is our Noise Terminator which incorporates a technology that employs a means and solution to noise reduction on pet dryers.

The Noise Terminator is a combination of many advanced technologies, including Noise Conduction Control, Noise Absorption, and Noise Transmitting in Air.

Now that this revolutionary technology is used on our Central Dryer System, the noise has been reduced to an amazingly low level of 65dB. This is a major achievement in the field of pet dryers.

The system very effectively does the work of two powerful dual motor dryers. This allows two groomers to work separately, drying two dogs simultaneously.

We could even incorporate 4 dryers into the housing, making the Noise Terminator more powerful. Four dryers inside with four individual wind outlets allows for 4 groomers to work at the same time.

We are not just reducing the ambient noise but also allowing you the option of removing the noise source from the grooming area.

Two 10m (over 32') long flexible hoses transfer the 90-95m/s high speed wind from the core unit to grooming table.


The hoses can be mounted to the ceiling of the room. This allows for ease of access, keeping them close at hand yet out of harm's way.

When the hoses are attached to our smartly designed spring balanced armature, they are available, yet out of the way. When you need use of the dryer, just pull the nozzle down gently, and then put it back when you finished.

To make this machine incredibly user friendly and convenient, you are able to control the dryer by pressing the button on the hand-held nozzle. This eliminates the need for you to walk to the unit to switch it on/off. There is a main One-Button-Start system.

You can place the unit in grooming room, which looks like a cabinet, or put it out of the room to cut the noise dramatically.

Ergonomic; people, pet and groomer friendly; this easy to use design brings you an extraordinary experience, while creating a healthier and more pleasant working environment.

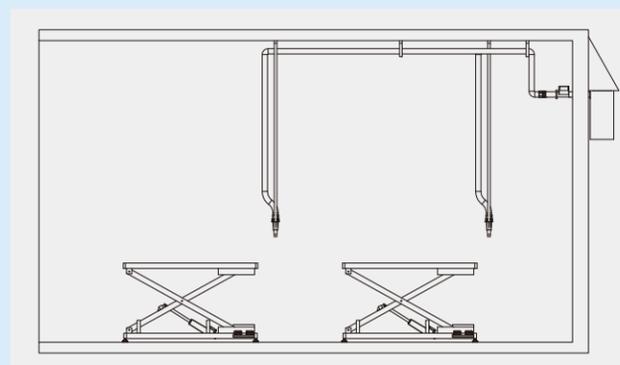

Prevent Noise Pollution

Restore Ear Health






GEMINI Parallel Dual Motor Dryer | TD-900CG

Don't call me "INSANE"!

After the first batch of Gemini was sold out, many reports from all over the world commented it has "insanely strong blow force"!
 Well, although it may be the strongest dual motor dryer in the world at this moment, please do not call it "insane". It dries large breeds fast!
 Incredible blow force of 1250g+ exceeds most similar dryers in the world.
 Variable speed control.
 Natural heat generated by two powerful motors.
 Conveniently mounts to the wall with specially developed wall mount bracket.
 10ft heavy duty flexible hose and 3 different types of nozzles included.

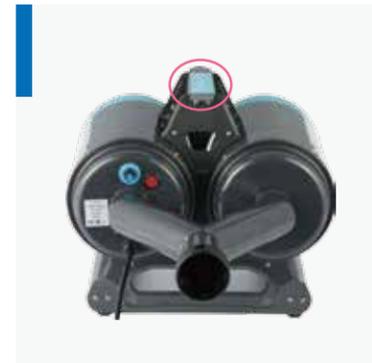

Specially Designed Handle for Safe Carry and Cable Organizing



Triple Rear Filtering System



Metallic Dark Grey Coating



Variable Speed Control with Natural Heating Output



Wide Nozzle Selection



Optional Stainless Steel Wall Mount Available



AEO LIAN

| TD-901GT

This is the most aptly named dryer in the pet world, AEO LIAN hosed dryer! Incredible power, incredible value. This single motor dryer outperforms most dual motor dryers in this market. Its wind speed and blow force has been experienced to be believed. This dryer is revolutionizing the world of single motor pet dryers. It's versatile too. Mount it on the wall or stand mount.



Hard extension pipe, for turning the hosed dryer into a finish dryer.

MANGO

| TD-901MT

Mango's single motor superior blow force and wind speed allows it to perform better than most double motor dryers. When equipped with a hard pipe extension, it turns from a forced air dryer into a finishing dryer instantly. Patented noise control technology brings the noise level down to about 80 decibels.





HERCULES TD-900XT

Hercules, The Hero in the name of Hero!
Super Blow Force

Incredible blow force of 955 grams and even more amazing wind speed of 95m/s, beats every recorded cyclone on earth, make TD-900XT the most powerful dual motor pet dryer available.

Motor Life Span

Longer carbon brush of 38mm (1.5"), is the longest currently on the market. Smaller sparks created by smaller carbon pellets and smoother rotor give the TD-900XT a longer lifespan (by 30%) compared to the previous model.

Energy Savings

A better air flow system and a smaller flow friction force, has allowed TD-900XT to perform at a higher energy efficient ratio, yielding a 30% savings in energy consumption over previous models.

Use our wall mount to get the dryer off the floor or countertop and away from the work space area. 10' hose provides enough length to allow you to walk around a grooming table. Variable heat and power controls.



CYCLONE TD-901T

There isn't a better value in the grooming realm than the Cyclone.

It is one of the best single motor pet dryers available.

Super Blow Force

New Launched 52mm (2.04") diameter hose system and highly improved motor, give the user a much higher blow force of 605 grams, an increase of 45% over previous model; unmatched wind speed of 60m/s, an increase over previous model by 40%.

Motor Life Span has been extended

Longer carbon brush of 38mm (1.5"), is the longest currently on the market. Smaller sparks created by smaller carbon pellets and smoother rotors give the TD-901T a longer lifespan (by 30%) compared to the previous model

Energy Savings

A better air flow system and a smaller flow friction force, has allowed TD-901T to perform at a highly energy efficient ratio. The TD-901T, 1200 watt motor, outperforms the 1600 watt dryer currently on the market.



SUPER CYCLONE TD-901ST

A very versatile dryer.

Performs the work of three dryers; forced air blower, finishing dryer and also an air massager; the most versatile dryer in our line.

Super blow force at 605 grams and fantastic wind speed at 60m/s.

Top grade SUS304 Stainless Steel promises years of use without rusting.

When this variable-speed dryer is detached from the workstation stand it can be used as a portable or floor dryer. The air nozzle swivels 360 degrees for accurate drying.

Weighing approximately 6kg (13-1/4lbs) makes it a portable dryer, easy and conveniently carried anywhere you need. It has dual mounted legs which allow both vertical and horizontal use.

Full steel workstation stand available.




Handheld Dryer TD-904N

Easy to maneuver in hand, lightweight and comfortable to handle.

High volume negative ion generator built-in, for closing cuticles and soothing hair.

Running quietly, good for long time closely use.

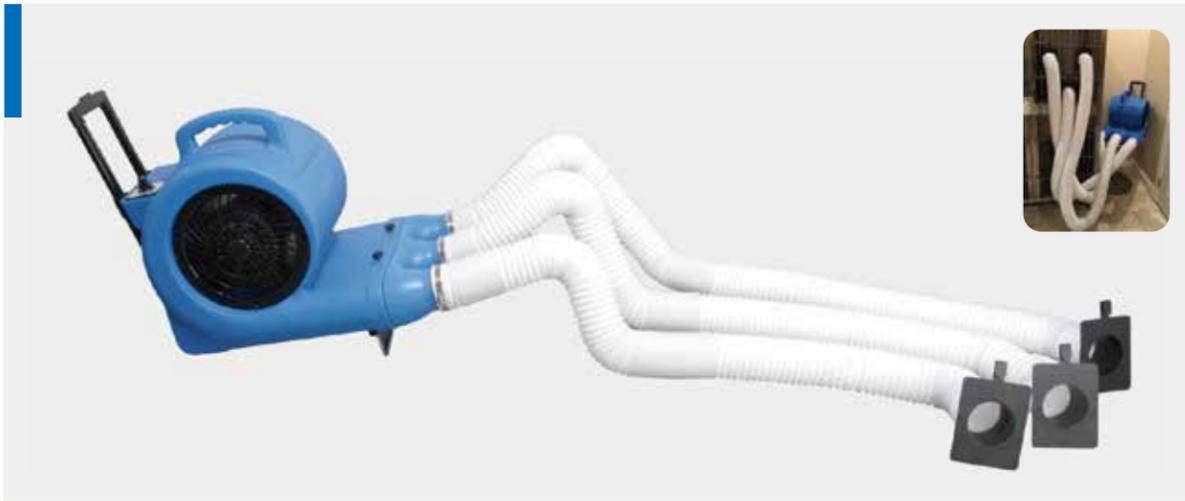

Cage Dryer

TD-903

900 Watts - powerful enough to make sure all nozzles have equal and strong wind pressure.
 3-speed control - wind speed from 15m/s to 24m/s.
 Delivers large volume and flow of fresh air into cage compartments at natural environmental temperature.
 Virtually maintenance free – no brushes to replace and no heating element to burn out.
 Both intake ports have filters that eliminate hair and debris from entering the motor.
 Includes a snout adapter, three long-lasting, heavy-duty 6' hoses, three cage attachments and clamps.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
TD-903	20	21.8	56×58×55	22-3/64×22-53/64×21-21/32	56×58×55H	0.18



TD-9013
Dry Stand for Single Motor Dryer



TD-9012X
Wall Bracket for Hercules



TD-9012
Wall Bracket for Single Motor Dryer



TD-9011
Groomer's Third Arm

Dryer Buying Guide

ITEM NO.	MAXIMUM POWER	HEAT POWER	BLOW FORCE	MAXIMUM SPEED	HIGHEST TEMPERATURE	NOISE	NW. (kg)	GW. (kg)	DIMENSION	PACKING MEASUREMENT	NOZZLE	ACCESSORY	ADVANTAGE	SUITABLE FOR
TD-901GT	2800W	600W	895g	90m/s	52 C	86dba	4.75	5.5	340×160×210mm 13.4"×6.3"×8.3"	450×300×240mm	1 Round 1 Flat Narrow	Dryer Stand Wall Mount Third Arm	Strongest Wind Speed And Blow Force Among Single Motor Dryers	Professional Groomers and Breeders
TD-901MT	2800W	--	895g	90m/s	52 C	86dba	5.25	6.0	340×160×210mm 13.4"×6.3"×8.3"	450×300×240mm	1 Round 1 Flat Narrow 1 Hard Pipe	Dryer Stand Wall Mount Third Arm	2-in-1 Combo Dryer	Professional Groomers and Breeders
TD-900CG	2400W	--	1250g	125m/s	65 C	85dba	7.00	10.0	500×350×300mm 19.7"×13.8"×11.8"	515×375×350mm	1 Round 1 Flat Narrow 1 Flat Wide	Wall Mount Third Arm	The Strongest Blow Force Performance	Professional Groomers and Breeders
TD-901T	1800W	600W	605g	60m/s	52 C	83dba	4.3	5.1	340×160×210mm 13.4"×6.3"×8.3"	450×300×240mm	1 Round 1 Flat Narrow	Dryer Stand Wall Mount Third Arm	Full Stainless Steel Housing	Professional Groomers and Breeders
TD-901ST	1800W	600W	605g	60m/s	52 C	83dba	4.3	5.1	340×160×210mm 13.4"×6.3"×8.3"	450×300×240mm	1 Round 1 Flat with Stainless Steel Long Extension	Dryer Stand Wall Mount Third Arm	Full Stainless Steel Housing	Professional Groomers and Breeders
TD-900XT	2400W	---	955g	95m/s	65 C	85dba	6.4	7.8	485×160×210mm 19.1"×6.3"×8.3"	600×300×240mm	1 Round 1 Flat Narrow and 1 Flat Wide	Wall Mount Third Arm	Superior Wind Speed and Blow Force	Professional Groomers and Breeders

ACCESSORY	UNIT	REMARK
Hose (D40)	Piece	For TD-901, TD-900X
2nd Generation Hose (D52)	Piece	For TD-901T, TD-900XT, Monster, TD-901ST
Carbon Brush	Pair	1 pair for single motor dryer and 2 pairs for dual motor dryer

ACCESSORY	UNIT	REMARK
2nd Generation Carbon Brush	Pair	1 pair for single motor dryer and 2 pairs for dual motor dryer
Dryer Motor	Piece	All accessories are suitable for AEOLUS brand dryers



Dog Show Aluminum Portable Trolley

KB-511 Series

Must-have star product for dog show of Europe and America.
 Multi-role —Cage, grooming table, work table, luggage cart, pet carrier, transportor, etc—for your multi-purpose at dog show.
 Heavy-duty aluminum parts, rust resistant, durable, strong, resilient in the roughest conditions.
 Light-weight, easy to clean. Simple to be loaded with your dogs and luggage or supplies.
 New and upgraded: rubber denoising protector—eliminate noise caused by metal bump, which often happened to similar products.
 Compact size fits in most standard cars and is easy for transportation and storage



Double door locking system



Fold down handle



Unique rubber denoising protector



Unique rubber denoising protector



Fixing slot of cage divider



Wavy, grooved skid-resistant surface



1.Folding state



2.Single-cage Trolley



3.Double-cage Trolley



4.Four-cage Trolley



5.Four-cage Trolley

Trolley Accessories




Trolley Skirt



Tool Organizer Cover



Anti-Slip Grooming Mat

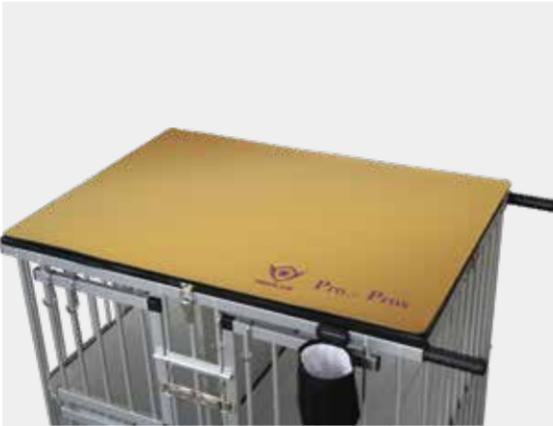

Adjustable Angle Umbrella with Clamp



Feeding Bowl with Mount



Hanging Water Bottle with Clamp



Cup Holder



Hanging Tool Box



Foldable Grooming Arm with Clamp



Attachable Carabineer Hooks




Tool Bag



Protective Covering for Trolley



Battery Powered Cooling Fan with Mini-USB port

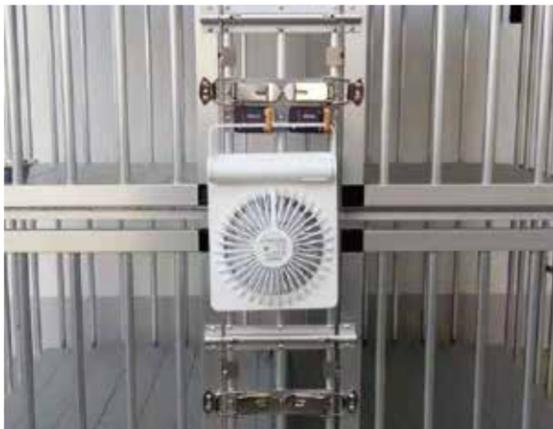

6" Heavy Duty Castor Set



Lightweight Collapsible Travel Chair



All-Terrain Foldable Travel Chair



Pink, Purple, Blue, Black and Silver Colors Available




Customized specialize colors available!

ITEM	Gross Weight (kg)	Dimension(cm)		Volume (cbm)
		Stand W×D×H(cm)	Package	
KA-511NAO	17.3	93×65×80H	95×84×16H	0.17
KA-511NAT	19.7	93×65×80H	95×84×16H	0.17
KA-511NAF	23.3	93×65×80H	95×84×16H	0.17

ITEM	Gross Weight (kg)	Dimension(cm)		Volume (cbm)
		Stand W×D×H(cm)	Package	
KA-513NAO	13.5	78×54×63H	80×65×16H	0.12
KA-513NAT	15.2	78×54×63H	80×65×16H	0.12
KA-514NAT	8.5	55×44×80H	58×84×16H	0.08

Aluminum Alloy Enclosure


Made of 100% adonic aluminum oxide alloy, strong and durable. Folding design and easy to carry and set up. Used as preparation area of dog show or exhibition area of dog sale.

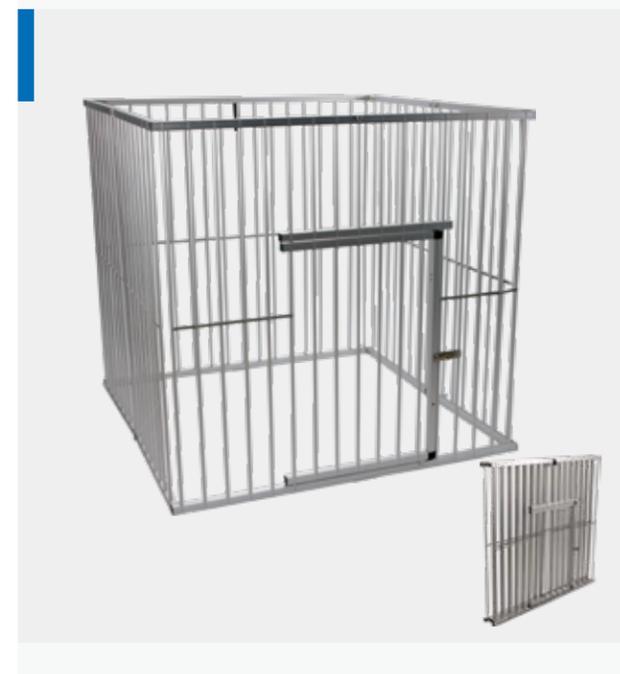

Item	Net weight (kg)	Gross weight (kg)	Unfolded dimension		Package (cm)	Volume (cbm)	
			One-piece width	Height			
Hexagon	EP-HEX-2824	6.0	7.0	70	60	71×23×62	0.10
	EP-HEX-2836	8.4	9.4	70	90	71×23×92	0.15
	EP-HEX-3624	7.8	8.8	90	60	91×23×62	0.13
	EP-HEX-3636	10.9	11.9	90	90	91×23×92	0.19
Octagon	EP-OCT-2824	8.0	9.0	70	60	71×28×62	0.12
	EP-OCT-2836	11.2	12.2	70	90	71×28×92	0.18
	EP-OCT-3624	10.4	11.4	90	60	91×28×62	0.16
	EP-OCT-3636	14.6	15.6	90	90	91×28×92	0.23

Item	Net weight (kg)	Gross weight (kg)	Size(cm)		Volume (cbm)	
			Unfolded dimension	Package		
Square	EP-4424	7.2	8.2	120×120×60H	124×15×62	0.12
	EP-4436	10.1	10.1	120×120×90H	124×15×92	0.17
	EP-6624	10.8	11.8	180×180×60H	94×27×62	0.16
	EP-6636	15.1	16.1	180×180×90H	94×27×92	0.23
Square	EP-8824	14.4	15.9	240×240×60H	124×27×62	0.21
	EP-8836	20.2	21.7	240×240×90H	124×27×92	0.31
	EPH-4448	13	15	120×120×120H	124×15×122	0.23
	EPH-6648	18.2	20.2	180×180×120H	94×28×122	0.32
EPH-8848	25.5	27.5	240×240×120H	124×28×122	0.42	

Dog Stacker | DS-600

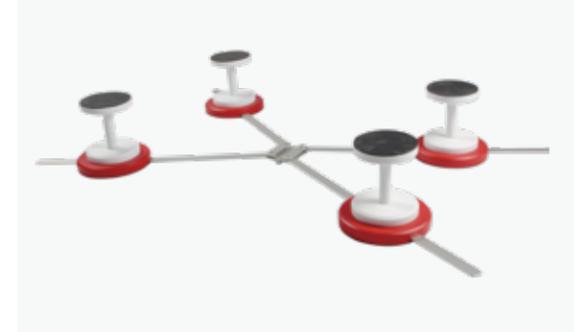

Snake Chain





ITDS-9060

Instrument Trolley

Constructed from 100% 304 stainless steel.
 Safety folds to all edges, with reinforced shelves.
 Smooth rolling 4" nylon wheels with S/S brackets.
 4 styles to choose, with or without rails, with multiple drawers too.
 Neat welding and polishing finish, ensures both decent looking and great durability.



ITFS-5050

ITFS-6050

ITFS-9050



ITFR-5050

ITFR-6050

ITFR-9050



ITDS-9050

ITDS-9060

ITDS-9050-D

ITDS-9090



ITDD-5050

ITDD-6060

ITDD-9090

ITDF-5050



Micro Bubble SPA Tub | BTS-150

This Professional Spa Tub offers unprecedented luxury and therapeutic experiences to the pet grooming industry. Patent pending technology generates billions of miniscule bubbles between 20 -40µm, immediately changing pure water into a creamy milky white liquid that soothes and caresses the skin and coat. The tiny bubbles carry negative ions that attach to dirt particles, lift them to the surface off the dog's skin and out of their coats. When dirt is saturated with micro bubbles, the particles will float up to water surface. The best part is you don't need any shampoo to clean the pet's undercoat, the bubbles do all the work.

When cleaning and conditioning solutions are implemented, the effect is even more pronounced, resulting in better moisturizing, cleaning and smoothing treatment. This is also the first micro bubble SPA tub that has high density Ozone in the bubbles. When air is absorbed into the water, the Ozone level is capable of reaching 14% -15% saturation; this aids in treating numerous skin issues and irritations problems. You could also connect this machine to an Oxygen tank to get 100% pure Ozone infused into the water. The therapeutic effects result is fantastic skin and great looking and feeling coat. There have been many cases where the effects on severe skin diseases, rashes and irritated conditions have been dramatically and positively altered and improved. Two countdown timers, one for SPA treatment and the other for Ozone function. The Ozone can also be switched on and turned off automatically, for maintaining daily sterilization conditions. This SPA Tub is not just a therapy treatment for pet, it's also a revenue generating machine. Treatments should be billed similarly to those charged for grooming; similar to what you would pay at a SPA. The money will flow in once your clients learn of this unique and wonderful SPA.

A petite round model (BTS-151) is suitable for smaller pets; available in red and white base colors.



Milky Delicate Micro Bubbles Spread in the Whole Tub



Micro Bubbles Permeate into Pet Fur and Clear out the Stubborn Dirt



New luxury stainless steel SPA bath tub is coming.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)		Volume (cbm)
			Stand	Package	
BTS-150	93	136	160×65×100	195×89×135H	2.34
BTS-151	49	82	ø66×100	77×68×120H	0.63

What Is "Micro Bubble SPA"? How Does It Work?


The definition of "Micro Bubble" is the bubble which is from 20~40µm. Bubbles that exceed this size range cannot be referred to as true "Micro Bubbles".

We manufacture our "Micro Bubble SPA" tub in its true sense. We added the unique function "Ozone Micro Bubble" that other similar products have not.

Observe that translucent water/liquid will immediately return to the milky white color 30 seconds after the SPA is switched on. This is why we also refer to it as the "Milk SPA".

There is no milk, milk powder any other additives. The white color is solely created by and due to billions of micro bubbles.

Look closely at a waterfall. You see white water instead of clear

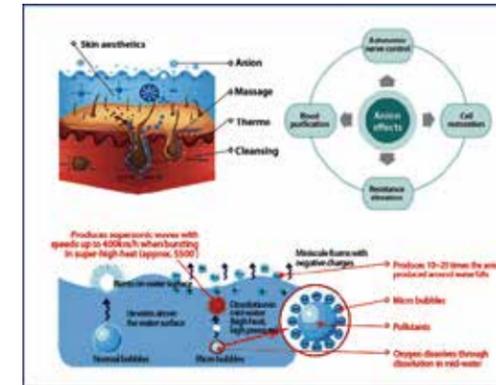

running water. Our SPA tub operates via the same theory, agminated micro bubbles make liquid seem almost opaque when saturated with the micro bubbles. When immersed in the SPA animals relax; their breathing is enhanced due to the soothing effects and increased negative ion release (because of the high concentration of anions in air and water). Their heart rates slows and they thoroughly enjoy the experience. Our SPA tubs dramatically increases the anion release through the high frequency crash of micro bubbles.

It is well known that anions (negative ions) are widely touted in the hair care field. The theory is anions can

keep hair soft and moisturized. The micro bubbles which have negative charges unite water molecules, make hair absorb more water, and deeply nourish hair.

If we put detergent or hair conditioner into then SPA liquid, they will combine with micro bubbles which have negative charges and cling to the hair root, thus improving damaged hair and improving overall feel and quality.

Micro bubbles massage the skin and coat. They can also deliver medications and conditioners to the desired areas needing treatment.

The effect will be immediate after one or two treatments, and cumulatively




better after regular ongoing treatments. Extensive tests have been conducted to ensure that our micro bubbles are effective. For allergies, skin diseases and more severe skin conditions and rashes, adding pure oxygen and Veterinarian recommended and prescribed serums and medications can enhance the SPA treatments.

Submerge your hand into the SPA for a minute, you'll realize immediately the benefits this wonderful tub will provide for your pets. There have been proven results of the positive dramatic effects the SPA treatment has on human skin as well. Therefore your pets will surely benefit from the treatment.



Glossy Stainless Steel Bathing Tub | BTS-161

Mirror looking fancy bathing tub.
 High polished smooth surface realizes both decent appearance and easy-to-shine function in one go.
 Outer dimension at 900mmLx550mmWx900mmH with inner tub at SPA version available with item number of BTS-161SPA, with micro-bubble pump built-in.

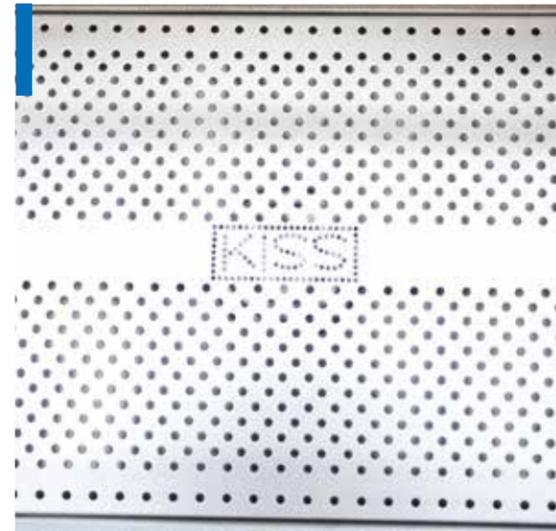



Stainless Steel Bathing Tubs BTS-130 / 131 / 130E / 131E

Made of SUS304 stainless steel, perfect manufacturing, classic constructure, occupying European and American market for years.
 Selected stainless steel from famous manufacturers meets or exceeds international standards.
 13 perfect manufacturing processes---shearing, cutting, welding, polishing.....
 Standard parts: Multi-filtering drain and exclusive drain hose, standard predrilled hot/cold faucet, ramp(standard version), water resistant height-control pedal(for E series only), etc.
 Optional up-down lifting door, left-right sliding door, electric lifting base, 2.0mm extra thick stainless steel material, ramp with stairs, ramp without stairs....multiple structure choices.



Classic Bathing Tub with Lifting Door | BTS-130

Made of SUS 304 stainless steel, up-down lifting door and ramp.



Classic Bathing Tub with Sliding Door | BTS-131

Made of SUS 304 stainless steel, left-right sliding door and ramp.



Electric Lifting Bathing Tub with Lifting Door | BTS-130E

Made of SUS 304 stainless steel, up-down lifting door and large-lifting-range base.



Electric Lifting Bathing Tub with Sliding Door | BTS-131E

Made of SUS 304 stainless steel, left-right sliding door and large-lifting-range base.



Deluxe Stainless Steel Tub | BTS-130T

Wholly made of 2.0mm extra thick SUS304 stainless steel, durable and stable.




Stainless Steel Bathing Tubs with Swing in/out Ramp | BTS-131P

Made of SUS304 stainless steel, including left-right sliding door and push-and-pull ramp with stairs for easy storage.



Stainless Steel Bathing Tubs With Swing in/out Ramp | BTS-130ECO

Made of SUS304 stainless steel, simple assembly, easy to operation, good value for money.

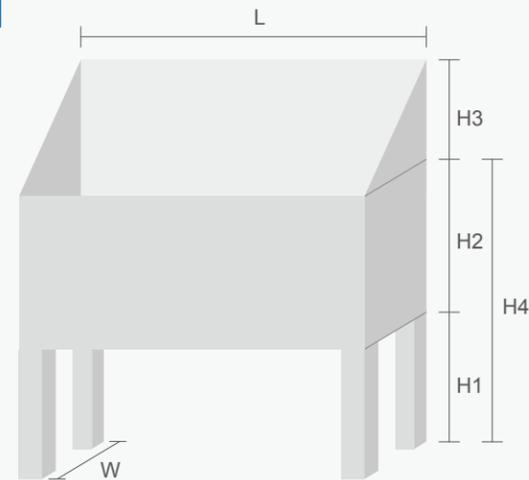

Dimension Figures

H1: Height of the base/leg.
 H2: Internal height of the tub
 H3: Height of the splash board.
 H4: Height from the opening of tub to the ground.



BTS-130T

Detachable and easy to pack



BTS-130E/131E



Tie-down ring on back splash board



Puppy booster



Black plastic floor panel is not equipped Standard which is an optional part



Hose holder



Water faucet reserved holes



Stainless steel drainage And hair strainer



Water faucet



Optional part: SR-401 shampoo rack



Simple operation



BTS-130		BTS-131		BTS-130E		BTS-131E	
L: 127cm	N.W: 53kg	L: 127cm	N.W: 53kg	L: 127cm	N.W: 92kg	L: 127cm	N.W: 92kg
W: 67cm	G.W: 90kg	W: 67cm	G.W: 90kg	W: 67cm	G.W: 125kg	W: 67cm	G.W: 125kg
H1: 50cm	Volume: 0.72m ³	H1: 50cm	Volume: 0.72m ³	H1: 25-95cm	Volume: 0.89m ³	H1: 25-95cm	Volume: 0.89m ³
H2: 38cm	Package: 140x75x69cmH	H2: 38cm	Package: 140x75x69cmH	H2: 41cm	Package: 140x76x84cmH	H2: 41cm	Package: 140x76x84cmH
H3: 52cm		H3: 52cm		H3: 40cm		H3: 40cm	
H4: 88cm		H4: 88cm		H4: 66-136cm		H4: 66-136cm	

BTS-130T		BTS-131P		BTS-130ECO	
L: 127cm	N.W: 118kg	L: 127cm	N.W: 53kg	L: 127cm	N.W: 53kg
W: 67cm	G.W: 155kg	W: 67cm	G.W: 90kg	W: 67cm	G.W: 90kg
H1: 50cm	Volume: 0.81m ³	H1: 50cm	Volume: 0.72m ³	H1: 50cm	Volume: 0.72m ³
H2: 38cm	Package: 140x78x74cmH	H2: 38cm	Package: 140x75x69cmH	H2: 38cm	Package: 140x75x69cmH
H3: 52cm		H3: 52cm		H3: 52cm	
H4: 88cm		H4: 88cm		H4: 88cm	



Super Ergonomic Bathing Tub | BTS-139

Instead of wrestling with large dogs, they can walk into the tub at almost floor height with great ease. Not bend-over and struggling with bad postures all day long, you can adjust both the tub and inside rack height to achieve the optimum working height you need.

The inner platform could lift animals during bathing and drying works, even when they lay down for temporary rest, you can still adjust their height with gentle start and stop actions.

Super big loading capacity up to over 300kg, ensures enough power to lift up any dog on this planet.

Water resistant mechanism, the linear actuator and power system can still perform well under all day long high humidity and even water splashing.

Fully opened front door with buffering system, allows every pet to find their most comfortable way to get in and out.

Constructed from 304 grade stainless steel, easy to clean, robust and up to the daily challenges of a busy grooming salon.

Personalized painting on splashing boards available, please consult your local distributor for more details.





Characterized painted BTS-136E

Stainless Steel Bathing Tub | BTS-136/136E

A budget tub with fully welded tank and assembly splash guards. Very good choice for grooming salons and boarding kennels; a terrific value. Puppy booster, bottom panels, hair strainer, drain, flexible drain hose and silicon caulking kit included.

It's time, to add lovely colors to your bathing room! Get away from the dull, cold and boring stainless steel looking, make it pretty, make it "warm", and make it yours! The painting on the inside of the splashing board of BTS-136 can be characterized to meet your personal demand. The painting is vivid with high defination.



BTS-136



BTS-136E



Water Sprayer Holder



Black plastic floor panel



Accessories



Grooming technology

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand	W×D×H (inch) Package	
BTS-136	35	45	125×65×138	49-7/32×25-19/32×54-21/64	129×69×44H	0.39
BTS-136E	62	104	125×65×103-173	49-7/32×25-19/32×40-35/64-68-7/64	130X68X70H	0.62



BTN-140

Bathing tub

| BTN-140/141

A new choice for bathing tub, with large capacity and attractive appearance.

The front entry allows convenient washing from both sides of the tub.

Very good option for large breeds, especially when the animal is not easy to turn around physically.

Two options for interior tub, either made from ivory color fiberglass or white color powder coated stainless steel.

Side panels could be painted to customer choice.

All parts available, including stainless steel restraint ring, spray nozzle holder, Moen shower faucet, and chrome plated sprayer hose or vinyl hose, etc. (they are not included in standard package, please let our sales representatives know the requirements).

Exterior measures 142cm long, 88cm wide, 91cm tall.

Interior size is 122cm long, 72cm wide, 50cm deep.




Entry step



Round Cornered Fiberglass Tub



Hair Strainer

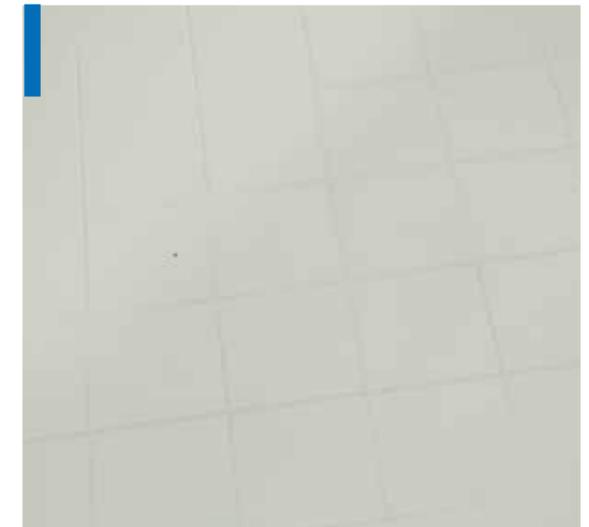

Anti-slippery Floor



BTN-141



BTN-141



Inline Stainless Steel Tub | BTS-137E

Great solution for limited space could be lined up in a row for crowded rooms. Easy and comfortable lift door with low-low stable electric base leaves easiest access for large breeds. Robust electrical base, lifts up to 250lbs as maximum capacity. Stainless steel drainage, hair strainer and plastic standing floor panels included. No assembly needed, shipped in solid one piece.



Black plastic floor panel



Stainless steel drainage and hair strainer

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume (cbm)
	(kg)	(kg)	Stand	Package	
BTS-137E	77	110	126.5×67×65-131	140×76×84H	0.89



BTS-132TE

European Style Stainless Steel Bathing Tub | BTS-132T/132TE

European Styling.

BTS-132T has a fixed base.

BTS-132TE has the same smart design with an electric base.

Constructed from 13 gauge stainless steel panels. This tub will perform for years to come.

Wide front door opens across the full width of the tub, allowing large breeds easy access.

Large breeds and older dogs can stroll in and out (of the BTS-132TE) when the tub is lowered.

Dog fixture rings on back splash guard for easy dog restraint.

Comfortably bath small breeds without bending by over raising the tub to a suitable level for your stature (and/or by adding a booster rack).

Water resistant foot pedal controls make operating the tub simple.

The gentle movement up/down does not startle pets.

Tubs also come with stainless steel hair strainer, flexible drain hose.

Tub bottom panel, puppy booster and faucet are optional parts.



Extra low access for dog to walk on



BTS-132T

BTS-132T free standing style



Stable electric lift base with great Loading capacity



Fully opening front door



Black plastic floor panel is not equipped standard, which is an optional part.



Accessories Included

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH(cm)	Stand WxDxH (inch)	Package	
BTS-132T	50	53	120x60x140	47-1/4x23-5/8x55-1/8	123x63x88H	0.35
BTS-132TE	77	83	118x64x105-175	47-1/4x23-5/8x41-11/32~68-57/64	123x63x88H	0.58
					129x69x25H	



Stainless Steel Hydro Bath | BTS-145

World Class Hydro Bath
(performs double duty as a therapeutic massage treatment tub).


High volume water pressure helps generate increased lather, which in turn saves time when washing pets, consuming less shampoo.

This tub employs top rated, Australian produced, Davey Pumps.

Powerful water stream makes the product not merely a tub but also a SPA machine, providing a gentle and thorough full body massage for your pets.

Two recirculation tanks offer options of solution treatment; fresh water rinsing available too. Advanced temperature and water level monitoring system ensure 100% safe usage.

Water circulation system

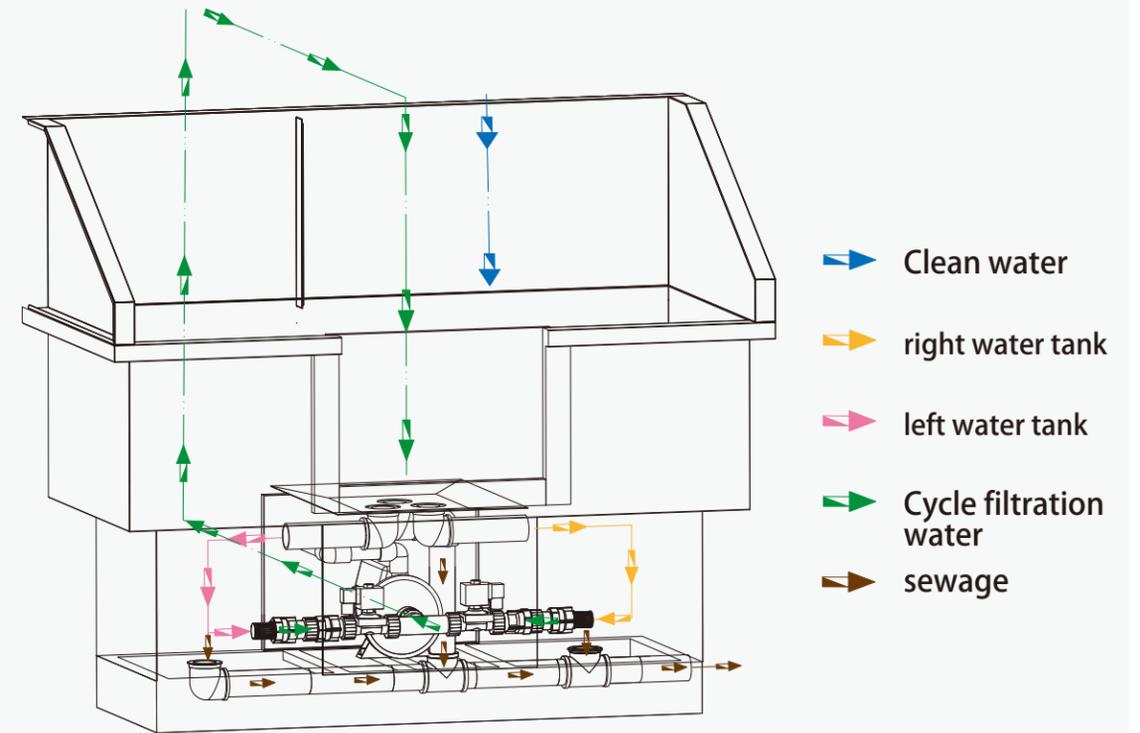

Industrial strength hinge



Recirculation plugs



Removable grids



Increased lather production



Water inlet filter



Instant on heater



BTS-142

Round Lip S/S Tub

| BTS-142/142E

Round corner front lip for comfortable body contact.
 Lift-up door realizes easy entrance for large sized animals.
 Built from 100% 304 stainless steel with whole body fully welded for great durability.
 Pre-drilled faucet holes adapts all the existing faucet center distances.
 Electric lifting version BTS-142E available too.






Stainless Steel Electric Lifting Bathing Tub

| **BTS-138E**

Made of SUS 304 Stainless Steel, anti-rust with lifelong durability.
 Fully welded bathtub main body, with round corners and smooth polishing, the best solution for cleaning and keeping hygienic.
 Good for containing water and providing medical treatment for pets with skin problems.
 Pre-drilled holes on each side allow for 4" to 8" faucet installations.
 Foldable Splash guard panels offer flexibility based on your needs.
 Easy height control by using step-on water resistant pedals.
 Simple management of small breeds (by raising the tub and adding puppy booster) and larger /older breeds (by lowering the tub for them to stroll in).
 The gentle movement up/down does not startle pets.
 Also comes with flexible drain hose, stainless steel hair strainer, two tub bottom panels and one puppy booster.

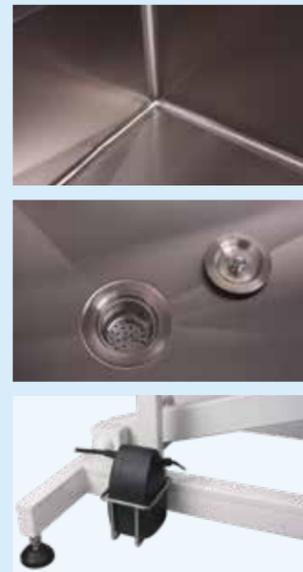

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume
	(kg)	(kg)	Stand	Package	(cbm)
BTS-138E	89	122	127×67×106-176	140×76×84	0.89



BTS-85ECO

Small Animal Bathing Tub | BTS-85/85ECO

Overall size:
 BTS-85: 840mm×680mm×1200mmH (33"×26-3/4"×47-1/4"H)
 BTS-85-ECO: 900mm×500mm×1150mmH (35-1/2"×19-3/4"×45-1/4"H)
 Small animal bathing tub, saving precious grooming room space. Budget solution for small sized stainless steel bathing tub. Stable freestanding base made from fully welded stainless steel tubes. Completely welded tub main body prevents leakage, working well for merging animals into water. Pre-drilled mounting holes for faucet installation, adaptable for both 6" centered or 4"/8" faucets, please clarify this point when ordering. Includes: Splash boards, drainage, flexible waste pipe, sprayer holder, and tie-down ring.



Stainless steel drainage and hair strainer

BTS-85

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume
	(kg)	(kg)	Stand	Package	(cbm)
BTS-85	84	96	68×84×120	88×66×135	0.92
BTS-85ECO	84	96	50×90×108	95×52×30	0.26



Small Fiberglass Tub

| BTF-90

Nautical grade fiberglass, the same material and thickness as used in yacht construction makes this product a top quality tub.

Fully sealed tub tank, made from injection mold, ensures great durability and reliability.

Softer touch and warmer feel than regular stainless steel tubs, better care for animals.

BTF-90 is equipped with firm and stable freestanding base, BTF-90E comes standard with low-drop electric base.

Includes: stainless steel drainage, hair strainer, and flexible waste pipe.



Skid-proof and decorative solid figure on the bottom



Simple and firm connection



Fashionable and beautiful streamline



ITEM	Net Weight		Gross Weight		Dimension (cm)		Volume (cbm)
	(kg)	(kg)	(kg)	(kg)	Stand	Package	
BTF-90	82	97	60×90×110	65×95×110			0.36



Round Cornered Economical Bathing Tub

| BTS-90-ECO

Measures at 900Wx500Dx900H(front lip), a really compact but feasible small size tub.

100% made from 304 stainless steel with fully welded and polished radius corners, it's really sturdy and durable.

Fully sealed tank, good for SPA and medical bathing treatment.

Knock-down package, cost effective for transportation.





Pet Washer

| PW-001

A perfect tool to provide a more effective bath in animal clinics, grooming salons, pet clubs, animal shelters and kennels.

It can help save 20% shampoo or conditioner and make groomers work easily.


The perfect combination of air and water flow, eliminating the cumbersome traditional manual dilution of shampoo; Automatic shampoo dilution without electricity and extra cost;

Labor-saving, time-saving, and "worry-saving"; Perfect dilution produces quick foaming and thorough cleaning;


Save shampoo by more than 20%, no need to worry about excessive waste of shampoo;

4 optional shampoos;


Easily setting dilution ratio for different shampoos.

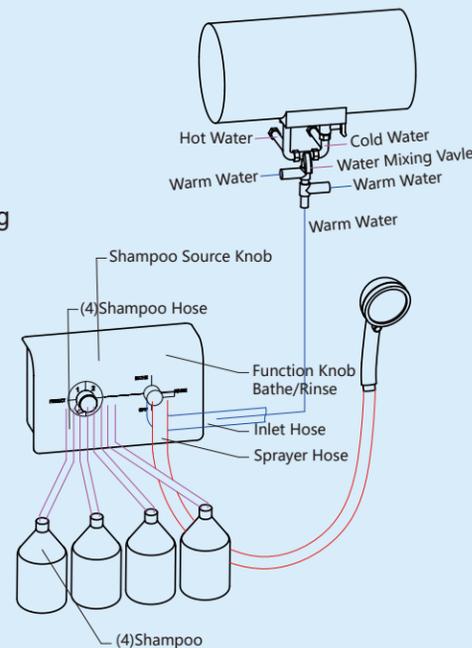

Optional Parts

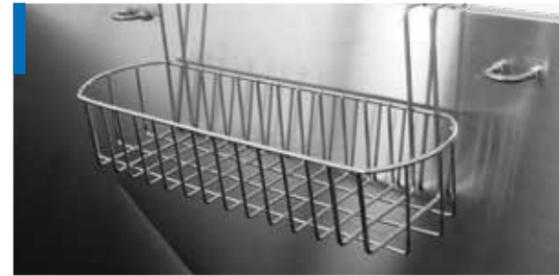

Shampoo Rack
SR-401

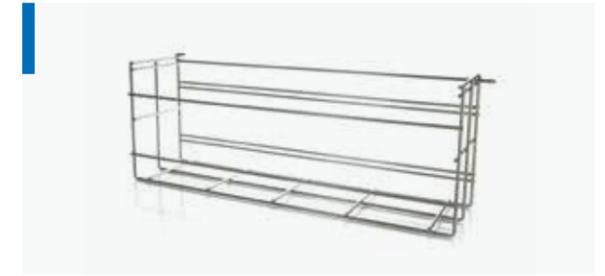

Powder Coated Shampoo Rack - Holds 4 Gallons
SR-402



Water Faucet



S/S Water Sprayer
WS-003



Overhead Dog Holding Arm

Standard Parts



Black Plastic Floor Panel



Stainless Steel Drainage and Hair Strainer



Hose/Sprayer Holder



Tie-down Ring On Back
Splash Board




Flexible Waste Pipe


Classification	Material		Standard parts and functions										
	Model	Stainless Steel	Other material	Fixed splash-board	Foldable splash-board	Ramp	Lifting	Puppy booster	Bottom panel	Drain assembly	Adjustable leveler	Standard predrilled threaded holes	Screws
BTS-150	☆	Acrylic Tub					☆	☆	☆				
BTS-151		Acrylic Tub						☆	☆				
BTF-180E		Nautical grade fiberglass	☆			☆			☆		☆	☆	
BTS-130	☆			☆	☆		☆	☆	☆	☆	☆	☆	
BTS-131	☆			☆	☆		☆	☆	☆	☆	☆	☆	
BTS-130E	☆			☆	☆	☆	☆	☆	☆		☆	☆	
BTS-131E	☆			☆	☆	☆	☆	☆	☆		☆	☆	
BTS-130T	☆			☆	☆		☆	☆	☆	☆	☆	☆	
BTS-131P	☆			☆	☆		☆	☆	☆	☆	☆	☆	
BTS-130ECO	☆			☆	☆		☆	☆	☆	☆		☆	
BTS-137E	☆						☆	☆	☆				☆
BTS-132T	☆		☆					☆	☆	☆			☆
BTS-132TE	☆		☆			☆		☆	☆				☆
BTS-136	☆		☆				☆	☆	☆	☆	☆	☆	
BTS-136E	☆		☆			☆	☆	☆	☆		☆	☆	
BTF-90		Nautical grade fiberglass	☆						☆	☆			☆
BTS-85	☆		☆						☆	☆	☆	☆	
BTS-85ECO	☆		☆						☆	☆	☆	☆	

Classification	Package				Optional parts					
	Model	Packless	Simple package	Wooden case	Carton	Optional size	Hot/cold faucet	Ramp	Grooming arm	Shampoo storage basket
BTS-150	☆			☆		☆		☆		
BTS-151	☆			☆						
BTF-180E			☆	☆			☆			
BTS-130			☆	☆		☆	☆	☆	☆	☆
BTS-131			☆	☆		☆	☆	☆	☆	☆
BTS-130E	☆			☆		☆	☆		☆	☆
BTS-131E	☆			☆		☆	☆		☆	☆
BTS-130T			☆	☆		☆	☆		☆	☆
BTS-131P			☆	☆		☆	☆		☆	☆
BTS-130ECO			☆	☆			☆		☆	☆
BTS-137E	☆			☆						☆
BTS-132T			☆		☆					☆
BTS-132TE			☆		☆					☆
BTS-136			☆		☆		☆		☆	☆
BTS-136E			☆		☆		☆		☆	☆
BTF-90			☆		☆		☆			
BTS-85	☆			☆		☆	☆			☆
BTS-85ECO			☆		☆		☆			☆


BTS-150

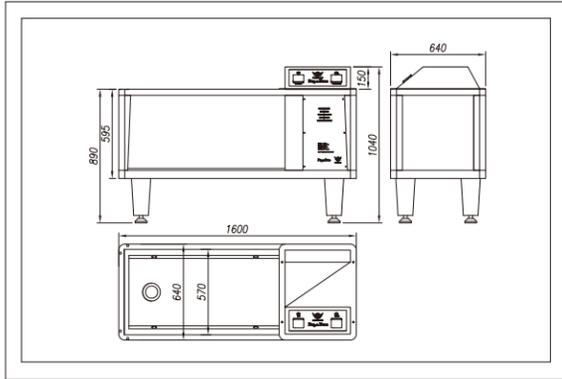


BTS-151

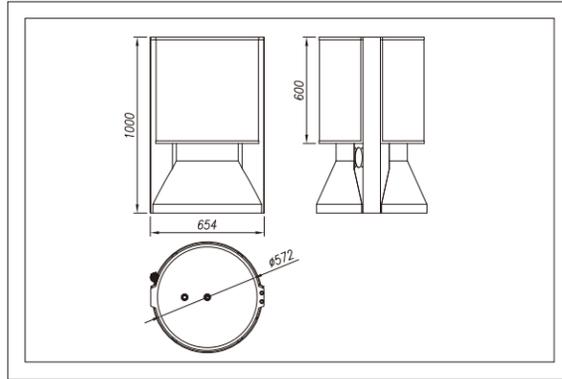


BTS-136E

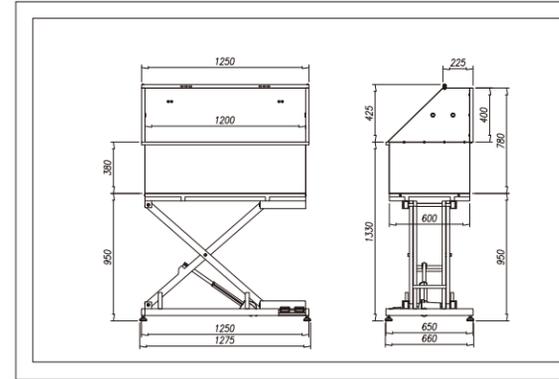


BTS-136E

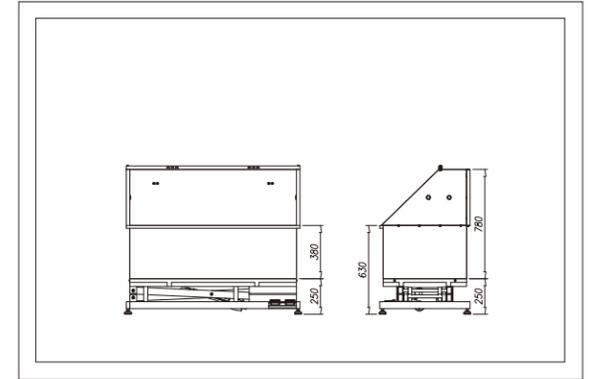


BTS-130/131

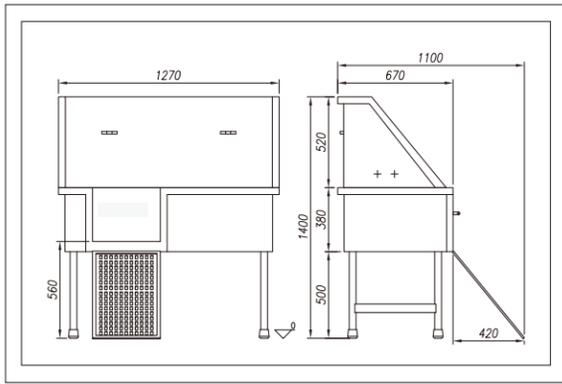


BTS-132T

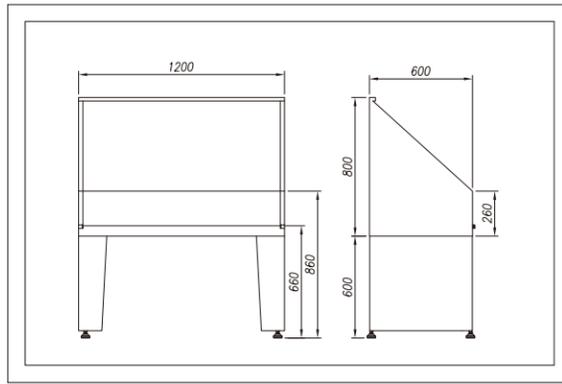


BTS-132TE

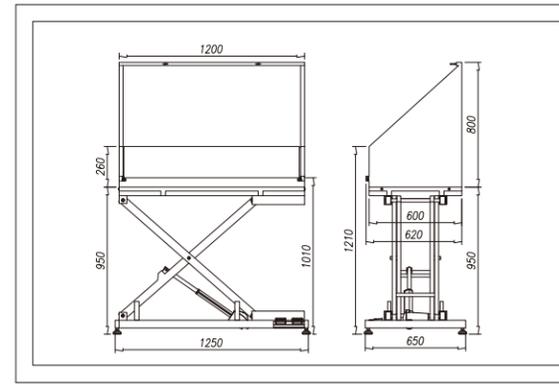


BTS-132TE

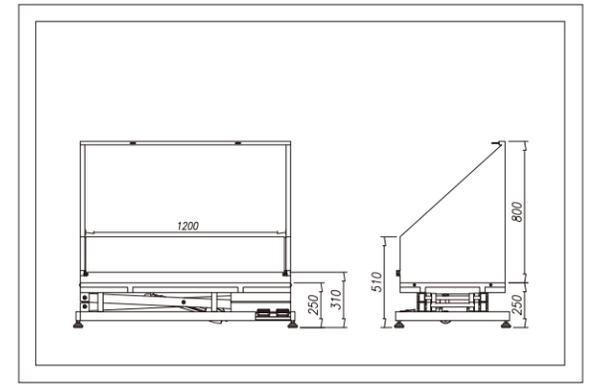


BTS-180E

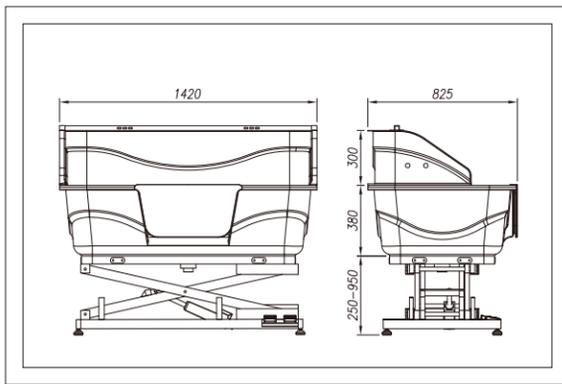


BTS-136

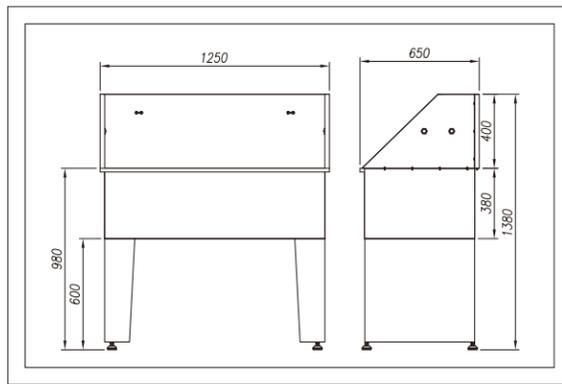


BTS-130/131E

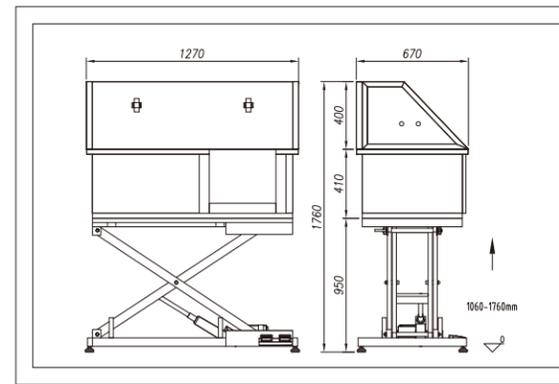


BTS-130/131E

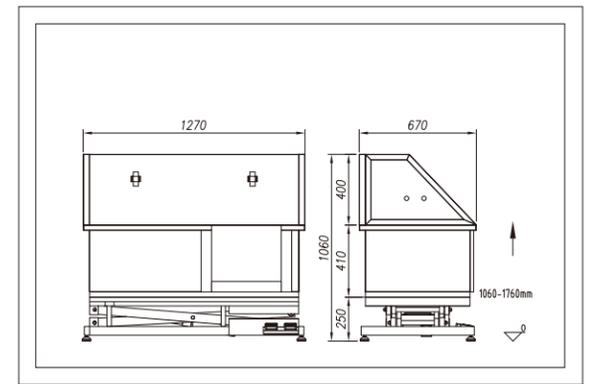


BTF-90

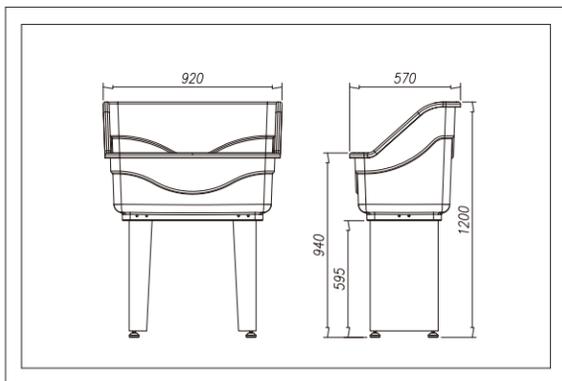


BTS-85ECO

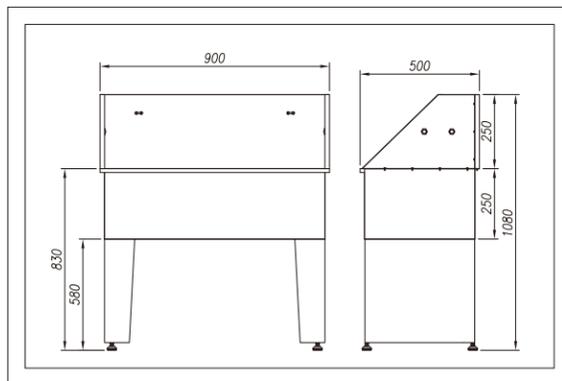

BTS-138E

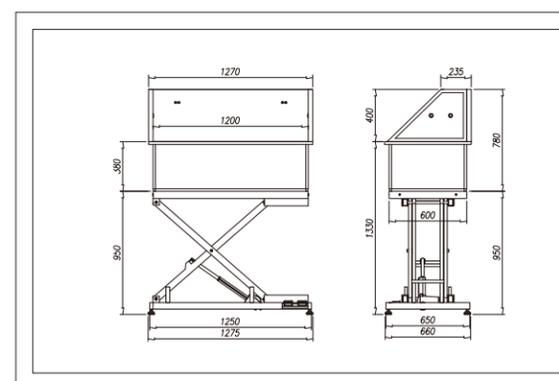

BTS-138E

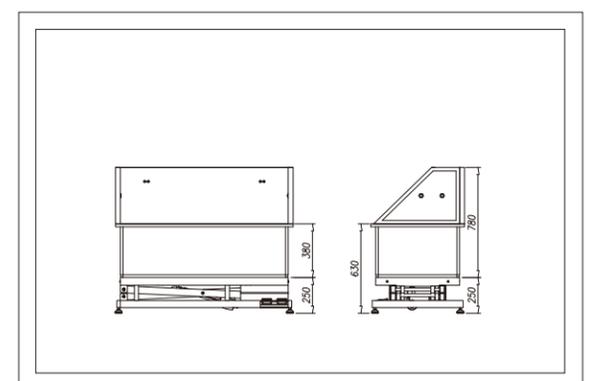



Rechargeable Pet Trimmer | MC-220

- Super lightweight, rechargeable with long running time.
- Safe and close-to-skin blade ensure smooth cutting at hard-to-reach areas.
- Ideal handy trimmer for ears, faces, paws and sanitary areas.
- T-liner blade available on this trimmer, good for creative carving too.
- Mini-USB charging port, could be charged by any standard power banks.



Rechargeable Pet Trimmer | MC-230

- Though slim and tiny, but as powerful as all other big trimmers.
- A perfect solution for partial trimming needs.
- Very light weighted and easy to carry around.
- Good in pocket, perfect in palm.
- Insanely low noise, not annoying pets at all.
- Safe, reliable, stable, durable and efficient.



LED display of battery capacity and RPM High-tech sense Detachable blade Good in pocket, perfect in palm.



Pet Clipper

MC-670

Heavy duty 12V motor built-in, born for toughest coats!
 Slim and lightweight mainframe fits any sized hands comfortably.
 5 detachable blades available, leaving hair length from 3 to 12mm.



Rechargeable Pet Clipper

MC-680

Four detachable blades and four comb attachments make easier to trim different pets' hair.
 LCD displays motor speed, battery percentage and out-of-oil.
 Quiet running comforts pets.
 2000mA high performance battery ensures 4-hour working time.





High Torque Rechargeable Clipper | MC-690

High torque, quiet running motor ensures very smooth cutting.
 Stainless steel blade, standard and ultra thin two sizes available.
 Metal combs available for whole body trimming finish.
 Lithium battery with 4-5 hour working time.



CL-SERIES

CL-530

Two Ceramic Blades
 One Battery
 Comb Set




CL-SERIES

CL-6110

Two Ceramic Blades
 Two Batteries
 Plastic Hard Case



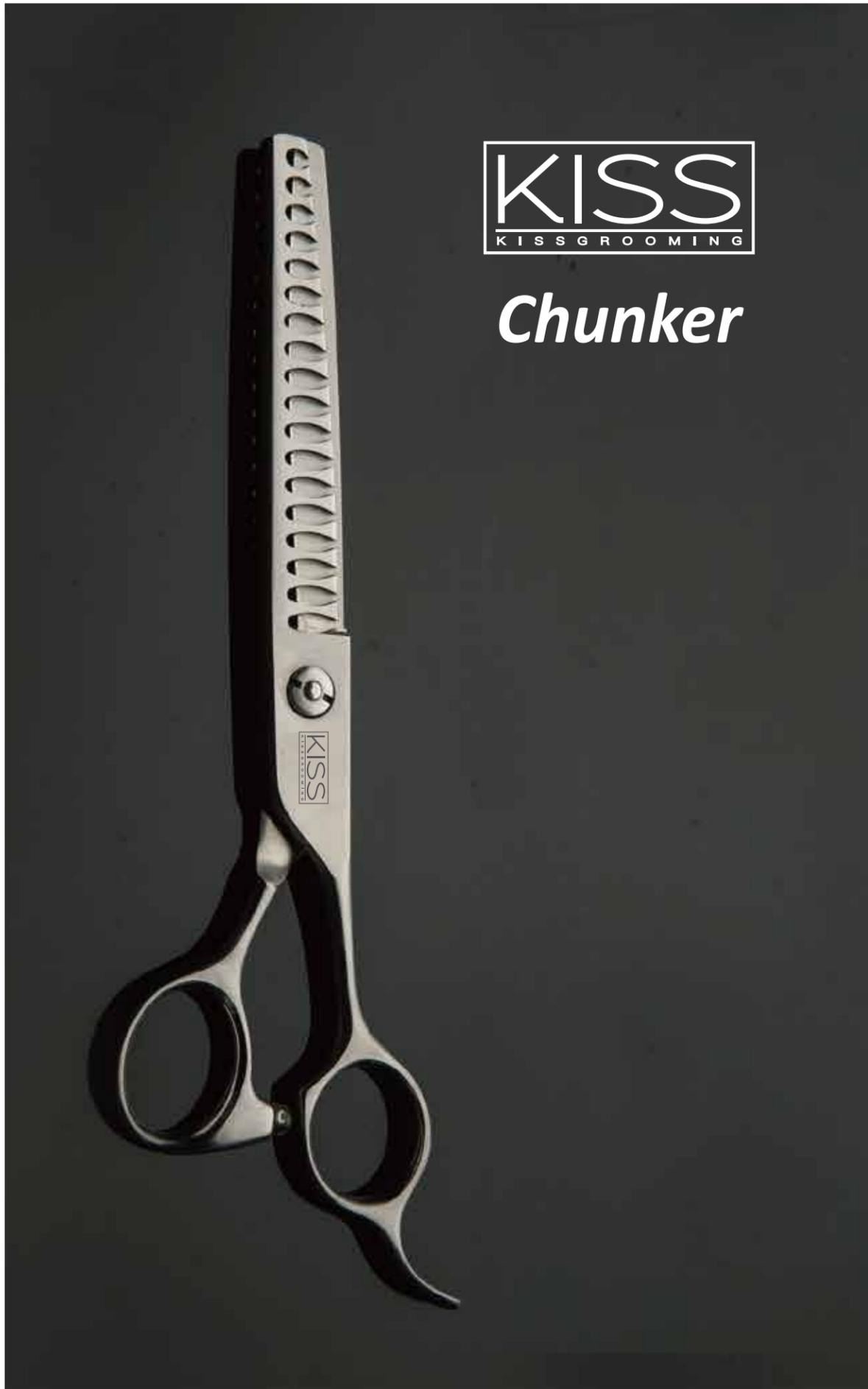

KISS
KISSGROOMING

Chunker



KISS Straight Scissor / Curved Scissor





The KISSGROOMING line offers a broad spectrum of shears. Ideal choice for busy grooming salons. Each shear is made of top performance handmade Japanese steel and ensuring long term durability and performance.



On a par with other top of the line shears, The Rainbow and Deep Space Black colored Prism series provides superior performance. It is made of titanium coated high quality Japanese 440C stainless steel. It offers a fantastic value and style. The popular choice for champion groomers the world over.





Velvet

The crème de la crème of shears.

The Velvet series is one of the highest quality lines in the market. Master groomers choose these shears over others prices considerably higher. This series is made of high quality Japanese 440C stainless steel with hand polished finish. The workmanship is exceptional. It holds it edge longer than most. The modern and ergonomic shear handle design features the “crane grip”, designed to help relieve stress on both shoulders and wrists.


Brand	Material	Surface Treatment	Type	Model	Color
	Japanese 440C stainless steel	Hand polished finish	6.5 inch teenth thinner shears	SV21265	Metallic natural luster
			7.5 inch teenth thinner shears	SV21275	
			4.5 inch facial straight shears	SV21145	
			5.5 inch facial straight shears	SV21155	
			7.0 inch straight shears	SV21170	
			7.5 inch straight shears	SV21175	
			8.0 inch straight shears	SV21180	
			8.5 inch straight shears	SV21185	



Swan

The Swan series is the shear for fashionable groomers who demand performance and value. Each shear is made of stainless steel and painted in 5 gorgeous colors; **Barbie Pink**, Lavender, **Emerald**, Sky Blue, **Sunflower Yellow** and Space Black.


Brand	Material	Surface Treatment	Type	Model	Color
	Stainless steel	Painted	6.5 inch teenth thinner shears	SW5 X 265	X: 1-Barbie pink 2-Lavender 3-Emerald green 4-Sky blue 5-Space dark
			7.5 inch teenth thinner shears	SW5 X 275	
			7.5 inch straight shears	SW5 X 175	
			8.0 inch straight shears	SW5 X 180	
			8.5 inch straight shears	SW5 X 185	
			7.5inch curved shears	SW5 X 375	
			8.5inch curved shears	SW5 X 385	
			4.5 inch facial straight shears	SW5 X 145	
			5.0 inch facial straight shears	SW5 X 150	

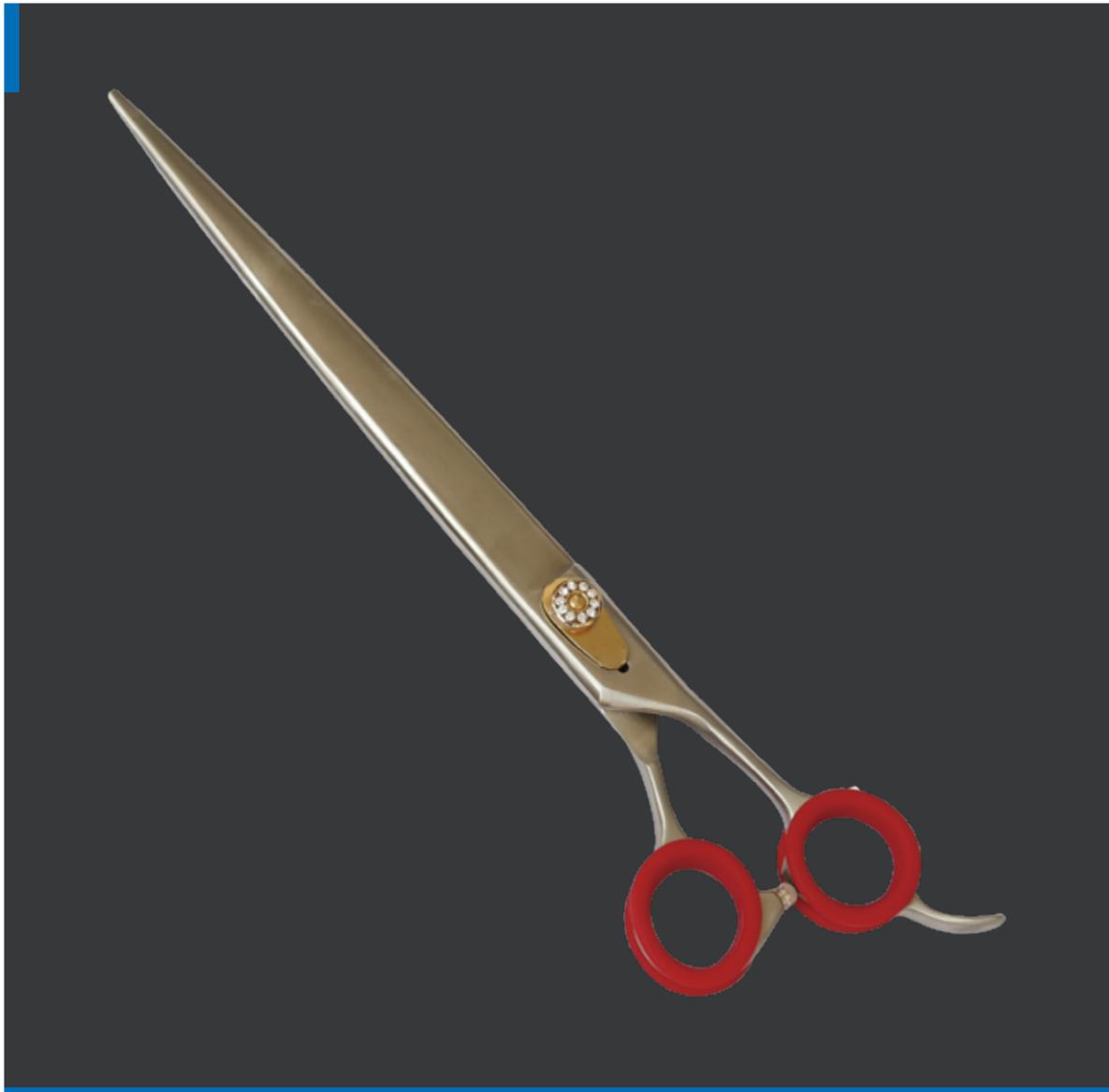

Scissors

Pro S

A combo set for one stop fast grooming solution!
 All tailor made from different portions on animal body to achieve systematic grooming operation with high efficiency.
 9" Straight, 8" Straight, 8.5" Curved, 6" Curved, 7" Double Sided Chunker, 6.5" Curved Chunker, 8" Chunker, forming a complete system for full body scissoring.
 For educational video guides, please search "Fast Jenny Grooming" on Youtube!

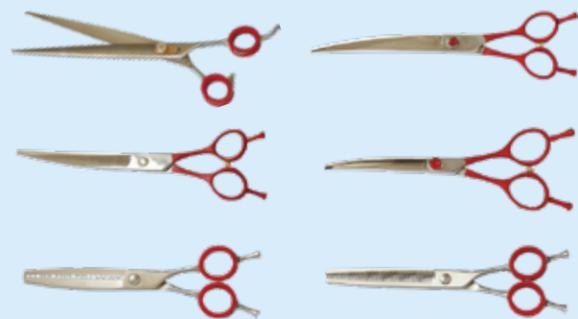

Basic Information

Series	Material	Coating	Type	Item No.	Color	
PRISM	440C from Japan	Titanium in Rainbow Color	6.5" Thinner	SP11265	Rainbow	
			7.5" Thinner	SP11275		
			4.5" Facial Straight	SP11145		
			5.5" Facial Straight	SP11155		
			7.0" Straight	SP11170		
			7.5" Straight	SP11175		
			8.0" Straight	SP11180		
	440C from Japan	Tatinium Coating in Black	6.5" Thinner	SP12265	Black	
			7.5" Thinner	SP12275		
			4.5" Facial Straight	SP12145		
			5.5" Facial Straight	SP12155		
			7.0" Straight	SP12170		
			7.5" Straight	SP12175		
			8.0" Straight	SP12180		
Velvet	440C from Japan	Hand Polished	6.5" Thinner	SV21265	Chrome	
			7.5" Thinner	SV21275		
			4.5" Facial Straight	SV21145		
			5.5" Facial Straight	SV21155		
			7.0" Straight	SV21170		
			7.5" Straight	SV21175		
			8.0" Straight	SV21180		
KISS <small>KISSGROOMING</small>	5 Star 9Cr13mov	Hand Polished	6.5" Thinner	SH35265	Chrome	
			7.5" Thinner	SH35275		
			4.5" Facial Straight	SH35145		
			5.5" Facial Straight	SH35155		
			7.0" Straight	SH35170		
			7.5" Straight	SH35175		
			8.0" Straight	SH35180		
	4 Star	Hand Polished	8.5" Straight	SH35185	Chrome	
			7.5" Curved	SH35375		
			8.5" Curved	SH35385		
			6.5" Thinner	SH34265		
			7.5" Thinner	SH34275		
			4.5" Facial Straight	SH34145		
			5.5" Facial Straight	SH34155		
7.0" Straight	SH34170					
Swan	Stainless Steel	Painted	7.5" Straight	SH34175	Chrome	
			8.0" Straight	SH34180		
			8.5" Straight	SH34185		
			7.5" Curved	SH34375		
			8.5" Curved	SH34385		
			6.5" Thinner	SW5 X 265		X: 1-Barbie Pink 2-Lavender 3-Emerald Green 4-Sky Blue 5-Space Black
			7.5" Thinner	SW5 X 275		
			7.5" Straight	SW5 X 175		
			8.0" Straight	SW5 X 180		
			8.5" Straight	SW5 X 185		
7.5" Curved	SW5 X 375					
8.5" Curved	SW5 X 385					
4.5" Facial Straight	SW5 X 145					
5.0" Facial Straight	SW5 X 150					



MATTI

As the first electric dematter in the market, it's a huge helper for stressful dematting work.

Compact, lightweight, easy to grip, reducing strain on groomers' hands, wrists and shoulders.

The metal safety guide rails prevent any discomfort for the dog by stopping the pins from coming into contacts with the dogs' skin.

Equipped with brushless motor, ensuring long lasting service life.

Adapting 100 -240V voltage worldwide.

Unique hair collector gathers all the loose hair together easily.

Replaceable cassette has S/S safe pins; it's easy to be replaced by few simple steps.



25KISS015



25KISS014

Kiss Slicker Brush

The Magic Tool for dematting!

Specially designed curved surface with heavy duty rubber mat mounted.

430 stainless steel pins implanted for comfortable and safe operation.

Removes mats, loose hair and debris easily.

Lightweight, eliminate the stress on your wrist.

Two sizes available.



Professional Brushes

Handles made from beautiful walnut, with wonderful texture, grain and natural colors. 3 different length pins, implanted tightly in durable rubber pad. Ionic brush, combined with fine cooper wires and nylon, generates rich ions that aid in helping to flatten coats and provide sheen.

Item Number	Description	Packing Size
BR-200	High Quality Thermo Ash Wood Brush Luxury Metal Pin Brush - 20mm Pin	23cm*24cm*7cm 0.5kg 5pcs/CTN
BR-270	High Quality Thermo Ash Wood Brush Luxury Metal Pin Brush - 27mm Pin	23cm*24cm*7cm 0.5kg 5pcs/CTN
BR-350	High Quality Thermo Ash Wood Brush Luxury Metal Pin Brush - 35mm Pin	23cm*24cm*7cm 0.6kg 5pcs/CTN
BR-B01	High Quality Beech Wood Brush Palisander Luxury Bronze Ionic Brush	26cm*24cm*10cm 0.6kg 5pcs/CTN
BR-B02	High Quality Beech Wood Brush Palisander Luxury Bronze Ionic Body Brush	10cm*16cm*7cm 0.24kg 1pcs/CTN



BR-200




BR-270



BR-350




BR-B01

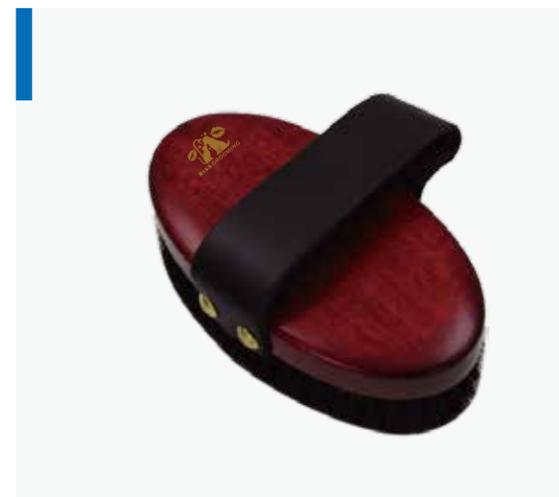

BR-B02



From Germany

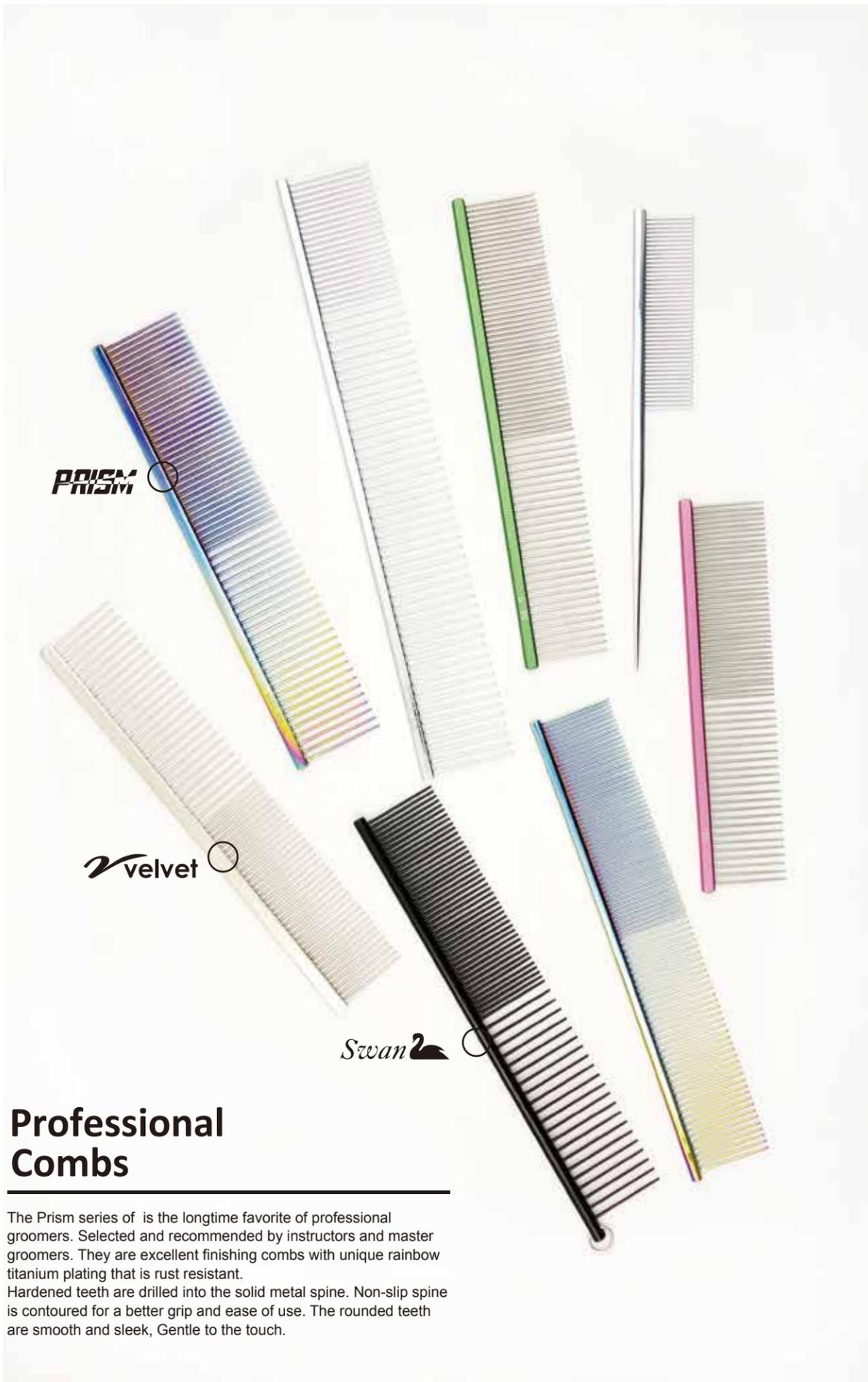

Professional Combs

The Prism series of is the longtime favorite of professional groomers. Selected and recommended by instructors and master groomers. They are excellent finishing combs with unique rainbow titanium plating that is rust resistant. Hardened teeth are drilled into the solid metal spine. Non-slip spine is contoured for a better grip and ease of use. The rounded teeth are smooth and sleek, Gentle to the touch.



Professional Combs

Basic Information				Size				Teeth/inch		Details	
Series	Coating	Feature	Item No.	Handle Length	Teeth Length	Medium	Coarse	Handle	Color		
PRISM	Rainbow Titanium Plating	Medium/Coarse	CMP-102	7-1/2"	190mm	1-1/16"	27mm	12	7	5mm Flat Spine	Rainbow
		Medium/Coarse Ultra Long Teeth	CMP-104	9-5/8"	245mm	1-1/4"	32mm	12	6,4	5mm Flat Spine	Rainbow
Velvet	Chrome Coating	Medium/Coarse Ultra Long Teeth	CMV-206	9-5/8"	245mm	1-1/4"	32mm	12	6,4	5mm Flat Spine	Chrome
Swan	Anti-static Teflon Coating	Medium/Coarse	CMS-401	7-1/2"	189mm	1-1/16"	28mm	12	7	5mm Round Spine	Black
		Medium/Coarse Long Teeth	CMS-402	7-1/2"	189mm	1-1/4"	32mm	12	7	5mm Round Spine	Black



Replica Model Dog

High simulation effect, easy to comb, nap, style and cut.
 Enough quantity and length of the hair to cut a professional competition hair style or cuddly Teddy style.
 Moderate intensity and hardness of the hair can avoid possible hand muscle strain, or any wear and tear of the scissors.
 Safety of the hair ensures no odor or allergen, and will not cause itchy or allergic skin, mouth and nose.
 Reasonably designed removable limbs ensure easy and convenient replacement of the hair from the model.
 Reasonable price of product exported to Europe will be offered to domestic groomers and groomer school.



Combing effect



Greatly matching head shape



Enough quantity and length of the hair



Cutting effect



New Headset Coming Soon



Item Number: RMD-23B1
Prouct Nanme: Head Hair(White)



Item Number: RMD-23B2
Prouct Nanme: Head Hair(Gray)



Item Number: RMD-22B1(White)
Prouct Nanme: Teddy Bear Full Body Coat



Item Number: RMD-22B3(Brown)
Prouct Nanme: Teddy Bear Full Body Coat



Item Number: RMD-23B3
Prouct Nanme: Head Hair(Brown)



Item Number: RMD-4101
Prouct Nanme: Head Set Mannequin



Item Number: RMD-22B2(Gray)
Prouct Nanme: Teddy Bear Full Body Coat



Item Number: RMD-12A1
Prouct Nanme: Teddy Bear Mannequin



Item Number: RMD-12C1
Prouct Nanme: Short Full Body Coat



Item Number: : RMD-12B1
Prouct Nanme: Standard Full Body Coat



Item Number: RMD-32A1
Prouct Nanme: Bichon Frise Full Body Coat



Item Number: RMD-3101
Prouct Nanme: Bichon Frise Mannequin



Item Number: RMD-1101
Prouct Nanme: Toy Poodle Mannequin



Item Number: RMD-12A1
Prouct Nanme: Long Hair Full Body Coat



Item Number: RMD-5101
Prouct Nanme: Torso Set Mannequin



Item Number: RMD-53B1
Prouct Nanme: Torso Coat

Grooming Apparels



 GA-030106



 GA-030601



 GA-030203



 GA-040106



 GA-040601



 GA-040203



 GA-030302



 GA-030405



 GA-030504



 GA-040302



 GA-040405



 GA-040504

Angelia Grooming Shirt

A trendy shirt with comfortable crew neck and short zipper.
Anti-static, breathable nylon fabric with water-resistant treatment
Adjustable back belt, set-in sleeve, 2 useful front pockets.

Eirene Wrap Shirt

Fashionable V-neck wrap shirt, constructed from ultra light nylon fabric, easy to clean, dries quickly. A practical tie belt (matching trim color), allows for easy size adjustment
Two large pockets with elastic are very functional for keeping tools close at hand.

Grooming Apparels



GA-050106



GA-050601



GA-050203



GA-060103



GA-060601



GA-060301

Pheme Grooming Jacket

Zippered front jackets features high collar and stylish fitted back. Fashionable classic styling, decorated with matching pocket and zipper trim. Reinforced at all stress points for long wear.



GA-050302



GA-050405



GA-050504



GA-080100



GA-080108



Pistis Grooming Jacket

Elegant fitted jacket with a high collar and a full length zipper. Long, for maximum protection and coverage. Adjustable back belt, set-in sleeve and covered front pockets.

Hebe

Unisex short cut styling. Water repellent fabric, made from 100% nylon, anti-static and easy to clean. Decoration choices, Schnauzer, Poodle, Bulldog, Scissor, Paw, etc.



GA-910103

GA-920101

Hemera

Hair Band

This useful hair band will keep you looking smart while working. Solid colors matches our apparels; with elastic on the back.

EOS

Tapered Bottom Pants

Comfortable pants made from 100% nylon, specially treated to repel hair, water and chemicals. Elastic waistband with drawstring.

Ceres

Straight Bottom Pants

Creates nicest style with any of our shirts or jacket. Very light weight, anti-static and easy to clean, dries quickly. 4 big pockets, two of them with elastic.

Color													
No.	00	01	02	03	04	05	06	07	08	09	10	More Colors for you	

Galene Footwear

Stylish, ventilated, comfortable, soft and light weight. Deep grooved outer sole is skid resistant, comfortable foot beds, swivel back-strap. Waterproof and easy to clean. Many color combinations.

00	01	05	06	08	10	11	12	13	14	More Colors for you
----	----	----	----	----	----	----	----	----	----	---------------------

GA-821108

GA-811200



Color					
European Size	36	37	38	39	40
US Size	5.5	6.5	7.5	8.5	9
Color					
European Size	40	41	42	43	44
US Size	7	8	9	10	11



GA-2113



GA-2214



GA-2111



GA-2112

Professional Grooming Apparel

Our specially coated material is waterproof, dustproof, germproof, antistatic and ventilate well.

3 quarter sleeves and sleeveless patterns. Custom decoration available on large orders, 100 pieces per color per style MOQ. Adjustable fit.

Item	Description
GA-2111 / 2112 / 2113 / 2114	Medium Size With Sleeves Pink / Silver / Blue / Purple
GA-2211 / 2212 / 2213 / 2214	Medium Size Without Sleeves Pink / Silver / Blue / Purple
GA-2121 / 2122 / 2123 / 2124	Large Size With Sleeves Pink / Silver / Blue / Purple
GA-2221 / 2222 / 2223 / 2224	Large Size Without Sleeves Pink / Silver / Blue / Purple



GC-003

Deluxe Saddle Stool

GC-003

This workhorse stool will help you get through all your seated work. Encouraging you to sit properly, it helps to avoid and relieve back strain and fatigue.

A beautiful, large ergonomic saddle shape can accommodate most people comfortably.

Long lasting top pad made from injected plastic, infused with shock absorbing soft material.

Rigid 5 pointed star shaped base made from zinc die-casting, solid, one piece construction, anti-rusting and durable.

Foot rest ring adds comfort.

Broad height adjusting ranges from 70cm to 94cm.



GC-003

Super big height adjusting ranges from 65cm to 89cm.




Durable top pad made from injection mold.




High-class hydraulic pump




Chrome plated feet rest ring

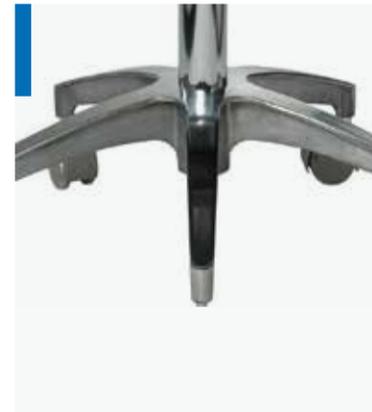

Decent overall finish

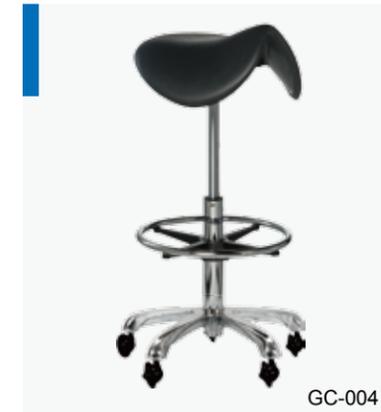

GC-004

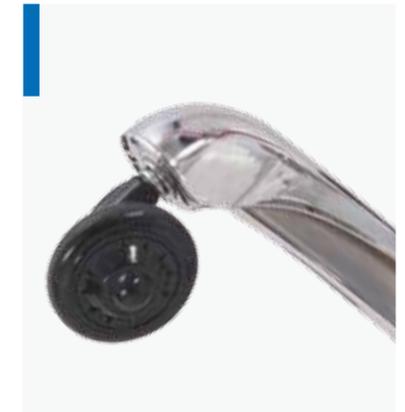

Top quality casters



Saddle Stool

GC-005

Fully plastic injection mold made stoop pad, ensures a super long service life without any concerns of delaminating. Ultra big sitting pad, a good choice for long termed operation, such as eye or dental surgeries. Zinc die-casted caster base, with super strength to support most operators. Smooth and comfortable foot rest ring, let your lower body rest in ease. Newest planetary casters rolls smoothly with little chance of hair jam inside a caster.

Ergonomic Saddle Chair

GC-002



Work chair with backrest

GC-001





TK-605-BK

TK-605-PK

Aluminum Tool Case

TK-605

Multiple combinations for using occasions.
 Free combinations: single-layer, double-layer, three-layer;
 hand-held, trolley-style. Reasonable space design, safe and efficient to organize tools.
 Ingeniously designed construction with stability and large bearing capacity.
 Multiple optional colors.



TK-608-BK

TK-608-PK

Tools Case

TK-608

Made from aluminum and plastic composite, ensuring both light weight and good strength.
 Collapsible dragging handles makes it easy to move around.
 Versatile usage realized by flexible snapping clips.
 Extending trays built-in for easy sorting of your small stuffs.
 Two standard colors available, pink and black.

size:340mm×240mm×418mm



TK-606-BK

TK-606-PK

Small Tools Case

TK-606

Measures 295mmLx190mmWx200mmH.
 Lightweight and practical solution for handy tools organizing.
 Lifting collection trays included inside.
 Black and Pink colors available.



Rolling Tools Bag doubling as Stool

TK-610

Very strong, very sturdy.
 Aluminum brackets cooperating with smooth castors ensure comfortable mobile using.
 Practical top pad gives a great working seat to groomers.
 Good for dog shows, grooming contests and all mobile groomers.
 Removable cover skirt for easy cleaning purpose.
 Telescopic dragging handle gives great mobility.
 Two zipper bags and one plastic tray included.



Medium Tools Case

TK-607

Medium size tools case.
 Sitting on 4 smooth rolling castors.
 Size: 360mmLx235mmWx510mmH.
 Collection trays and organizing pockets inside.
 Black and pink available.



TK-607-BK

TK-607-PK



Aluminum Tool Case

TK-603



Flat tool case

TK-609



SO-605
Pink

SO-608
Gray

SO-604
Purple



SO-504
Purple

SO-505
Pink

SO-503
Green

SO-501
Blue



SO-402
Orange

SO-404
Purple


SO-405
Pink

SO-403
Green

SO-407
Black

SO-401
Blue

SO-406
Red

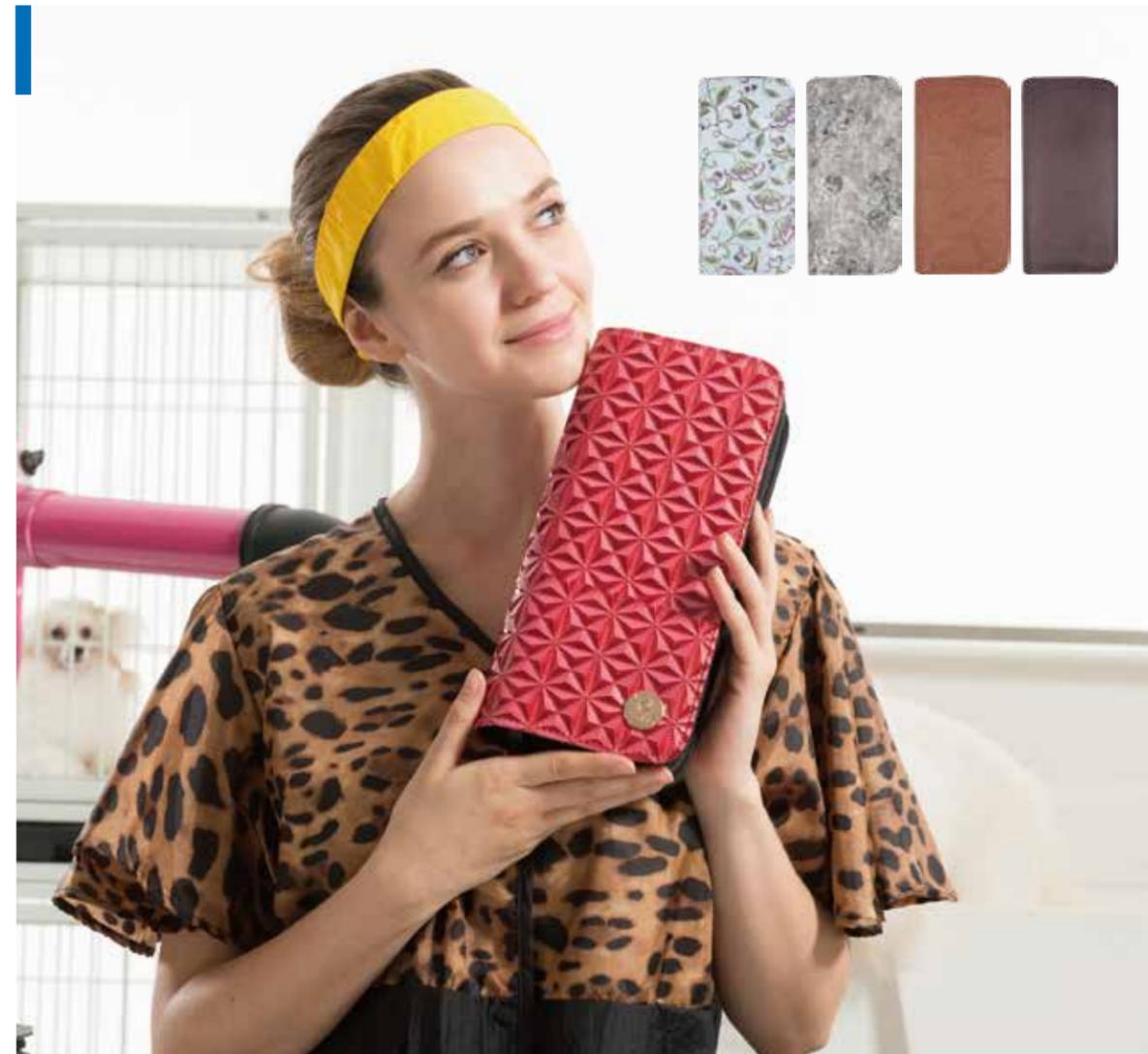

Professional Scissor Cases

AA	Pattern	Demonstrate	BB	Color
01	Diamond		01	Black
			02	Purple
			03	Pink
			04	Red
			05	White
			06	Yellow
02	Chinese Phoenix		01	White Background, Gold Pattern
			02	White Background, Silver Pattern
			03	Black Background, Bronze Pattern
			04	White Background, Purple Pattern
05	Luxury Alligator		01	Pink
			02	Grass Green
			03	Yellow
06	Python Skin		01	Sapphire Blue
			02	Earth Yellow
			03	Rose Red
			04	Purple
			05	Black
07	Steady Cow Fur		01	Black
			02	Blue
			03	Brown
			04	Deep Blue
08	Leopard		01	Leopard
			02	Black Leopard
			03	Purple Leopard

Blade & Shear Organizer



BO-201
16-cell blade organizer



BO-202
8-cell blade organizer



SO-301
12-cell shear organizer



SO-302
6-cell shear organizer



SO-303
5-cell shear organizer


Basic Information				Packaging Dimension (cm)			Details	Description
Category	Material	Size	Item No.	Length	Width	Height	Holes	
Blade Organizer	Acrylic	Large	BO-201	45	6	8	16	This professional blade organizer is made of clear Acrylic. Break-resistant holder makes it easy to store up to 16 clipper blades (16 holes in 2 rows). Large size holds up to 16 blades; small size holds up to 8 blades. They protect and keep your clipper blades safe, organized, and close at hand. It can be mounted to walls with ease. It is a great way to keep clipper blades off surfaces and out of drawers and where they could be damaged.
		Small	BO-202	23	6	8	8	
Shear Organizer	Acrylic	Large	SO-301	42	13	6.5	12	This convenient wall-mounted shear organizer is made of high-quality, break-resistant Acrylic. It keeps shears well-organized and well-protected. This top performance organizer is easy to install. There sizes available, large holds 12 shears, small holds 6 shears and stand holds 5 shears.
		Small	SO-302	24.5	13	6.5	6	
		Small	SO-303	10	10	18	5	



Blade Storage Box

BB-101/102


Smartly designed for blades storage, with sponge on the bottom to contain oil and prevent blades from oxidizing.
Made from high quality acrylic material, the supplier for UK pavilion in Shanghai World Expo.
Two sizes available, small holds 8 blades and large holds 16 blades.

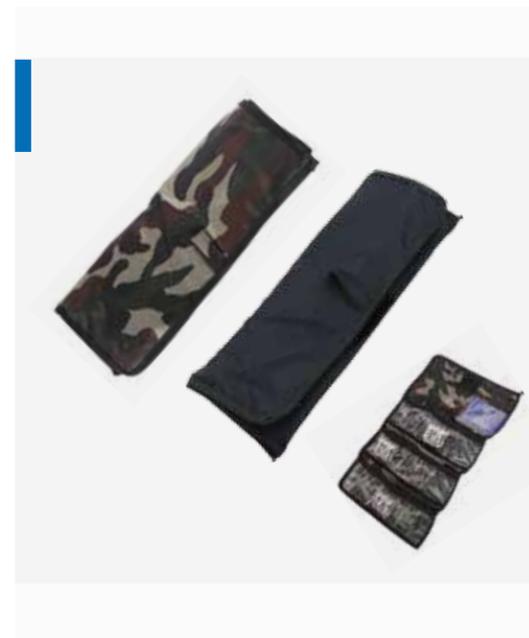

Blade Vest

BV-505

Portable "blade apartment", organizes and protects your blades, keep them close at hand and travel ready.

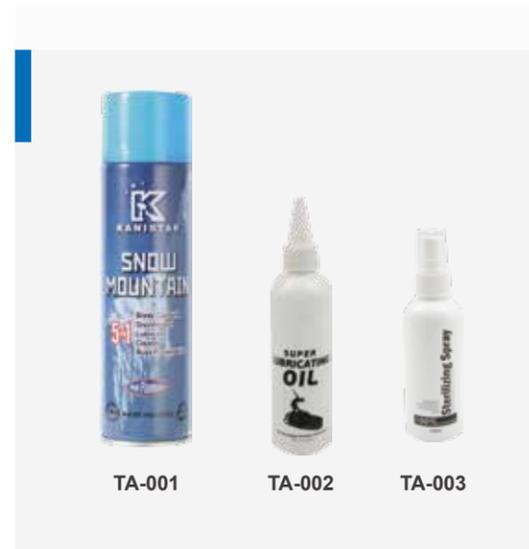

Cooler Spray TA-001

Blade coolant
Disinfectant
Lubricant
Cleaner
Rust Preventative
Non-flammable and CFC free.

Lubrication Oil TA-002

Specifically for blade lubrication. Does not contain impurities. Help maintain and extend the life of blades. Can be used for soaking blades.

Sterilizing Spray TA-003

Disinfection and sterilization spray. Applicable to shears, blades, combs, etc. Help avoid cross contamination



Kitten Chest of Drawers

DF-801

Modelling vogue is lovely, the work is exquisite and exquisite, the space design is reasonable.

The finishing touch of fashion shop.




Series Motor Nail Grinder

NG-003

Strong motor, dynamic, safe and reliable, save effort.


There is a replacement file for the purchase.

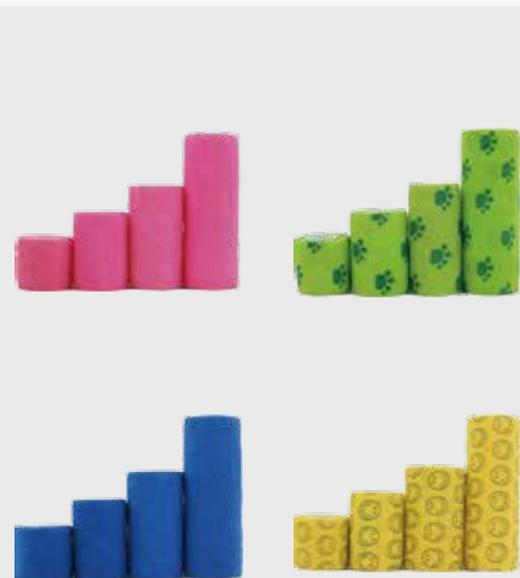

Elastic Bandage

Colorful, lightweight, clings to itself, not to the coat or skin, non irritating.

Series	Color	Width	Code	Description
EB 02 06	2	2→2"	0101	Pure Color/Neon Pink
		3→3"	0102	Pure Color/Purple
		4→4"	0103	Pure Color/Dark Blue
		6→6"	0104	Pure Color/Caramel
			0105	Pure Color/Neon Green
			0206	Printed Smile Face/Yellow
			0306	Printed Paw/Yellow
			0307	Printed Paw/Grass Green
			0400	Printed Camo/Camo
			0508	Printed Love /Pink
			0609	Printed Horse /White
			0708	Printed Cat /Pink

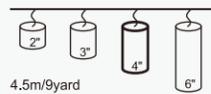

4.5m/9yard

AC-DS-100 DISPOSABLE MASK



DT-65 ABSORPTION TOWEL



AC-DC-100 DISPOSABLE CAP



KIDO BUNGEE



Mobile Van Solution Provider

Mobile grooming service is getting more and more popular all over the world. Without annoying overheads, based on inexpensive investment, with much more flexible working time, mobile service has shown its beauty to many groomers.

For the mobile medical care and emergency service, a fully equipped van with all the necessary equipment seem to be more important, which could affect the animals' living or death.

In different area of the world, you will meet varied van chassis and quite different requirement from clients. Based on over 8 years of customization service experiences to over 10 different typed vans, we are able to provide very detailed solutions to almost all the potential requests.



Our engineering team knows how to play the game with compact space, acknowledged of how to arrange plumbing, water supply, illumination and air ventilation in the small room, also familiar with all the tricks of how take good advantage of a wall, like making it as a storage rack and also a TV display wall by using a sliding TV mount, etc.



After dealing with the very big RVs and very tiny small transportation vans, we found a balancing point that how to make both "elephant" and "ant" happy. Being requested from all the aspects in a limited space, we start to know

how to make all the components and every corner "functional". Have been "tortured" by many extreme challenges, our designers and workers are very calmed and easy with all the potential tasks.

We have made special grooming tables with extremely big legroom for groomers to be able to sit down along a very narrow aisle; worked out a solution for one bathing tub could contain SPA machine, force blaster, powerful heater, boosting pump and water recirculation in a 0.25cbm compact box; maximized the use of a 200L water tank by an efficient water purification and recirculation system, etc. Our team organized all the medical monitors, oxygen generators, endoscopes, ultrasound machine, and other devices into one integrated hanging pendant; put all needed lights, cameras, X-ray machines into one 50cm wide by 100cm long roof; invented a brand new convertible operation table for the van extension part, etc.



After all the solution tested by rattling country in South Africa, challenged by the Australian desert, inspected by the US 110V voltage, now we are ready to give you, the integrated van equipment R&D service! Take your van, we are the one!



Mobile Pet-Grooming Van

Mobile grooming service is getting more and more popular all over the world. Without annoying overheads, based on inexpensive investment, with much more flexible working time, mobile service has shown its beauty to many groomers. Based on over 8 years of customization service experiences to over 10 different typed vans, we are able to provide very detailed solutions to almost all the potential requests. Our engineering team knows how to play the game with compact space, acknowledged of how to arrange plumbing, water supply, illumination and air ventilation in the small room, also familiar with all the tricks of how to take good advantage of a wall, like making it as a storage rack and also a TV display wall by using a sliding TV mount, etc. We have made special grooming tables with extremely big legroom for groomers to be able to sit down along a very narrow aisle; worked out a solution for one bathing tub could contain SPA machine, force blaster, powerful heater, boosting pump and water recirculation in a 0.25cbm compact box; maximized the use of a 200L water tank by an efficient water purification and recirculation system, etc.





Mobile Pet Medical Van



For the mobile medical care and emergency service, a fully equipped van with all the necessary equipment seem to be more important, which could affect the animals' living or death.

Based on over 8 years of customization service experiences to over 10 different typed vans, we are able to provide very detailed solutions to almost all the potential requests. After dealing with the very big RVs and very tiny small transportation vans, we found a balancing point that how to make both "elephant" and "ant" happy.

Being requested from all the aspects in a limited space, we start to know how to make all the components and every corner "functional".

Have been "tortured" by many extreme challenges, our designers and workers are very calmed and easy with all the potential tasks.





Our team organized all the medical monitors, oxygen generators, endoscopes, ultra-sound machine, and other devices into one integrated hanging pendant; put all needed lights, cameras, X-ray machines into one 50cm wide by 100cm long roof; invented a brand new convertible operation table for the van extension part, etc.

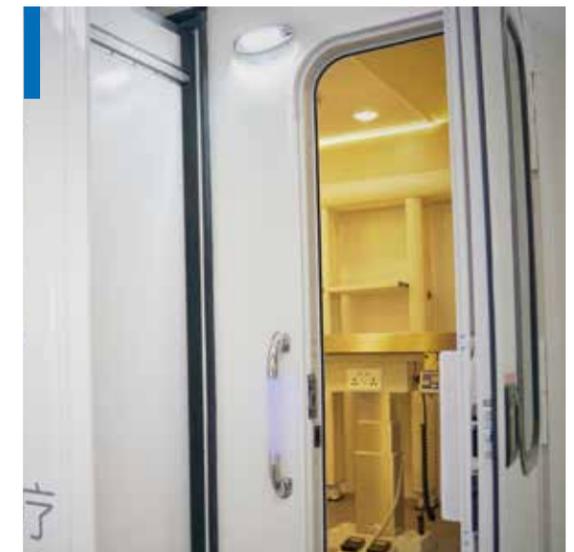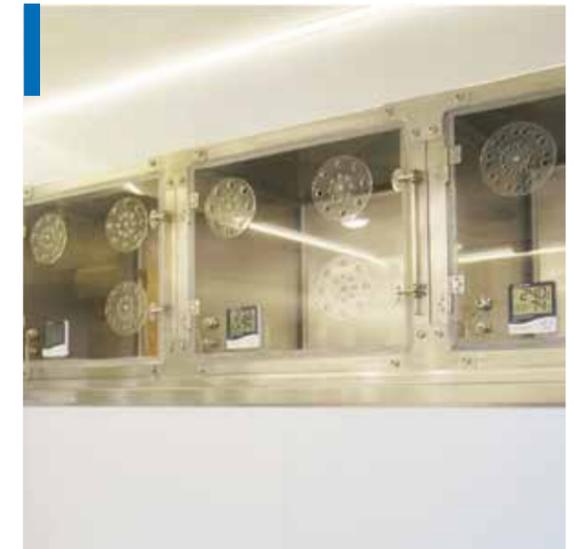



Pet Medical and Grooming Van

A well-designed multi-purpose pet van allows for flexible service choices and more potential customers. Our engineers make extra use of compact van room. Multi-purpose lifting table can easily transform between grooming table and operation table by tabletop change. Intelligent electric and water supply makes for bath and electrical equipments. Drying Cabinet saves your labor and time. Just put one pet in and start drying, and then your can turn to another. Pets have no change to scratch or escape. Unbelievable cleaning effect brought by our Professional Micro Bubble SPA Tub! Storage shelves makes van work part-time as a mobile shop. Seeing is believing!

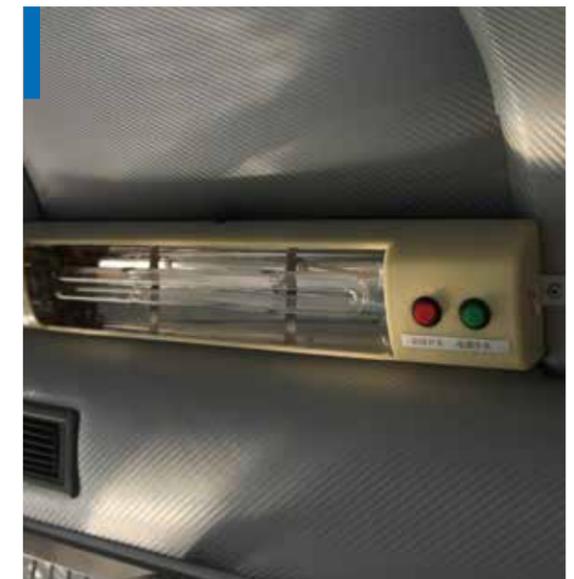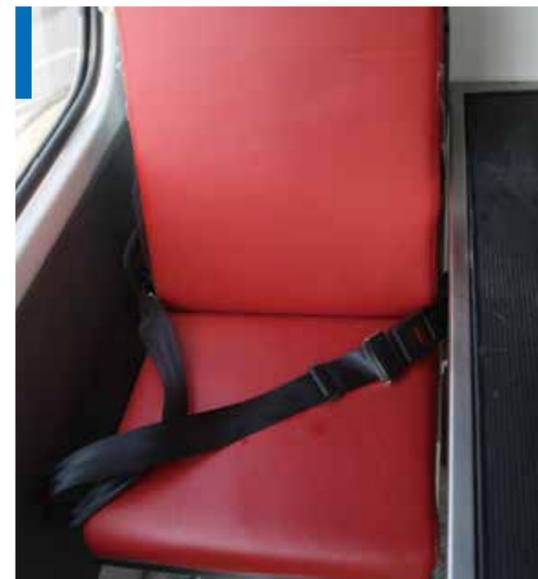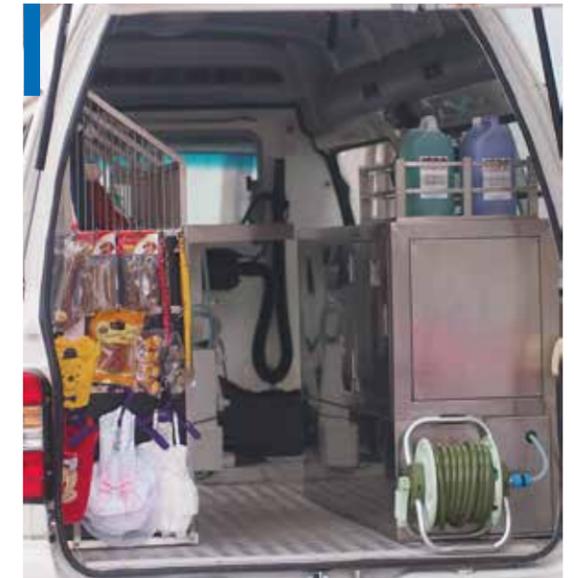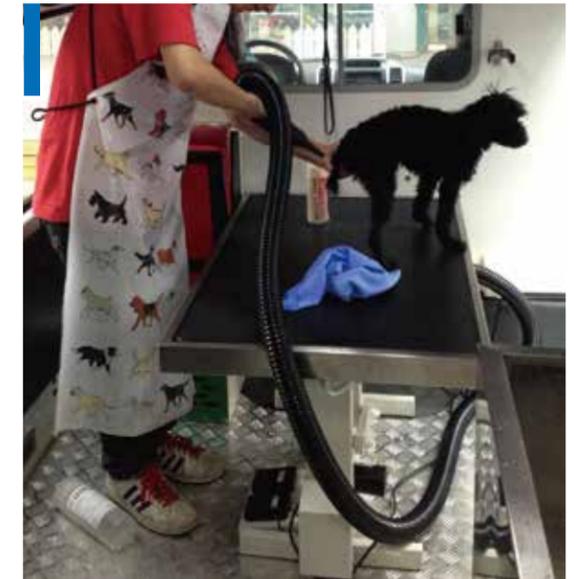




Animal Positioning bracket

Veterinary MRI system

Magnetic Resonance Imaging (MRI) is a noninvasive imaging technology that is known for its superior soft tissue contrast and ability to provide morphological, microstructural, functional, and metabolic information. MRI does not involve X-rays and the use of ionizing radiation, which distinguishes it from CT or CAT scans. MRI is the investigative tool of choice for neurological cancers, as it has better resolution than CT and offers better visualization of the posterior fossa. The contrast provided between grey and white matter makes MRI the best choice for many conditions of the central nervous system, including demyelinating diseases, dementia, cerebrovascular disease, infectious diseases, and epilepsy. MRI's Cardiac applications include assessment of myocardial ischemia and viability, cardiomyopathies, myocarditis, iron overload, vascular diseases, and congenital heart disease. MRI's applications in the musculoskeletal system include spinal imaging, assessment of joint disease, and soft tissue tumors. Hepatobiliary MRI is used to detect and characterize lesions of the liver, pancreas, and bile ducts. MRI is also used to evaluate the arteries of the neck and brain, the thoracic and abdominal aorta, the renal arteries, and the legs. Since MRI does not use any ionizing radiation, its use is generally favored in preference to CT when either modality could yield the same information. In certain cases, MRI is not preferred as it may be more expensive, time-consuming, and claustrophobia-exacerbating. Our veterinary MRI system is your most cost-effective choice.

Technical characteristics:

1. Powerful permanent magnet material N52, open magnet design to reduce the weight. Field strength: 0.3T, 0.35T, 0.4T (customizable)
2. High performance self-shielding water-cooled gradient coils, eddy current suppression and **active shielding technology** make better imaging performance.
3. Well-designed RF shielding system makes installation easy and fast.
4. **Customized RF coils** for all kinds of animals, suitable for small animal like mouse to large dogs.

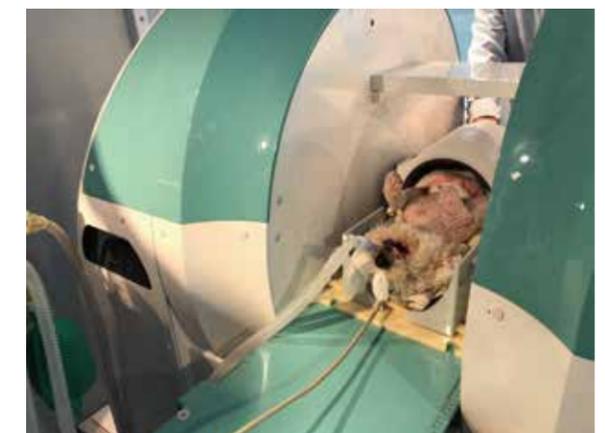



Large-size RF coil




Middle-size RF coil




Small-size RF coil



Large-size U-shaped coil

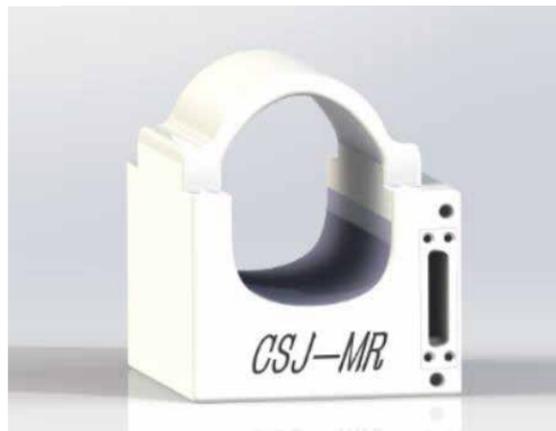

Small-size U-shaped coil

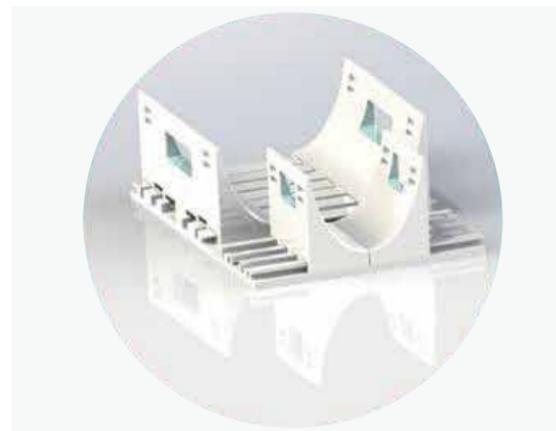

Pet restraint

Technical parameters:

Magnet

Field strength: 0.3T,0.35T,0.4T (customizable)

Homogeneity: <2PPM @20cm DSV

Gradient Type: High performance self-shielding water-cooled gradient coils

Gradient Amplitude: 18-25mT/m

Slew rate: 45-75mT/m/ms

Pulse sequences

Spin echo (SE)

Turbo spin echo (TSE)

Standard IR (IR)

Short inversion time IR (STIR)

Fluid Attenuation IR (FLAIR)

2D Gradient echo (GRE2D)

3D Gradient echo (GRE3D)

3D turbo spin echo (TSE3D)

Imaging capabilities


field of view: ≤100mm

Acquire Matrix: 64×64-1024×1024


Space resolution: ≤1mm



Canine Knee

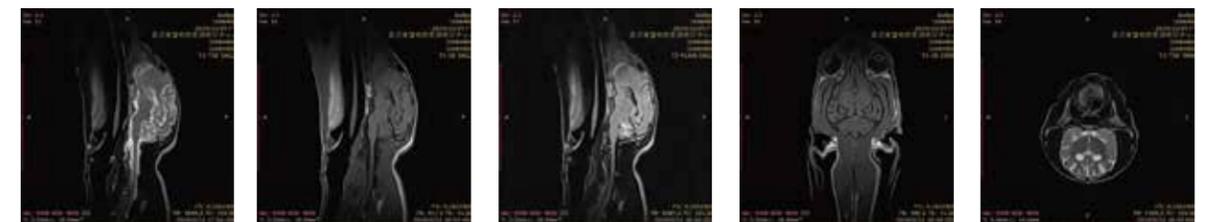

Canine Head

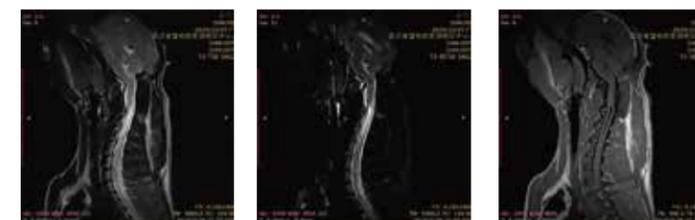


Canine Spine

5. Special designed magnetic resonance compatible pet anesthesia monitoring system, much safer.
6. Adjustable height of animal table, makes it possible to position different size of animal to the central of the MRI system.
7. Special designed pet placement tools, makes it easy to positioning the animals in seconds.
8. Easy using software. Abundant 2D and 3D imaging sequences and technologies
9. One phase AC required and low power consumption.




Preclinical Molecular Imaging System SmartView-3000B Small Animal PET/SPECT/CT


SmartView dual-layer LYSO detector achieves good DOI resolution.
 The smartView electronic system uses proprietary ASIC chips to read in the large number of channels in high-density SiPM arrays.
 SmartView SPECT uses a multi-pinhole collimator design with pinholes embedded on a cylinder. Rats and mice can be imaged respectively by axial movement of the cylinder.
 List-mode data acquisition preserves maximum original information.
 Spatial resolution recovery is enabled through accurate point spread function modeling.
 DOI response model of dual-layer crystal enhances resolution uniformity in the full FOV.
 Attenuation correction using high-resolution CT images ensures the quantitative accuracy of PET and SPECT images.
 Dynamic PET imaging assists pharmacokinetic studies.
 SmartView support units for animal study: Anesthesia system, video monitoring system, body temperature maintenance system.

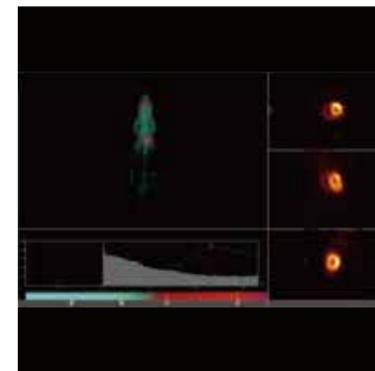


PET Mouse Myocardial Imaging

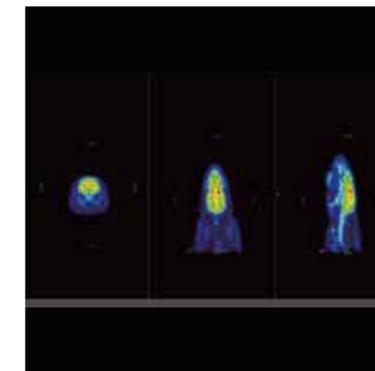


PET Rat Brain Imaging

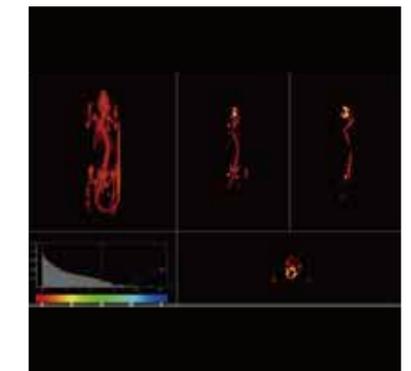


PET Rat Bone Imaging

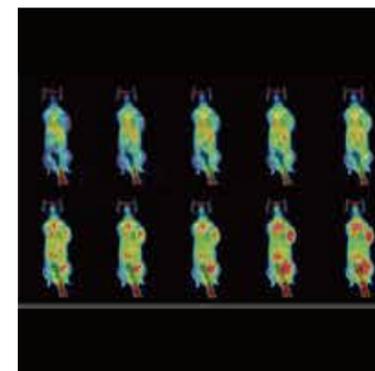


PET Mouse Tumor Dynamic Imaging

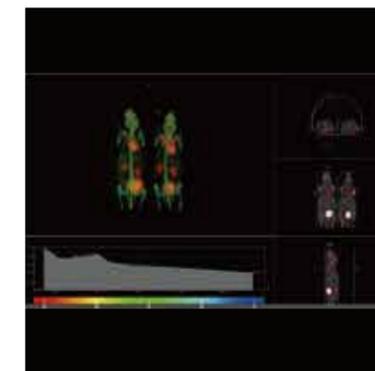


PET Double-mice FDG Imaging

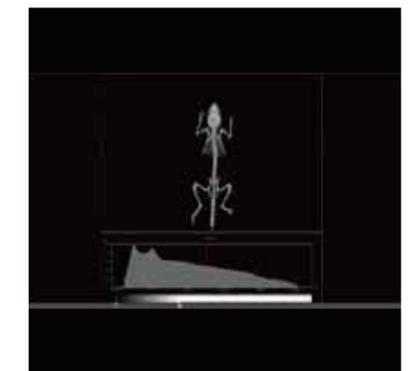


Mouse CT Bone Imaging

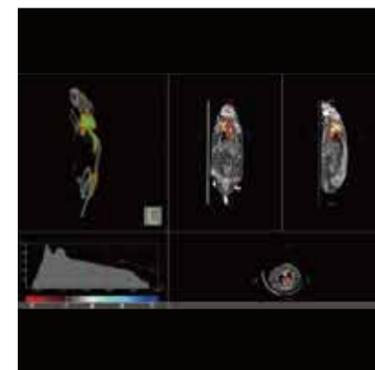


PET/CT Mouse Brown Fat Imaging

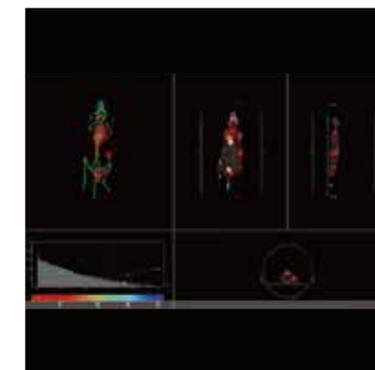


Mouse Tumor PRT/CT Imaging

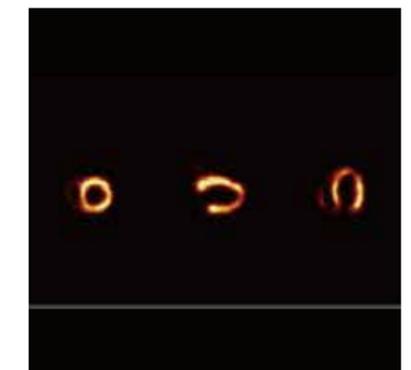


SPECT Rat Myocardial Imaging

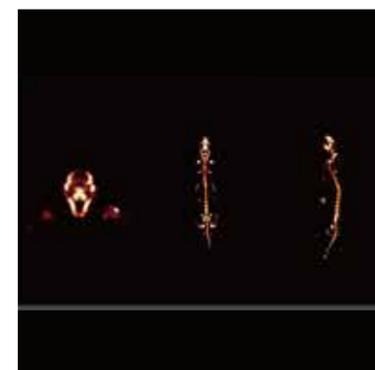


SPECT Rat Bone Imaging

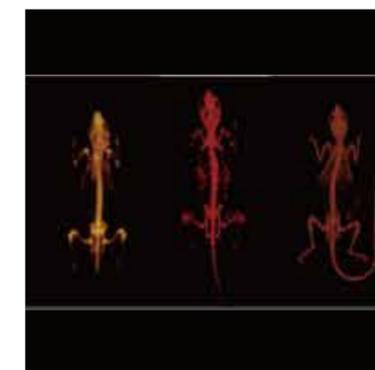

Mouse PET/SPECT/CT Bone Imaging

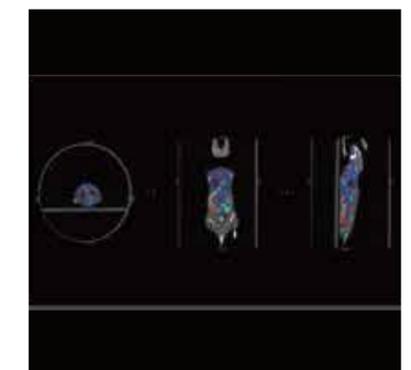


SPECT Rat Bone Imaging

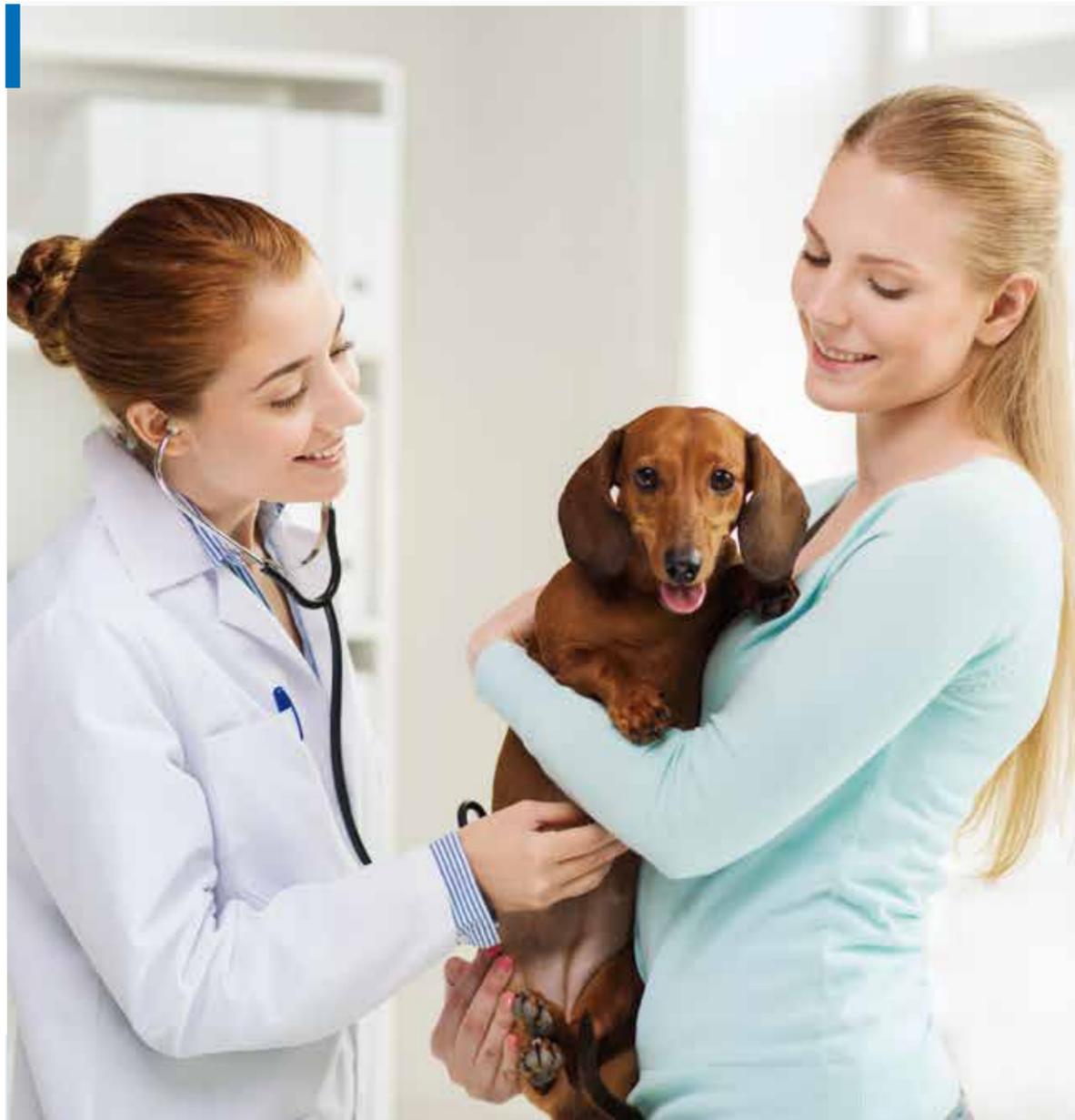

UltraBright 7 Vet

UltraBright 7 Vet Ultrasound is a wonderful all-in-one system, which covers all the vet specialties for daily use. With the most advanced patented technologies, UltraSmart 7 Vet brings excellent image quality, unbeatable color sensitivities and user smart operation. UltraBright 7 Vet empowers precise and confident diagnosis like never experienced before.



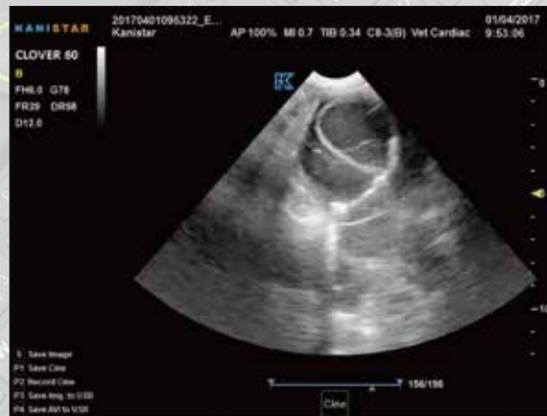

Exceptional Imaging Quality

Uncompromising whole field best in class image quality

General
Abdominal
Reproductive
Musculoskeletal

Cardiac
Urology
Small Parts

Scanning Mode

- B mode
- M mode
- Tissue Harmonic Imaging
- Color Doppler Imaging
- Directional Power Doppler Imaging
- Pulse Wave Doppler (PW)
- Continuous Wave Doppler (CW)
- Spatial compounding imaging (SCI)
- Time-spatial speckle reduction imaging

Transducers



Convex



Micro-convex



Linear



Phased



High-frequency Phased





UltraBright 50Vet

UltraEasy 50Vet Full Digital Ultrasound Diagnostic System

Specifications: Professional Vet use software

3 Probe connectors

15 Inch Medical LCD Screen

8 segments TGC adjustment

Standard PC keyboard

Original Windows Operation System

1000 frames Cine loop

1000 frames image storage

320GB hard disc

DVD Read and write optional

Optional Digital printers available

Optional video printers available

Free hand 3D function optional

Internal Lithium battery working continuously up to 3 hours

5 times real time local zoom

Multiple connection ports including Video, XGA, 4 USB 2.0, Dicom 3.0, RS-232

95 types of body mark.

Probes available:

Convex

Linear


R20 Micro-convex

R15 Micro-convex

Animal rectal linear.

Full Digital Technology

Full digital engine ensures high definition images.

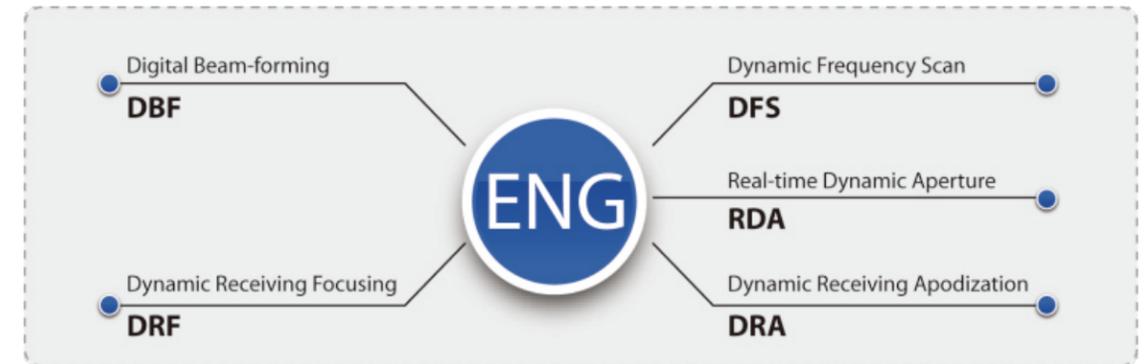


Compact and Portable

Outstanding Image Quality

Large Capacity for Data Storage

Tissue Harmonic Imaging

High-caliber harmonic signal processing technology

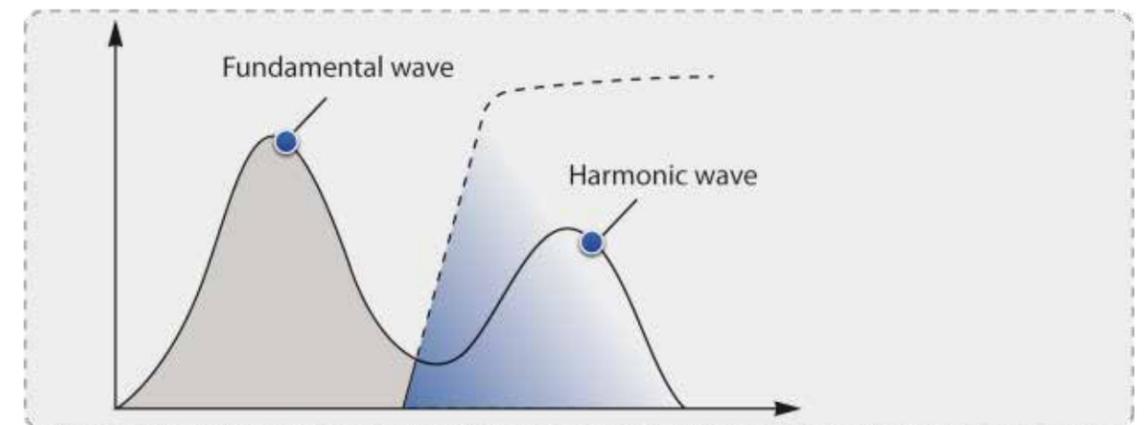

Technical Specification

15 inch LCD Display

B, B+B, M, B+M, 4B

320G Hard Disk

B/W, Colorized Images

Tissue Harmonic Imaging

Image Enhancement

Noise Rejection

Zoom In Real and Frozen

OB Measurement for Vet

Dual USB Ports

VGA Image Output



Transducers:




3.5MHz
Convex Probe



5.0MHz
Micro-convex Probe



7.5MHz
Linear Probe

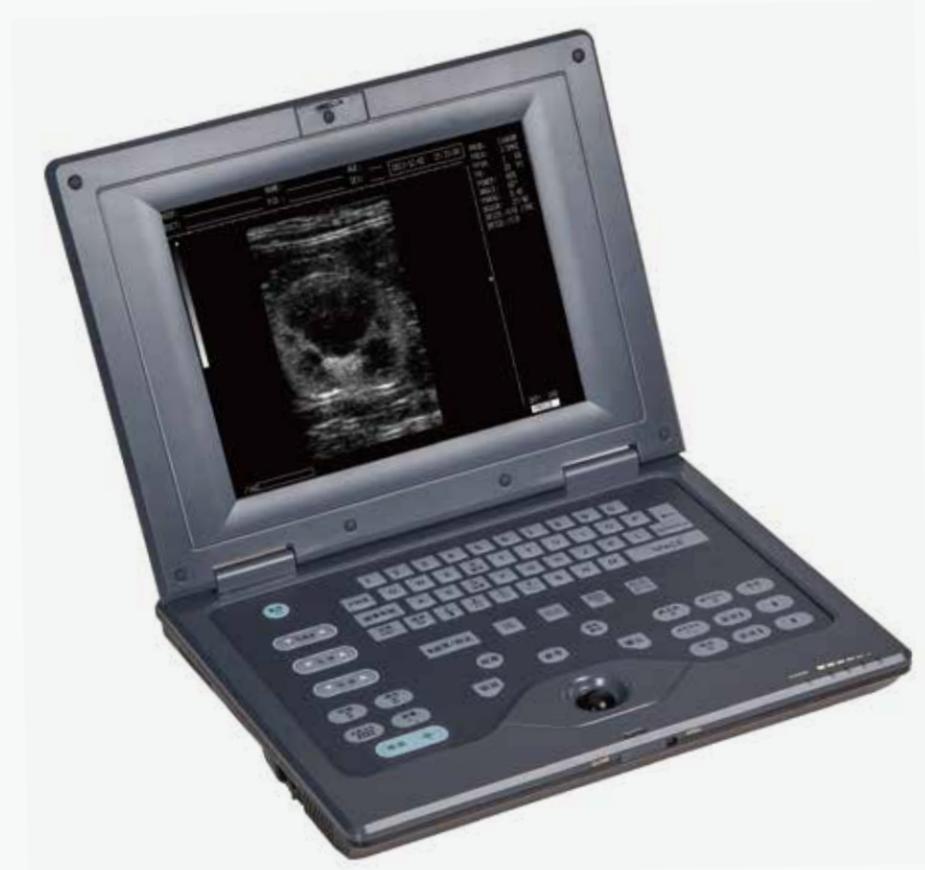

UltraEasy 20Vet

UltraEasy 20Vet Full Digital Ultrasound Diagnostic System Specifications:

Small and compact for easy carrying

Notebook design

10.4 Inch Medical LCD Screen

Internal Lithium battery working continuously upto 3 hours

Optional Car charger

400 frames Cine loop

64 frames image storage

Optional video printer available

Probe automatic identification.

Probes available:

Convex

Linear

R20 Micro-convex

Animal rectal linear.



Ultralight 9 Vet

A brand new tablet type ultrasound adapting all the probes with great mobility!
 Extremely lightweight case, smart design and built-in customized workflow;
 More efficient and convenient experience for daily operation, especially musculoskeletal, pain management, physical therapy and vascular.

Best in Detail Tissue Imaging: 64-beam processing and achieves up to 5000 frames per second;
 Able to provide not only the whole field best in detail image quality but also highlight the display of bones, muscles and nerves;

13.3-inch touch screen, horizontal or vertical display, height adjustable, support surface liquid disinfection, horizontal or vertical transformable design;

A wide range of transducers (optional) covers full range of general applications, from large dogs to small cats.

An innovative middle line design along with project probe direction indicator, provide users best experience and accuracy;

3 Probe buttons: Corresponding function name display, Remote control;

With a Built-in Storage Tray: A dedicated probe holder for easy disinfection without any need of assistance from others.



13.3-inch touch screen, horizontal or vertical display, height adjustable




A wide range of transducers




3 Probe Buttons



Hockey Stick Probe--small joint and superficial

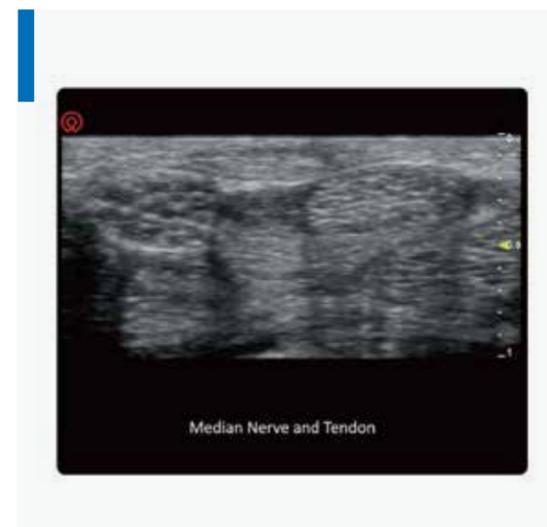

Best in Detail Tissue Imaging

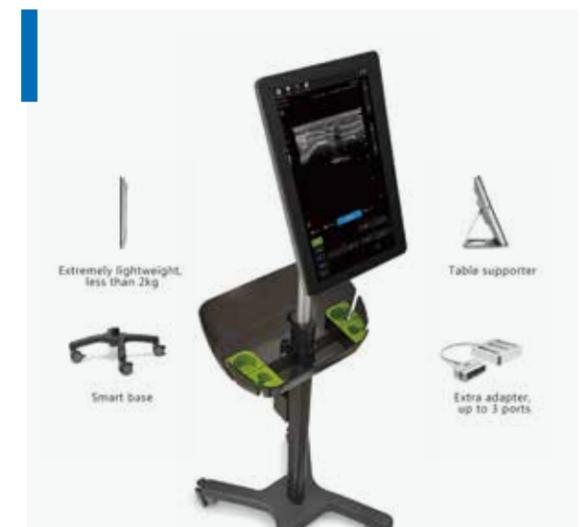

Flexible Style, More Possibility



Endoscope

| ES-801

Endoscope made and assembled in France.
 Long Life LED light source is made in USA.
 Japanese made processor and motherboard.
 First class accuracy and sharpness.
 Captures image immediately without delay and stays in focus during movement and repositioning.
 True 720P resolution, dedicated medical monitor.
 Four diameters of endoscopy available, 3mm, 6mm, 9mm and 12mm.



Medical Class Monitor



720P High Resolution Image Processor



US Made LED Long Lasting Light Source



Lens Made and Assembled in France



Ergonomically, Easy to Operate Handle



Sturdy and Functional Trolley

Endoscope Instruments



TX-EDS-104-037
Reusable maryland dissector
φ3×280mm



TX-EDS-104-033A
Reusable curved scissor
φ3×280mm



TX-EDS-101-047-250
Reusable babcock forceps
φ5×250mm



TX-EDS-301-044C,
Reusable bipolar with cable,
curved jaw



TX-F4311
Trocar and Trocar Sleeve



TX-YM270
Telescope 0°, angled 45°,
instrument, channel
φ6.2×370




TX-101-151
Suction and irrigation,
with trumpet valve
φ 5/10×330mm

TX-EDS-101-034-450
Reusable hook scissor
φ5×450mm



TX-EDS-101-046A-450
Reusable claw forceps(2+3 teeth)
φ 5×450mm

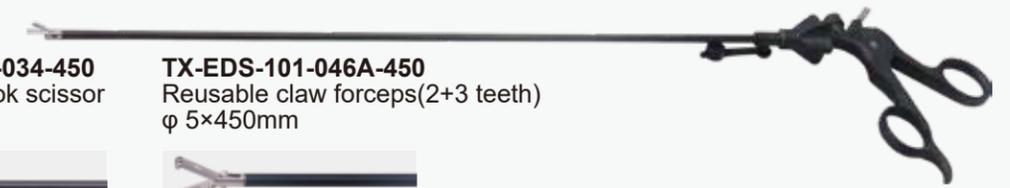

TX-EDS-101-033A
Reusable curved scissor(Dual action),
φ5×330mm



TX-EDS-101-037C
Reusable blunt dissector
φ5×330mm, length of jaws:18mm



TX-EDS-101-049A
Reusable long fenestrated grasping forceps
φ 5×330mm



TX-EDS-101-033B-250
Reusable micro curved scissor(Dual action)
φ5×250mm



TX-EDS-053-250
Reusable biopsy forceps
φ5×250mm



TX-EDS-101-113B
Clip applicator for polymer ligation clips, ML
φ 5×330mm



These are only some typical examples of our big instrument range, for more instruments, please contact local distributor and refer to KaniStar official website.



TX-BM-01


Biological Microscope

TX-BM series biological microscope is best choice for college teaching and hospital clinical examination, and laboratory;
 Infinite Optical System, plano-achromatic objectives, wide-field eyepieces, anti-fungal cameras;
 High resolution and perfect definition from infinity correct optical system;
 100X water objective instead of 100X oil objective, more convenient for using;
 Drop-style and elegant looking;
 Perfect Ergonomics Design of low position operation, convenient for long-term using;

Technical Specifications

Magnification: 40X-1000X
 Eyepiece: 10X
 Objective: 4X, 10X, 40X, 100X
 Mechanical Tube Length: 160mm
 Head : Inclined 30° Seidentopf Binocular (Trinocular) , Diopter adjustable range $\pm 5^\circ$, anti-fungal .
 Interpupillary Distance: 50-75mm
 Nosepiece: Quadplex Nosepiece
 Mechanical stage: Size 125mm \times 125mm , Y-X travel 35mm \times 75mm
 Focus Knob: Coaxial Coarse and Fine Focus Knobs, Fine division 2 μ m.
 Condenser: Abbe, N.A. 1.25, Adjustable aperture
 LED Light Source:
 Input voltage AC100V \sim 240V , 50/60Hz Output voltage DC0.4 \sim 3.4V ;
 LED 3.4V/3W Brightness adjustable Life time: above 10000 hours.

Optional Cameras Installation

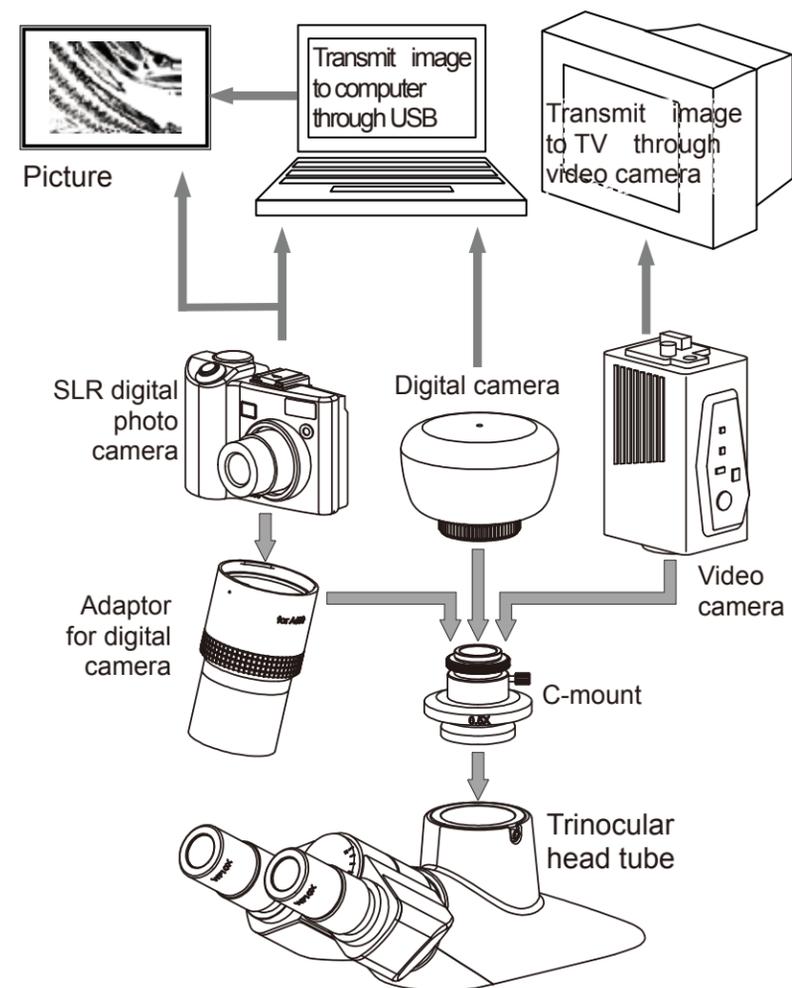



Biological Microscope

TX-BM-01

Binocular without accessory camera



Biological Microscope

TX-BM-02

Trinocular with Digital Camera
The digital camera can provide microscopy images, which have 2952×1944 active pixels in size of 2.2um×2.2um, and output images to the software in a computer via USB2.0 cable.



Biological Microscope

TX-BM-03


Trinocular with Video Camera
The video camera can provide microscopy images of the resolution of 1920*1080 and display them on TV or monitor via HDMI cable.



Biological Microscope

TX-BM-04

Trinocular with built-in Digital Camera and Monitor
The microscopy images can be displayed on the monitor.

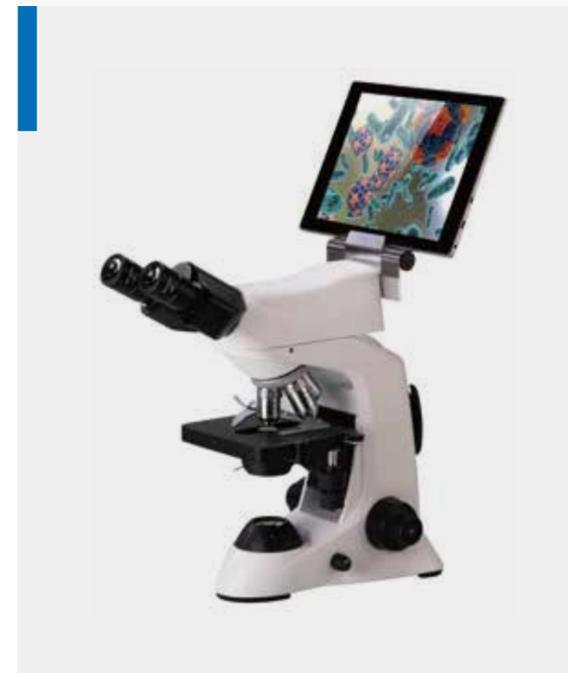

Biological Microscope

TX-BM-05

Trinocular with Tablet Computer
The microscopy images can be displayed on the tablet computer.



KANISTAR Veterinary Anesthesia Machine

| ANS-1701

Used for pet anesthesia before and during surgery, especially thoracotomy, lung tumor resection, cardiac surgery, thoracic and esophageal surgery, diaphragmatic hernia repair surgery and so on.

Safe, reliable, and controllable. Isoflurane's pulmonary (lung) metabolism is fast, low-load with few side-effects. Especially suitable for pets who are old, weak, or have metabolic disorders in liver or kidney.

Minor allergic reaction. (During traditional injection anesthesia, once allergy happens, it's too late to stop anesthesia and rescue. During gas anesthesia, veterinarians can immediately cut off the anesthesia gas source).

Good anesthesia effect and analgesic effect, particularly suitable for long time surgery. Tracheal intubation is convenient for first-aid oxygen supply.



Pressure gauge
Left side is airway pressure, right is input oxygen pressure.



Anesthetic vaporizer



CO₂ absorber



Flowmeter




Reservoir bag



Durable casters

Advantages of the Gas Anesthesia

- Safe;
- Good anesthetic and analgesic effect ;
- Easy to control; adjustable anesthetic depth at any time during surgical anesthesia state;
- Little influence on heart and lung;
- No influence on liver or kidney for isoflurane is exhaled through lung;
- Keep airway patency and good for oxygen supply;
- No allergy;
- Fast recovery;
- Conducive to rescue in case of emergency.

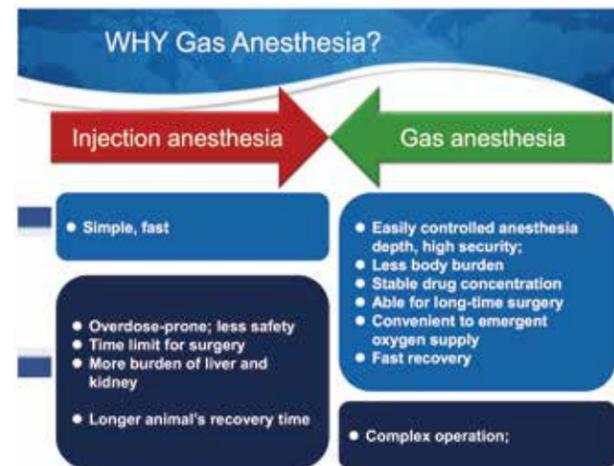

Overview of Gas Anesthesia System



Veterinary Anesthesia System

Friendly Powerful Reliable

The Anesthesia machine makes a good performance in Intensive Care Units (ICU), Operation room, Anesthesiology Department and other departments. Good performance based on High Security, High Accuracy, High Stability and accurate monitoring.

Easy to use, easy to move with 4 wheels;

Patient Centered Ventilation: Precision in an anesthesia ventilator, from conventional ventilation to advanced modes, including 5 modes(IPPV; A/C; PCV; SIMV; SIGH)/ 3 modes(IPPV; A/C; SIMV.) Flexible configurations to suit your needs;

International standard and advanced technology suitable for wide range use;

Compact interface and big screen give you better operating experience;

10.4" TFT LCD screen displays the Ventilation parameters, Alarming information and Waveform;

High precision Mechanical/ Electronic flow meter, instantly know the fresh gas flow to your patient;

Integration breathing circuit design, ensure easy operating and keep tidy;

Multiple working modes such as volume control and pressure limit, adapt to wide range patient;

Vaporizer with temperature, flow compensation and self-lock function, keep safety anytime;

Multiple parameters monitoring interface, make every parameter clear, let users know the patient conditions in all aspects;

Real time pressure-time, flow-time loop graphics and high precision O₂/ETCO₂&O₂ concentration detection function included;

Vital sign monitor and Anesthetic gas monitor are optional;

Wide use range including animals within 200KG.

Safety

Three level alarming system, visual and sound alarm information;

With lots of alarming, reminding and protection functions;

Advanced power management control technology;

With built-in backup power source, when outside power source goes off, back-up power source starts to work;

Low O₂ pressure alarm.

Parameters Comparison

Model	ANS-1801A	ANS-1801B	ANS-1701
Picture			
Display Screen	10.4" Touch Screen	10.4" Touch Screen	N/A
Gas source	O ₂ , N ₂ O	O ₂ , N ₂ O	O ₂
Ventilation Modes	IPPV, A/C, PCV, SIMV,	IPPV, A/C, SIMV	MANUAL
Tidal volume (Vt)	0/20~1500ml	0/20~1500ml	N/A
I/E	4:1 ~ 1:8	2:1 ~ 1:6	N/A
Pressure Limit	20 ~ 100 cmH ₂ O	20 ~ 100 cmH ₂ O	20 ~ 100 cmH ₂ O
Vt, PEEP Control	Electric Control	Electric Control	N/A
Oscillographs	P-T, F-T, P-V loop, V-T, ETCO ₂ -T	P-T, F-T	N/A
Flow meter	Electronic	Mechanical	Mechanical
Alarms	14 Kinds of alarms	14 Kinds of alarms	N/A



KaniKare ICU

| UC-1801

2018 Newly Launched ICU.

Improved Humidification System.

The new water reservoir is placed in the middle of the right sided chamber, easy to refill water and check running status.

The new heavy duty ultrasound humidifier ensures an ultra high humidification performance, now the humidity could be increased to 85RH (at 28 Celsius Degrees), which is required by many exotic species.

The visible system also provide easier maintenance access.

Improved Nebulization System

Double air compressor installed in the nebulization system, to ensure a more secured working stability. Even one line get kinked during operation, the other compressor will still provide enough air stream to complete the needed operation.

Accurate Temperature Control—Eliminates hot spots and provides overall consistent temperature throughout.

Accurate Humidity Control— Provides a healthy and balanced environment.

Negative-ion Generation—Augments the Medical Effects of the ICU.





Sterilization Function—Greatly diminishes the possibility of cross-contamination.

Nebulization/Medical Atomization Treatment Function—Making it more convenient to administer full strength medicine.

Carbon Dioxide Concentration Monitoring System—A Defense Line to safe guard patients' lives.

ICU Illumination Function—Creating a comfortable therapeutic environment.

Setting of Security Mechanism—No worries in use.

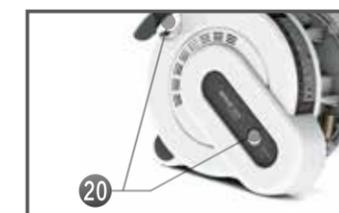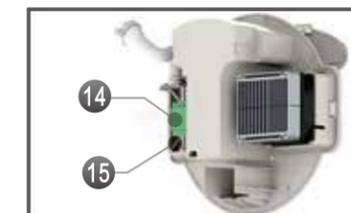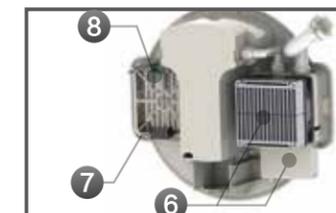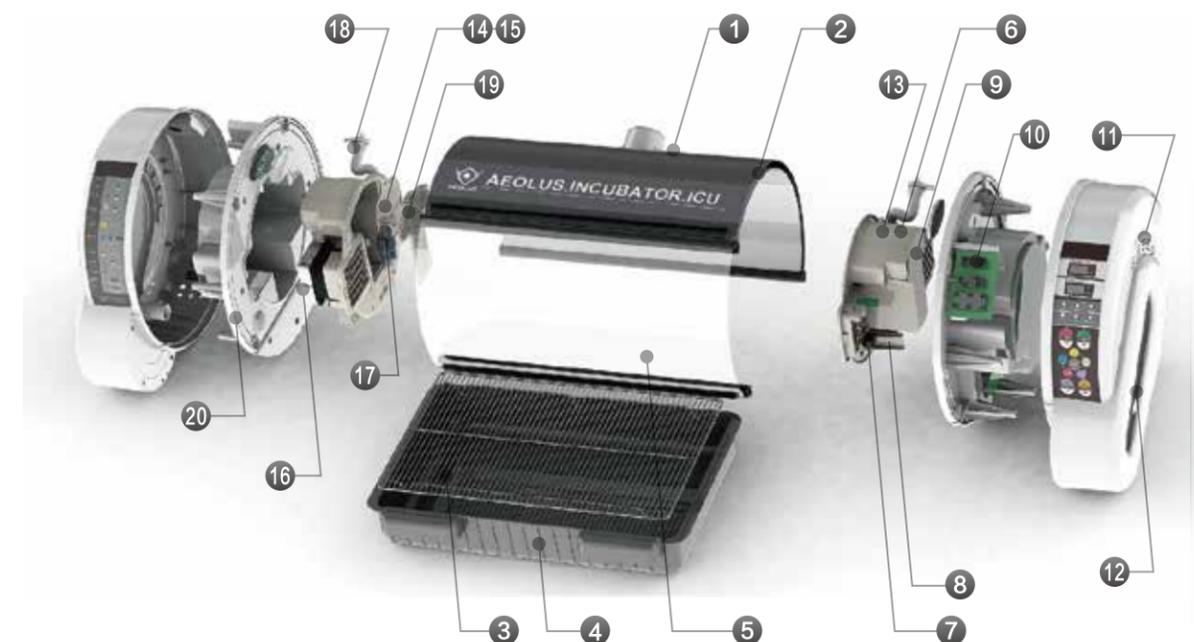

- | | |
|---|-------------------------------------|
| 1 Air Conditioning Port | 11 Air Vents |
| 2 Double Layer Insulation | 12 Lighted Locking On-off Switch |
| 3 S/S Floor Grid | 13 10 Level Lighting System |
| 4 Pull-out Waste Tray | 14 Carbon Dioxide Detector |
| 5 Heat Resistant Translucent Door | 15 High Density Filters |
| 6 Dual Heating System | 16 Consistent Dual Flow Circulation |
| 7 Negative Ion Generator | 17 Forced Air Nebulizer |
| 8 Precision Temperature and Humidity Detector | 18 Oxygen Port |
| 9 Humidifier Water Tank | 19 Ultraviolet Sterilization Light |
| 10 Gold-plated Contacts on Circuit Board | 20 Infusion & Oxygen Therapy Port |

Risk Free Operation

Our exclusive dual flow circulation system prevents any concerns regarding overheating by consistently and gently ventilating the ICU unit. Based on unique feature of PTC (porcelain temperature control) heating element, overheating is eliminated.

While in operation mode, the dual port ventilation system constantly refreshes the air inside. The CO₂ monitoring sensor alarm will be triggered immediately once the concentration of carbon dioxide rises above the set level (Default at 1000ppm).

Therapeutic Benefits

Dramatically raises the survival rate and recovery rate of patients. Should be used in various conditions, such as trauma, organ failure and transplantation, recovery from medical procedures and anesthesia, severe shock, septicemia, intoxication, severe skin disease, nurturing of newborns and Intensive care and treatment for numerous diseases.

Financial Benefits

Low maintenance, energy efficient, economical to operate. Revenue generating –billing should be based on time of use construed hourly. Most machines pay for themselves in the first quarter of use (according to recent customer reports).



ICU actual application examples





UC-1803

The new economical version of the fabulous ICU unit will be launched in 2018! Keeping most of the technical advantages in the last successful generation, only minimizing the control system and wiring to lower the manufacturing cost. Providing much more cost efficiency for all the veterinarian practices and breeding kennels. Internal Nebulizer has been removed, extra external nebulizer is available as optional part.



ICU/Incubator Organizing Rack

UC-1801-RCK


Double decker design, well for two units of KaniKare UC-1801 ICUs. Specialized top shelf for air conditioner or oxygen generator locating. Stable and Sturdy, made from Gauge 17 steel panel, ensures long service life. Comes standard with casters, for easy moving and sterilization operation. Trendy looking but functional too, a good mate for ICU accessories organizing.



Air-Conditioning Unit

UC-1801-AC

Perfect match for the KaniStar ICU. Gives the ICU the ability to set the desired stable temperature for the patient by producing chilled airflow. Remotely controlled by the ICU. Amazingly quiet without startling pets.

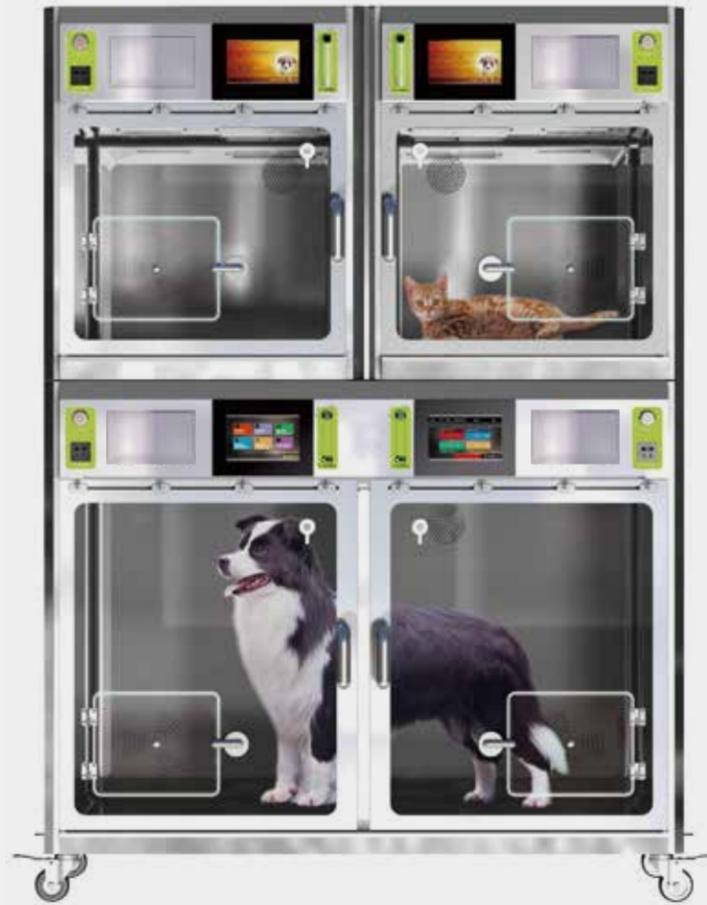

Stainless Steel ICU Unit | UC-1901

Widely used for Normothermia maintenance, fluid transfusion, oxygen supplement, blue light therapy, nebulization treatment, emergency rescue, and in house observation. Convenient self-diagnostic and all the components are easy to be replaced. Safer as a result of power supply failure protection and ventilation fan control. Creating a friendly medical circumstance by Ultra Violet Sterilization, High Concentration Negative Ion Air Purification, Temperature and Humidity Control, Carbon Dioxide Monitoring and Removal, Oxygen Concentration and Medical Nebulizer. Four separately controlled compartments provide specialized and individual cooling and heating. Two medium compartments on top tier, one large compartment on bottom tier which could be divided into two medium compartments when needed. Operating History Documenting Temperature, Humidity, CO₂ Density and O₂ Density, displayed in intuitive curves, providing up to 72 hours maximum tracking. Designed for the mounting all necessary accessories, such as IV pole, infusion pump, syringe pump, monitor, and door shade curtain.



5'W×3'D×6'H Approx.

- | | | |
|------------------------------------|---|-----------------------------------|
| 1 Heating & Cooling Equipment | 7 3 Way Illumination System | 12 Filter (Vent) |
| 2 Ultraviolet Sterilization Light | 8 Four Directional Heating & Ventilation System | 13 Built-in Printer |
| 3 Vent Door | 9 Filter (Vent) | 14 Nebulization Compartment |
| 4 Functional Door | 10 Negative Ion Generator | 15 Smart Control LCD Touch Screen |
| 5 Carbon Dioxide Monitoring System | 11 Emergency Vent | 16 Oxygen Flow Meter |
| 6 Central Divider | | |



Setting Interface



Trouble Shooting




Operation Guide




Running Status



Setting Parameters

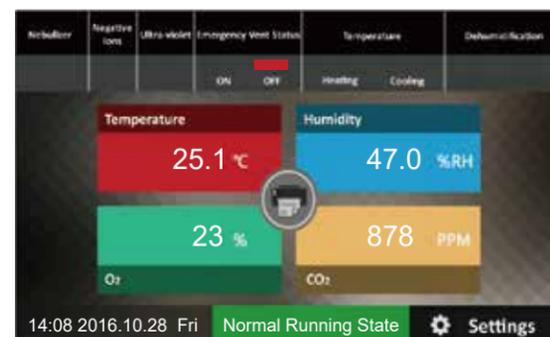

Main Interface

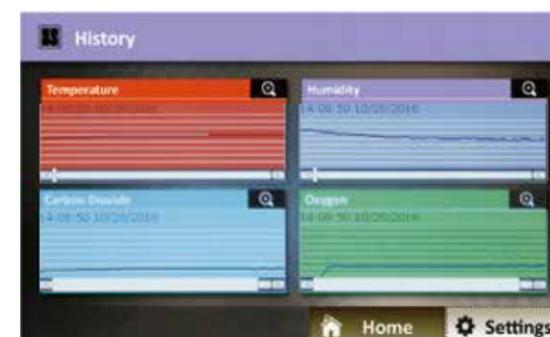

History



Fault Navigation



10-stage illumination light



Blue light therapy



Examination white light



Ultraviolet sterilization light



Air conditioning evaporation tank



Safe printer



Built-in nebulizer



Emergency vent



Ventilation intake



LED touch screen control



Stainless Steel Intensive Care Unit

An inexpensive but still professional solution for animal body temperature management or for intensive cares service.

Precision temperature control with only 0.1 degrees tolerance provides a safe and comfortable warm environment for sick animals.

Full sealed front door with oxygen inlet access, very helpful for oxygen treatment.

Controllable ventilation on front door too, allowing fast fresh air ventilation when needed.

Digital real time temperature and humidity display on front door too.

Two sizes available now, medium measures 65cmL-66cmD-79cmH and large is 130cmL-66cmD-100cmH, with varied combinations.




Veterinary Warming System | WMS-1501

Anesthetized patients are unable to regulate their own temperatures.

Hypothermia could be caused by core body temperature drops, which could increase the risk of complications. Including but not limited to uncontrolled shivering, wound infections, heart problems, the need for blood transfusion; all of these problems could lead to longer hospital stays. Even worse, it could lead to higher risk of death.

KANISTAR Warmology will provide a consistent sterilized warm atmosphere, surrounding the patient reducing the chance of Hypothermia (resulting from the administering of anaesthetics during surgeries).


WMS-1501 WARMOLOGY HEAT CORE

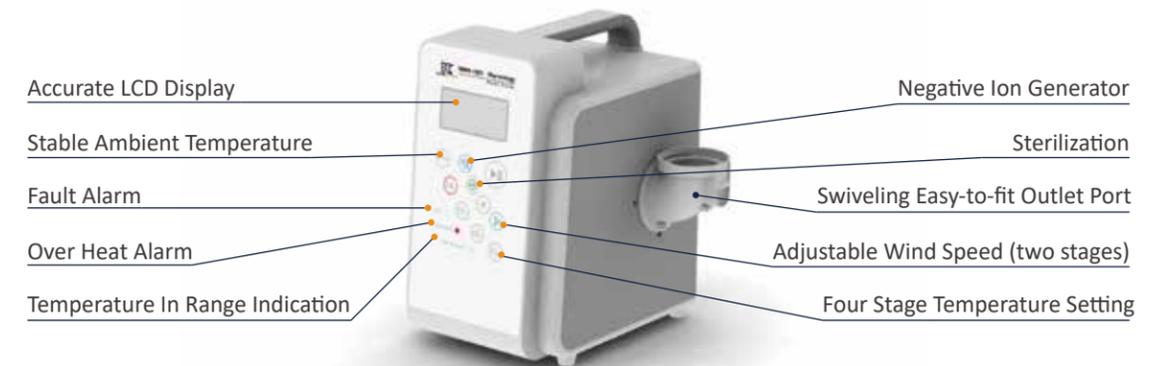

- Accurate LCD Display
- Stable Ambient Temperature
- Fault Alarm
- Over Heat Alarm
- Temperature In Range Indication
- Negative Ion Generator
- Sterilization
- Swiveling Easy-to-fit Outlet Port
- Adjustable Wind Speed (two stages)
- Four Stage Temperature Setting


The main heating unit offers a multitude of beneficial functions and features.

1. The operating noise level is very low. (53DBA).
2. Three tiered sterilization system, 254nm ultra violet light, high concentration negative ion and 0.2µm HEPA filtering.
3. Dual detecting temperature control with 1.5 Celsius degree minus tolerance ensured. Two sensors are located at the out let from the heat unit also at the end of the hose, detecting and compensating for heat loss from the hose.
4. Adjustable wind speed (two stages) with stable temperature control.
5. Easy to handle and maintain, smart trouble shooting system installed.
6. Accurate LCD display, real time temperature of blanket monitoring.
7. The unit can be freestanding or easily attached to our specialized stand.
8. Compact and light weighted, using little space, easy to transport.
9. Four stage temperature settings, 93F°(34C°), 98F°(37C°), 104F° (40C°), 109 F° (43C°).
10. Swivelling easy to fit outlet port, makes attaching the hose effortless.

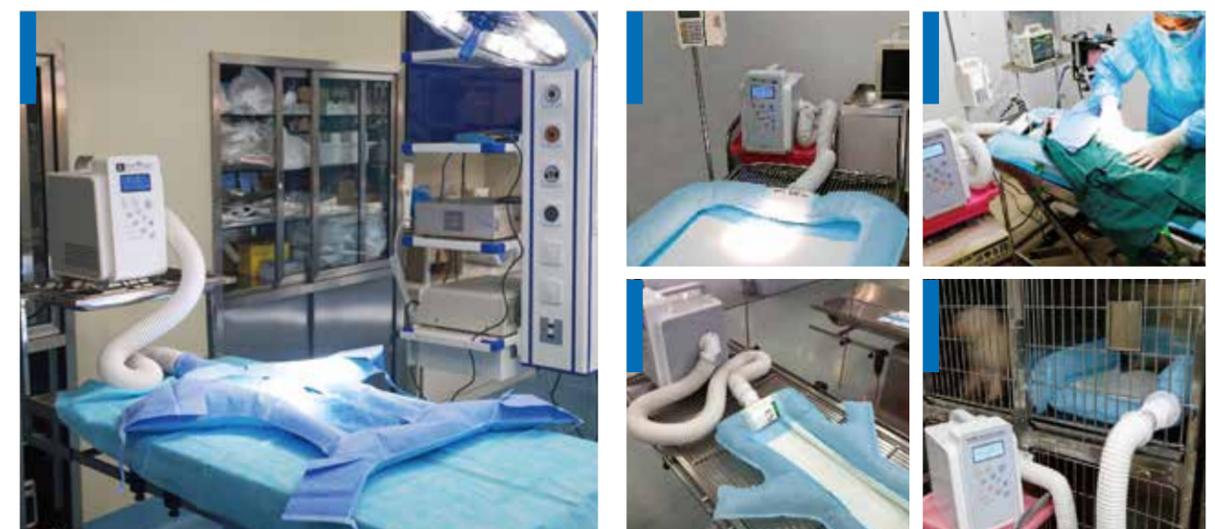

The system could be used before, during and after surgical or dentistry operations.

WMS-1501-ST WARMOLOGY BODYSUIT

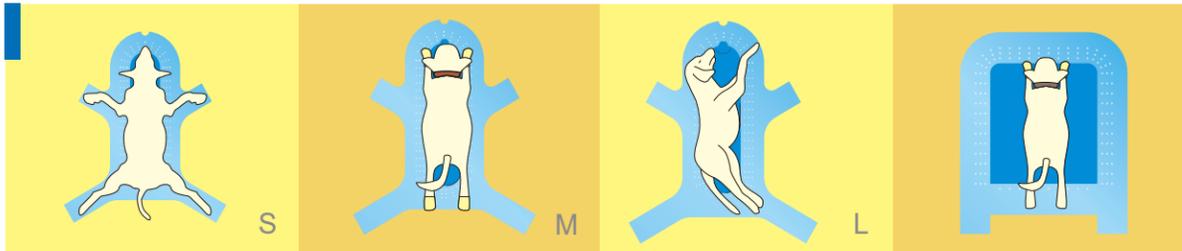

1. The body suit fully envelops the torso and limbs. (The limb wrap is our exclusive patented system).
2. Disposable nonwoven material, very cost effective.
3. Evenly distributes heat through 3D outlets system on wrapping "suit".
4. The best choice for anesthesia temperature management.
5. Different sizes available for small, medium and large animals.
6. Can be used with WMS-1501-MT mat to form a complete warming system.
7. Non-inflatable paw guards against thermal injury to the heat-sensitive paws and lower legs.
8. Durable and flexible materials that provide the ability to warm patients in various positions.



WMS-1501-MT WARMOLOGY MAT

1. The best solution available for convalescent animal normothermia maintenance.
2. Ventilation and fluid transfusion, allows the fabric to remain dry and warm for extended periods of treatment. They allow the patients to avoid the pain and discomfort of eczema and bedsores (pressure ulcers). Urine, feces and other bodily fluids always irritate the patient's skin and cause skin problems, especially when patient's activities have been restricted in a cage where the mat or floor is wet for longer than normal. This mat is an ideal choice to help create a dry and warm environment for animals. This significantly helps the healing process.
3. Made from patent pending fabric with cubic honeycomb structure, provides 3D support to the patient's body (40 points per square centimeter), filtering dirt in dense fibers.



1. Flexible, easy moulded to different sized animals.
2. The sterilized warm air will provide much better ventilation also significantly lower the risk of infection.
3. A very good choice for young animal incubation also.
4. Machine washable, easy to clean, eco-friendly (SGS approved).
5. No electric hazard risk, (no electric wires are used internally), preventing possible injury from exposure to liquids and chewing.





Headlight

TX-HDL-01

3W international brand LED with high brightness and long life.
Perfect optical structure, uniform illumination, clear and ultrabright spot.
High capacity battery heat dissipation design, and the silica gel soft headband.
Come with safety, comfort and ability to be disinfected.
Ultra-light and decent design.
Suitable for cardiothoracic surgery, brain surgery, ENT, dental and plastic surgery.



Headlight

TX-HDL-04

3W LED bulb with high brightness.
Perfect optical structure, uniform illumination, clear and ultrabright spot.
The dual high capacity battery supply easy and quick to replace.
Stable and reliable performance long lifetime.
Perfect optical structure, clear and bright light spot.
Widely used for local lighting in ENT, dental micro plastic operations.



Headlight

TX-HDL-02

5W world-famous LED with high brightness and CRI, nearly sunlight color temperature.
Reliable dual high capacity battery supply easy and quick to replace.
Perfect optical structure, uniform illumination, clear and ultrabright spot.
Equipped with a pendant of magnifying glass under the lamp holder.
Widely used for cardiothoracic, cardiovascular surgery neurological surgery, ENT, dental, plastic surgery brain surgery.



Headlight

TX-HDL-05

World famous longlife 5W LED bulb with high brightness and CRI.
High capacity battery with stable and reliable performance.
Perfect optical structure, uniform illumination, clear and bright light spot heat dissipation design.
Easy and comfort with the silica gel headband.
Ability to be disinfected.
No low frequency interference.
Widely used for ENT dental and other operations which has high requirements for local lighting.



Headlight

TX-HDL-03

5W world-famous LED with high brightness.
Perfect optical structure uniform illumination, clear and bright light spot
High capacity lithium battery power supply.
Workable while charging.
Light spot size adjustable by the ring on the lamp holder.
Made of high quality lightweight material comfortable and meeting medical disinfection requirements.
2200mAh battery with a work time up to 4 hours.
Very cost-effective and widely used for examination and operation of local lighting in ENT dental, general surgery micro plastic.



TX03-500C-LED

TX03-700C-LED

Multi-functional Surgical Operation Light

TX03-500-LED
TX03-700-LED

New light temperature adjusting technique, with 3-stage setting at 3200k, 4000k and 5000k.

Cold light source powered by long lasting LED components.

Unique “Endoscope” mode, which could lower the brightness significantly and quickly to adapt the endoscope operation needs.

Easy access left for installing the video camera onto this light by replacing the existing holding handle with camera assembly, as a good way for upgrading.

Trendier design, more powerful functions.

Two diameters available, 500mm and 700mm.



Double light combination



Ceiling mount version



Wall mount and free standing version available




Optional battery pack for free standing version



Easy upgrading by installing HD video camera



Adjustable brightness, 3 stage color temperature
Special Endoscope mode
Video camera control panel

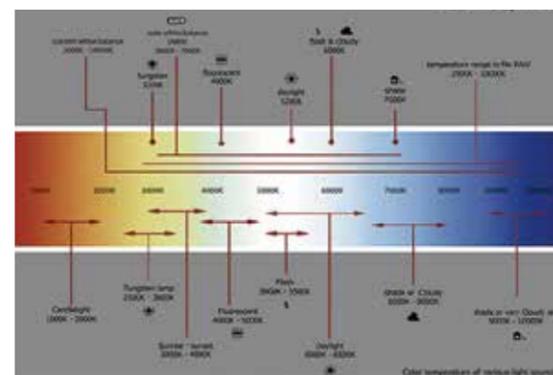

3-stage Light Temperature Adjusting, at 3200k, 4000k and 5000k.

Electrical Specification	
Description	700/500 LED
Supply Voltage	AC 220V±10%
Supply Frequency	50Hz±1Hz
Power Input	140VA
Bulb Power	3.2VIW
Main Fuse	F2AL-250V(quick-acting fuse)
Performance Specification	
Description	700/500 LED
Rated Illumination	100,000-160,000lx/80,000-130,000lx
Color Temperature	3800±200K/4400±200K/5000±200K



TX02-LED3+5





Shadowless Operation Lamp | TX02-LED

TX02-LED shadowless operation light is well suited for medical facilities. The lamp features fully closed streamline single or double lamp block structures.


The computerized integrated reflective optic system provides deep illumination and brings correct color temperature to the area of interests.

Microcomputer controlled, 6 brightness levels, memory function.


Operates in wide voltage range with power supply negative feedback alarm.

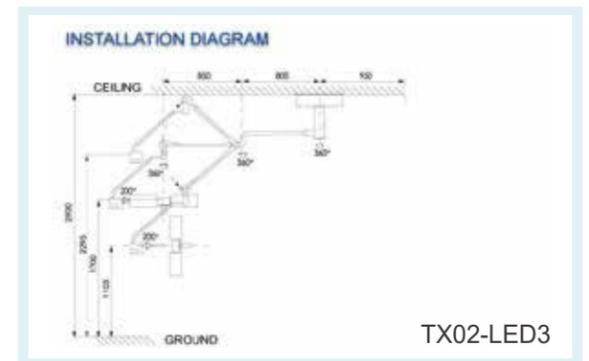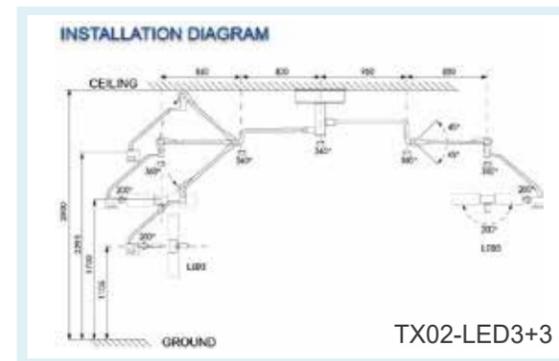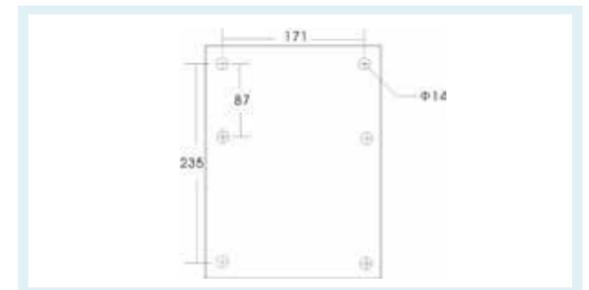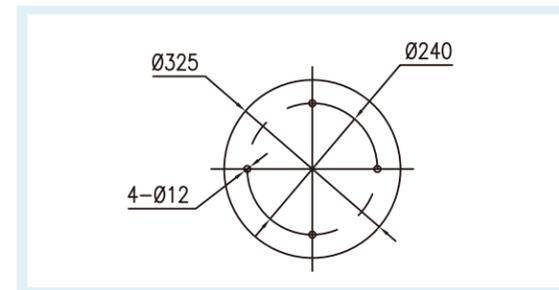

TX02-LED3 has two versions, wall mount and ceiling mount, distinguished by the suffix of W or C. The number next to LED indicates the quantity of lamp blocks in each group; where 3 stands for 3 blocks, 4 is for 4 blocks, and so on. TX02-LED3-W refers to "wall mount TX02-LED" shadowless light with 3 lamp blocks.

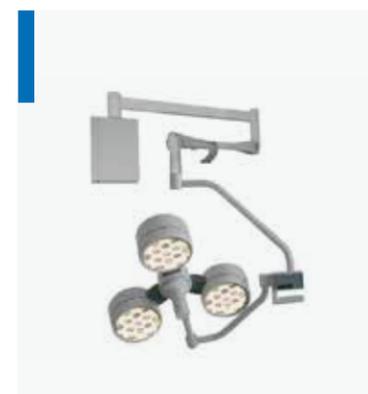

TX02-LED3

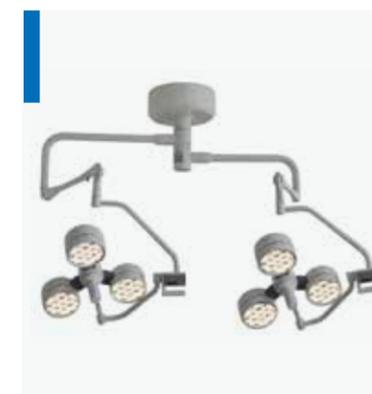

TX02-LED3+5



TX02-LED5-TV

TX02-LED3	
Rated illumination	80,000lux
Color temperature	4800±200K
Temperature rise	310W/m ²
Lighting depth	Not less than 600mm
Rotation of rotary body around central axis	360°
Rotation of balancer around rotary body	360°
Upward/downward turning	80°
Rotation of holder around balancer	360°
Rotation of holder block around connection tube	180°
Rotation of lamp around connection tube	180°

ITEM	Net Weight	Gross Weight	Dimension	Volume
	(kg)	(kg)	W×D×H(cm)	(cbm)
TX02-LED3	35	46	108×42×59H	0.52
	9	24	85×82×36H	



TX01-5LED



TX01-5LED-E

Cold Light Operation Lamp | TX01-5LED

TX01-5LED has three versions, free standing, wall mount and ceiling mount, with the suffix of W and C to indicate Wall and Ceiling Mount.

Adopting LED as light source, without causing surrounding temperature increase.


Excellent daylight quality and CRI.

Homogenous light spot, high-definition.

Energy saving and environment-friendly.

Flexible and stable balanced-arm support frame.

Suffix "E" refers to battery pack.

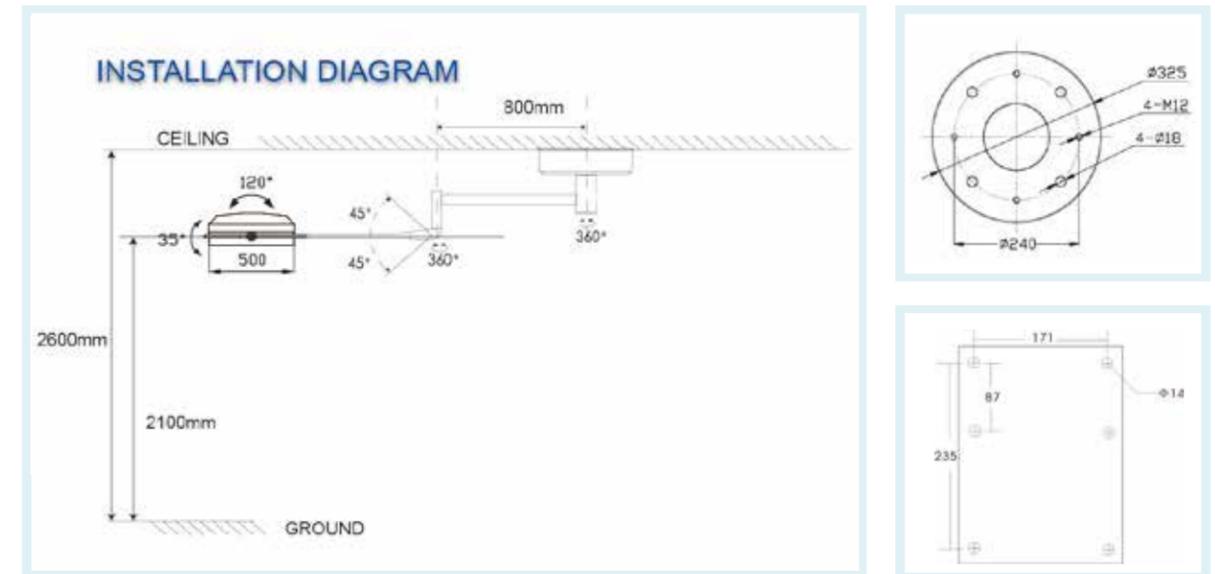

TX01-5LED-C



TX01-5LED-W




TX01-5LED	
Light intensity	≥100,000Lux/
Color temperature	4800±200K
Color-rendering index	92
Irradiance	310W/M2/420W/M2
Light field size	160mm~220mm
Lighting depth	500mm
Diameter of lamp holder	500mm
LED parameter	3.2V/W

ITEM	Net Weight	Gross Weight	Dimension	Volume
	(kg)	(kg)	WxDxH(cm)	(cbm)
TX01-5LED	42	53	114×64×33	0.24
TX01-5LED-W	37	47	111×66×39	0.28
TX01-5LED-C	26	42	90×70×47	0.29



Vertical Examination Lamp | TX01A-LED

TX01A-LED Stand Mobile Examination Lamp.
 A convenient, energy saving, mobile examination light.
 LED light source does not cause any unnecessary temperature rise.
 High luminance, good for mobile exam needs.

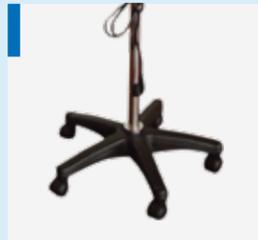

TX01A-LED				
Voltage	AC 110/220V,50-60Hz			
Property	AC			
Lamp power	3.3V/3W			
Color temperature	4800±200K			
Illuminance	30000lux 0.5m/10000lux 1m			
LED Life	20000hours			
ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension WxDxH(cm)	Volume (cbm)
TX01A-LED	6	8	96×55×17H	0.09




Spot Exam Light | TX01B-LED

Easily maneuverable arm allows comfortable up/down control.
 Positioning handle mounted on the body of the light.
 LED cold light source, long service life.
 Offers exceptional performance and versatility.



TX01B-LED				
Voltage	AC 110/220V,50-60Hz			
Power Supply	AC			
Lamp power	3.2V/1W			
Color temperature	3500±500K			
Illuminance	≥70,000 Lux			
Lamp head diameter	160-220mm			
ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension WxDxH(cm)	Volume (cbm)
TX01B-LED	7	9	96×65×32H	0.20

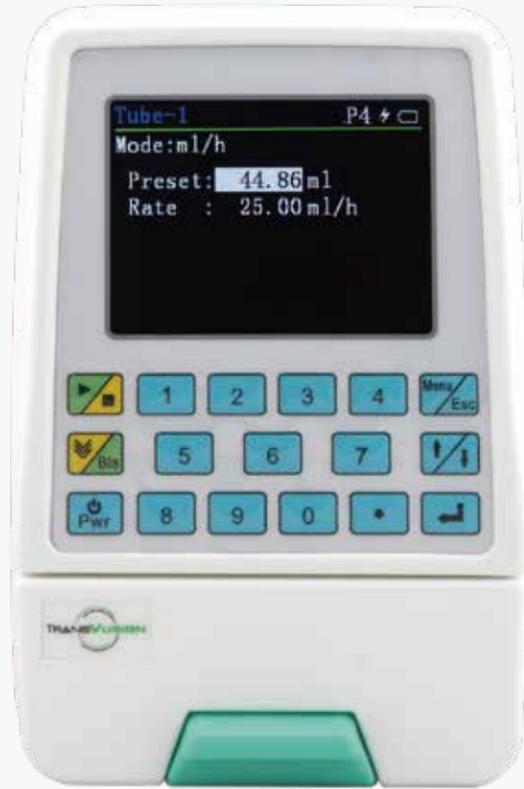

Intelligent Infusion Pump with Heater | TX-LP-2100H/2100HW

1. LCD color display, vivid and intelligent.
2. Build-in warming module to preheat medicine fluid, softening infusion line and prevent hypothermia (low temperature) during infusion.
3. Unique semi-extrusion technique, good for all kinds of fluid transfusion, including medical solution, blood, plasma and also nutritional liquid.
4. New industrial design, making this unit lighter weighted (only 1.3kg), more compact, achieving a real handheld machine.
5. Available with optional Wi-Fi supervision function. With Wi-Fi receiver module installed on computer, up to 60 units of pump could be monitored real-timely.
6. Safer bottom-hinged swinging door to avoid any core-part damage caused by possible liquid leakage.



Functions and Parameters

1	Pump type	Finger-type peristaltic pump
2	Applicable infusion sets	Special or quality PVC,TPE infusion sets with an external diameter of 3.4mm~4.5mm.
3	Infusion mode	Flow-rate mode, time mode, drip-rate mode, body-weight mode.
4	Flow-rate mode	Flow-rate setting: Between (1.0~1200.0)ml/h or (1.0~1200.0)ml/h (calculation result) at 0.1ml/h per step
5	Preset time mode	Preset time setting: Set between (1s~99h59m59s) at 1s per step.
6	Drip-rate mode	Drip-rate setting: Set between 1~333d/min at 1d/min per step.
7	Body-weight mode	Body-weight setting: Set between (1~300)kg at 1kg per step.
		Mg dosage setting : Set between (0.01~999.99)mg/kg/h at 0.01mg/kg/h per step
		Infusion dose setting : Set between (0.1~9999.9)mg at 0.1mg per step
8	Flow-rate setting for KVO function	Set between (1.0~5.0)ml/h at 0.1ml/h per step.
9	Alarm	1.Air in tube, graded sensitivity adjustable
		2.Bubble sensor fault
		3.1.Upper-end occlusion,-20~-60kPa, graded sensitivity adjustable
		3.2.Upper-end base-pressure fault(pressure sensor broken)
		4.1.Lower-end occlusion, -20~-60kPa, graded sensitivity adjustable
		4.2.Lower-end base-pressure fault(pressure sensor broken)
		5.1.Pump rotation fault,
		5.2.Inverse-infusion,
		6.Abnormal door-opening,
		7.Low battery power, battery fault
		8.Overtime pause,
		9.Completed infusion
10	Data storage	Memory calculation: Auto-record infusion set data for late usage
		Setting record: Auto-record last infusion setting for late usage
		History(SOE): up to 1000 latest infusion events can be recorded(start time, end time, flow-rate, stop reason and the like)
11	Safety classification	Class I, CF device,IPX4
12	Power supply	AC power:100-240V,50/60Hz,built in rechargeable NI-MH battery
		Power consumption: ≤20VA
13	Shape and weight	W100×H134×D80mm; ±5mm, Net weight≤1.3kg
14	Heating function	Built-in heating module
		Built-in heating function
		Optional external heating belt.
15	Wifi module	Built-in wireless data transmitter
		Wireless networking(optional signal receiver installed on computer)
16	Screen	LCD color screen, resolution:320×240
17	Special function	Alarm stops when occlusion disappears
4. 5. 6. 7 : Preset volume setting: Between (0.1~9999.9)ml at 0.1ml per step Display of accumulated volume : (0~9999.9)ml, "zero clearing" only and modification is unavailable		



Infusion Pump | TX-LP-1900

Reliable Patency-Restore function, working very well for our hair patients. Breakthrough semiextrusion technique, do not damage molecules in fluid, good for medical solution, blood, plasma and also nutritional liquid.

Patented "Double Pulse" compensation technique ensures very good linearity even at low delivery rates. Double CPU system, makes the running status super safe. Delivery will be cut off even control circuit is powered off and the machine will alarm. Rechargeable back-up battery supports at 4 hours continuous operation.

Different modes: Delivery Rate, Delivery Time, Patient Weight and Dripping Rate.

Alarm for door open, low battery, occlusion, end program, malfunction, etc.

Light weighted compact size, measures only 128mmW*15mmD*300mmH.

Historical date memory function, keeping up to 1,000 records that can be saved for over 5 years without external power or built-in battery.

Available with optional fluid warmer and Wi-Fi supervision function. With Wi-Fi receiver module installed on computer, up to 60 units of pump could be monitored real-timely.



Patency Restore Function



Semi-extrusion function makes it not merely an infusion pump, but also good for plasma/full blood and nutritional liquid transfusion.



Optional preheat module



Optional Wifi Supervision



Wifi Supervision Interface



Historical data recording function, up to 1000 records.



Wifi Receiver



<https://youtu.be/CGTj8yGbcL4>
Watch our AD video on Youtube!





KO2 Kanoxy Oxygen Concentrator

| KO2-2F

LCD display, allowing easy and accurate air flow management.
 Steady 2L flow rate at 95%+ concentration.
 Extremely quiet; running at only 43 decibels!
 USA made molecular sieve and long lasting compressor.



Control Panel



Oxygen Outlet Humidifier



Filter and Power Switch



Filter Cleaning and Replacement



Remote Control



Well-designed Package



KO2 Kanoxy Oxygen Concentrator

| KO2-10F

A true 10L capacity oxygen generator!
 Extracts medical-grade oxygen from the air and delivers it on demand.
 Steady and safe 10L oxygen supply, an economical and long lasting solution for most medical care needs.
 Capable enough for anesthesia machine running.
 Very quiet compressor; utilizing exceptionally fine quality German technology; combined with USA made molecular sieve, will provide long term service.
 Features standard fittings.
 UL Tested.



Humidifier



Display and Control Panel



I-stage Filter



Omni-directional Casters



TX-OTS-01



TX-OTS-02

Diagnostic Set - Otoscope and Ophthalmoscope with Rechargeable Handle

3.5V Fiber Optic Oscopes

Offers you clear images of the ear canal and tympanic membrane;
Wide-angle viewing lens allows for instrumentation under magnification;
Sealed system for pneumatic otoscopy and accurate diagnosis of otitis media;
Fiber optics project cool light to the distal tip with no reflections or obstructions for superior views of the tympanic membrane;
Lightweight, reliable and durable Reusable Ear Specula (2.7 and 4.2mm).

3.5V Coaxial Ophthalmoscope

Coaxial optics produce a shadow-free spot, easier entry into undilated pupils, and a larger field of view versus standard ophthalmoscopes;
18 unique aperture/filter combinations for greater versatility;
28 focusing lenses with a range of -25 to +40 diopters.



Diagnostic Set with Lightweight Rechargeable Handle

TX-OTS-01

Half the weight of traditional 3.5 V handles, and fits more comfortably in your hand.

Fits all our 3.5 V instrument heads;
Made of sturdy heat-resistant engineering resin ABS;
Stepless brightness regulation.



Diagnostic Set with Metal-case Rechargeable Handle

TX-OTS-02

Replacement Lithium Ion Battery included;

Fits all our 3.5 V instrument heads;
Aluminum construction for durability and reliable long life;
Smooth, nonabrasive knurled finish ensures a no-slip grip;
Stepless brightness regulation.



Cannulated Bone Drill | TX-BDS-01

Ideal tool for K-wire, intramedullary nail, and endoscopic surgery.
 Full aluminum alloy shell, high-temperature disinfection available at 135 degrees Celsius and 1.6Mpa.
 Only 20 minutes' disinfection will kill viruses and bacteria.
 Brushless motor with forward and reverse stepless variable speed, maximum speed 1000rpm.
 High coaxial degree, the radial offset is within 0.1mm.
 Low-noise operation. The noise of no-load drill is 65dB(A).
 Most prevailing drill bit of 0.8mm-8.1mm available.
 Supply voltage: 13.2V. Two battery cases included and each can last more than 2 hours.
 Lighted weighted and easy to carry. Safe and simple to operate.
 Full line of bone drill tools, including common bone drill, oscillating saw, reciprocating saw, overhaul saw, grinding drill and complete accessories for sale.



Aluminum alloy toolcase



Widely adaptive collet



Charger



Long-lasting battery



Aluminum alloy disinfecting channel



Exquisite craftsmanship



Cannulated Bone Drill

TX-BDS-02

Internal fixation



Swing Saw

TX-BDS-04

Mainly used for joint surgery in knee surgery



TPLO Saw

TX-BDS-03

Horizontal osteotomy of tibial plateau in animals



The Universal Swing Saw

TX-BDS-05

Limb joint surgery





The Multi-Functional Drill Saw | TX-BDS-06

Limb trauma joint surgery



Portable Intelligent Electrosurgical Unit | ESU-1000CI

Applicable to Cutting and Coagulation of all kind of operations;
 Automatic Boot Check-up and Error Detection;
 Output power control by using the adjustment button of the ESU for a faster operation;
 Memory function with a capacity of 10 settings, restart with the previous power setting;
 Two levels of coagulation: pinpoint and fulguration. Pinpoint coagulation provides precise control of bleeding in localized areas. Fulguration (spray) coagulation provides greater control of bleeding in highly vascular tissues over broad surface areas.
 When the power is 10w or less, you can adjust in the amplitude of 0. 1W for spray coagulation;
 Work with the Argon system to be a complete set of Argo Plasma System;
 For your convenience, the power output is calibrated in watts with large, illuminated digital displays;
 The output terminals are independent.



Dental Delivery System | DTL-681

Dental delivery system DTL-681

An integrated dental tools driving and delivery system for you.

Providing a clean, dry and aseptic work environment.

Designed for delivering water, compressed air, and electricity to dental instruments.

Allowing four input lines, good for one high speed drill, one low speed drill, polisher and scaler, etc.

Quality components. Safe, reliable, convenient and comfortable.

Quiet running, without irritating pets.

Ergonomical design, simple installation and operation.

Easy to maintain, disinfect and clean.



Complete Dental Operation Station (FT-851 dental table and DTL-681 delivery system)



Side view



Quality components



Tools tray



Handle



Smooth movement

Dental handpiece

Imported ball bearing;
 Safe, reliable, high precision;
 Low noise (Noise: $\leq 65\text{dB}$), big torque and vibration free;
 Sterilization by autoclave at $(136\pm 2)^\circ\text{C}$ & 0.22MPa for 15 minutes is recommended;
 Easy for lubrication and sterilization;
 Convenient to change burs.
 Water/mist cooling sprayer.

Low speed handpiece

Low speed handpiece is applied to trim, model trim and so on.
 Contra angle head and straight head;
 Inverting function of motor rotation direction;
 Torque $\geq 3\text{N}\cdot\text{cm}$.



Inner cooling waterway



Outer cooling waterway

High speed handpiece

During dental treatment, high speed handpiece is used for drilling and trimming while matched the high speed bur.
 Push button structure convenient for changing bur;
 Type of connection: Quick coupling or non quick coupling;
 Anti-retraction system preventing from cross infection;
 Air pressure: $0.25\text{Mpa}\sim 0.27\text{Mpa}$ Airflow speed: $30\sim 36\text{L}/\text{min}$
 Cooling waterflow rate: $> 55\text{mL}/\text{min}$ at water pressure of 0.2MPa
 Rotation speed: $\geq 310,000\sim 410,000\text{r}/\text{min}$
 Cutting force: radial $\geq 1.0\text{kgf}$ axial $\geq 2.5\text{kgf}$
 Dynamic balance for cartridge $\leq 180\mu\text{g}$
 Radial runout $\leq 0.01\text{mm}$
 Noise: $\leq 65\text{dB}$
 Torque: $0.12\sim 0.13\text{N}\cdot\text{cm}$
 Cooling waterway: inner



Mini head for small pet.
 Head size: $\phi 10.2 \times 13.1$
 German made ball bearing
 360° swivel quick coupling.



German made ball bearing
 360° swivel quick coupling.



Head size: $\phi 12.2 \times 13.7$;
 Imported ball bearing.

Slow handpiece

Used for dental implants.
 6 sets of ball bearings;
 Speed reduction ratio: 20:1;
 Head size: $\phi 2.35\text{mm}$;
 Rotation speed: $0\sim 2000\text{rpm}$;
 Cooling waterway: outer.



Silent oil-free air compressor

The excellent silent oil-free air compressor can provide clean, powerful and silent compressed air source. The compressor can be customized according to the number of the pneumatic dental units, up to 60 units, which are connected to the compressor.
 For all models, air dryer system and silent cabinet are as options.



Advantages

High quality air compressor pump, with pure copper coil;
 Tank inside has done anti-rust and antibacterial coated treatment;
 Longer service life: 7-8 years lifespan;
 Working noise is 5-10 dB less than normal oil free compressor supplier;
 Easy operation, no more maintenance fee;
 Meticulous workmanship;
 Low vibration: rubber feet limit its vibration to a lower level;
 Fashionable design and durable material;
 Low energy consumption air compressor;
 High precision filtration.



Compressor with air dryer




Compressor with silence cabinet



Mobile Dental Cart | DTL-P302S

Electric lifting and easy movement with silent medical casters and storage space for foot pedal;
 Decent looking with Acrylic top;
 Ergonomic design with rotatable tool hanger and arm;
 Medical grade air and water pipe, water bottle, drainage bottle, aluminum alloy bracket, and high quality pure copper valve;
 Standard accessories:
 1pc 3-way syringe, 1pc high volume ejector, 2pcs handpiece tube (4 holes), 1pc LED curing light, 1pc ultrasonic scaler, Self-contained water bottle, Self-contained drainage bottle, Electric height adjustment system (560-800mm), Rotatable operating hanger (Radius of rotation: 300mm), Foot control system, 110~240V, 50/60Hz, adaptive plug to fit your country.
 Optional accessories:
 High and slow speed handpiece, Fiber optic handpiece tube, LED woodpecker ultrasonic scaler, Silent dental air compressor.

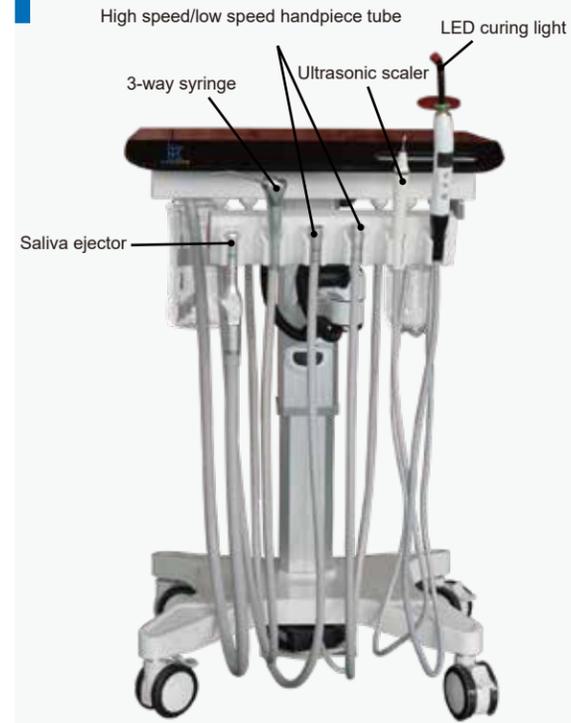

Accessories



Portable silent dental air compressor



Built-in Ultrasonic Scaler

UDS-N2 LED

Built in ultrasonic scaler used for dental calculus elimination and root canal treatment;
 Special design for dental unit;
 Detachable handpiece with light, better visibility and autoclavable ;
 Automatic frequency tracking.
 No cracking or deformation under the pressure for shell of the handpiece is made of special imported plastic from America.
 Output power: 3W to 20W
 Output primary tip vibration excursion: $\leq 100\mu\text{m}$
 Output half-excursion force: $< 2\text{N}$
 Output tip Vibration frequency: $28\text{kHz} \pm 3\text{kHz}$



Built-in unit



Tips



Handpiece


Intelligent power output;
 Experience painless treatment;
 Integrate subgingival scaling, reconstruction of minimal damnification, implant maintenance and root canal therapy into a whole;
 Intelligent negative feedback--- More effective scaling :
 In accordance with different tooth surface, intelligent degenerative feedback will automatically output accurate power to the tip. Stable and sustained scaling in the calcification removing process brings you professional operating experience.
 Output power: 3W to 20W
 Output primary tip vibration excursion: $\leq 90\mu\text{m}$
 Output half-excursion force: $< 2\text{N}$
 Output tip Vibration frequency: $28\text{kHz} \pm 3\text{kHz}$



Intelligent touching system: Clear control panel, Ultra-sensitive and Waterproof design.



Endo: Fast & Complete: After ultrasonic irrigation, there is no smear layer surface but a few amount of dentin debris in the root canal.

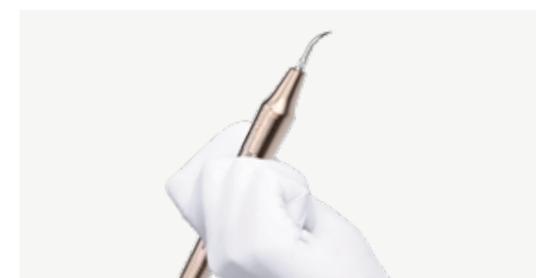

Sine wave drive: Avoiding stimulation caused by instantaneous current output; ensuring comfortable treatment.



Automatic water Supply System.
 Easy to change water, Increase work efficiency.
 Hydrogen peroxide, sodium hypochlorite, and chlorhexidine is available.
 Clear 50ml scale mark makes it convenient for formulating liquid.

Ultrasonic Scaler

U600 LED

Used for dental calculus elimination and root canal treatment;
 Detachable and autoclaveable;
 Automatic frequency tracking ensures that the machine always works on the best frequency and more steadily;
 Digital control, easy operation and more efficient for scaling.
 Scientific strength with multiple functionalities





UV Sterilization Trolley | TX-UDV-01

Designed to reduce viruses, mold spores, yeast, and harmful bacteria.

Using two high-powered UV-C lamps with 253.7nm rays to emit concentrated rays the machine is able to sterilize rooms with various surfaces.

Entirely based on ultraviolet light, it eliminates harmful microbes without using harmful chemicals.

915mm lamp arm length for maximum ray emission and optimal reach over surfaces

The lamp arm with 180° range of motion can be adjusted to rest at whatever angle and best suit the specifications of your room or surface.

Simple and Safe to Use, easy for storage too.



Front view



Back view



Inside view



Lamp holder and lamp protective cover




MEDICAL DISINFECTANT SPRAY

EPA HOSPITAL-GRADE DISINFECTANT, imported from USA and distributed in China, eliminates bacteria, viruses, mold and mildew safely and effectively.

OUR SOLUTION IS YOUR SOLUTION FOR INFECTION PREVENTION & CONTROL, A FRAGRANCE-FREE ODOR ELIMINATOR AND ECOLOGICALLY SAFE.

20+ kinds of bacteria and viruses are killed, including 99.999% of those causing MRSA, Norovirus, H1N1, E. Coli, Salmonella and Listeria.

CLEANS, KILLS & CONTROLS PARVO*, DISTEMPER*, BORDETELLA*, RINGWORM* & MRSA

Kills and prevents mold and mildew for up to 7 months.

Eliminates odors at their source.

Proven allergen eliminator.

Antimicrobial control for HVAC, carpet and more. BACKED BY THE CARPET AND RUG INSTITUTE.

Natural and environmentally sustainable.

NSF certified food surface sanitizer—no rinse required.

Non-irritating to the skin; no gloves required, for use on colorfast fabrics.

Ready to use, no mixing required—just spray, wipe or fog.

Air dries quickly and leaves no odor or sticky residue.

Product lasts for 12-months after opening; 24-month shelf-life when stored at room-temperature.

SUPER STERILIZATION

Proven effective against Methicillin-resistant Staphylococcus aureus (MRSA), Norovirus, Legionella pneumophila, E. coli and the H1N1 virus (Swine Flu) in testing conducted by certified independent laboratories under GLP conditions submitted to the USEPA.

SUPER ODOR ELIMINATOR

Highly effective at eliminating pet urine and fecal odors as well as skunk, smoke, trash and more. The disinfectant solution contains no masking agents or fragrances and leaves your facility smelling clean and fresh.

MICROBE ELIMINATOR

The disinfectant solution is a broad-spectrum antimicrobial that performs against a wide variety of bacteria, yeasts, viruses, fungi, spores, molds, mildews, and other microbes. It is effective at low concentrations, over a wide pH range.

ALLERGEN ELIMINATOR



Proven effective on dust mite and cockroach allergens and pet dander. Can be used directly on bedding.

SOFT SURFACE SANITIZER

Great for synthetic pet bedding, carpet, upholstery and curtains—no soaps or detergents in the product, no residue, and doesn't damage fibers. Kills odor-causing bacteria on contact by eliminating the source.

LABOR PROTECTION

For your operating room, clinic, lab, shelter, boarding facilities, breeding kennels, doggie day care, grooming salon or pet store. Eliminates bacteria, mold, mildew and other fungi including Ringworm.* Disinfects, eliminates odors and ensures air quality in HVAC systems.

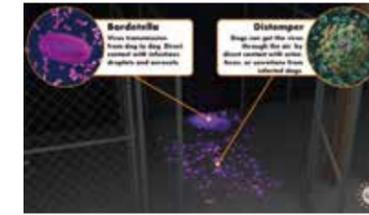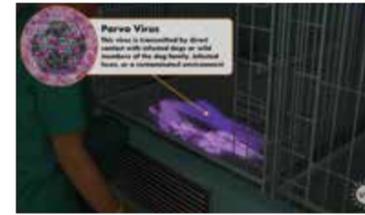

KILL LIST

Hard Surface Sanitization: 30 sec. contact time

Escherichia coli
Staphylococcus aureus

Carpet Sanitization: 10 min. contact time

Enterobacter aerogenes ATCC
Staphylococcus aureus ATCC

Soft Surface Sanitization: 5 min. contact time

Enterobacter aerogenes ATCC
Staphylococcus aureus ATCC

Disinfection Bacteria: 10 min. contact time

Pseudomonas Aeruginosa
Acinetobacter baumannii
Staphylococcus aureus MRSA
Staphylococcus aureus ATCC
Listeria monocytogenes
Legionella pneumophila
Salmonella enterica
Klebsiella pneumoniae (NDM-1)
Escherichia coli
Bordetella bronchiseptica m (Kennel Cough)*

Mildewstat (fungi): 10 min. contact time

Aspergillus niger – 1 wk. protection hard surfaces
4 wks. protection soft surfaces

Stachybotrys chartarum
Aspergillus fumigatus
Alternaria alternata

Penicillium sp
T. rubrum ATCC (Ringworm)*
T. mentagrophytes (Athlete's foot)*

Disinfection Viruses: 5 min. contact time

Respiratory syncytial virus (RSV)
Human immunodeficiency virus (HIV Type 1)

Rotavirus, Strain WA
Hepatitis A virus
Hepatitis B virus
Hepatitis C virus

Influenza B virus

Murine Norovirus (MNV-1)
Norovirus feline calicivirus
Swine influenza (H1N1)

Canine distemper virus*
Canine Parvovirus ATCC*
(10 min. contact time)
Adenovirus* (10 min. contact time)
Hantavirus*





OG-502

OG-501

Ozone Generator | OG-501/502

Grooming salon, pet shop, vet clinic, kennel, pet hotel, pet training school, etc. will all benefit from having this wonderful machine on hand. It removes odor while helping sterilize and disinfect the environment.

This advanced ozone machine is designed to destroy, not just cover up pet odors, cleaning chemicals odors, bacteria, fungus, viruses and other indoor air pollutants. Generating high level of ozone, 14g/h, it is one of the most efficient machines in current market. Purify the air in a 325 square ft. room in approximately 5 minutes!

Count-down timer and automatic shut-down function.

Capable of operating for extended periods, it enables you to maintain an odor free and clean environment.

Requires no supplementary additives and is relatively maintenance free.

120W, low carbon production.

User friendly.

OG -502

Single powder unit with simplified interface, easy to use and cost efficient.



Sterilization Chest | US-503

Utilizes Ultra violet sterilizing lights which operate at an incredibly efficient wavelength of 254nm. Double compartment design, lets you clean and store twice tools.



Pet Hospital Buying Guide

Introduction of Veterinary Operation Tables



Preparation Table

Your preparation and operating tables are the starting point and focal point around which you can build your clinic. It should welcome and draw patients and their owners in.

AEOLUS FT-851 and FT-851-AS series are ideal preparation table choices. The attractive white bottom cabinet offers storage for tools and supplies. It is also warm and inviting. You can choose between top quality Stainless Steel and the beautiful artificial stone basins. Both offer easy flow of waste liquid. The pre-drilled holes are ready to accept water taps (WK-301/302) installation.

FT-852 and **FT-861** electric tables are terrific choices that offer a different esthetic, perform quite well and have a smaller footprint. These tables both drop to very low height of less than 30cm height (less than 12"). This decreases the work load and effort needed by hospital/clinic staff, making it a bit easier to move big pets.



Operation Table

Veterinary operation tables encompass a broad spectrum of styles and uses, with many varied classifications.

Comprehensive Operation Table

Representing this classification are FT-826, FT-827 and FT-828.

Our flagship model **FT-828** is classic representation of an operation table. The angle of the large table top (145cm H / [54"]) is adjustable. The two table panels can be adjusted to V-type (with a center channel for fluid run off). The design aids in the prevention of asphyxiation due to regurgitation and choking. It enables a pet's head to rest in a safe position. The ergonomic width (23cm [10"]) of each table panel is deliberately based upon the length of an average forearm.

FT-826 is our unique design. This extension table was developed about 10 years ago by our team. The maximum width of table top can extend to 73cm (78") accommodating larger animals.

FT-827 also has a very distinctive feature; the table top rotates on the pedestal, turning smoothly on oversized bearings. The table top can revolve freely, allowing the doctor to remain comfortable and secure while repositioning the patient to complete the procedure.

Dental Operation Table

These operation tables are "recommended" for teeth cleaning and dental procedures. Most animals are under an anesthetic during the operation. Proper drainage is essential. A basin table is the preferred choice for oral care.

FT-851, **TT-202** and **FT-852** are good choices.

The FT-851's Dental Delivery System provides automatic control for handpieces, a scaler, a saliva ejector and an air-water syringe. The system includes air and water delivery plumbing system, rotatable handpiece hangers and an air brake system to ensure proper positioning and stability.

Optional accessories: 1pc air compressor, 1pc saliva ejector; 1 pc high speed handpiece, 1 pc low speed handpiece, 1pc ultrasonic scaler.

Ophthalmology Operation Table

Our reverse top dip tables provide doctors a more convenient surface to properly examine the details of animals' eyes. The trough collects waste liquid during the operation, then slides out (for easy cleaning). The V-top helps to guide waste liquid flow away from the operation area.

We offer many tables suitable for ophthalmology operations.

The **FT-853** is a good one to consider.

Extended Operation Table

New structural designs in some tables are bringing comfort while easing the workload. A vertical lifting column as implemented in **FT-88X series** provides multiple use options. The resulting extra leg space underneath allows for freedom of movement.

Let us help you fulfill all your clinic and hospital table needs.

Exam Table/Transfusion Table

These tables are used for basic operations; bandaging and changing dressing, small operations and simple suturing. Sturdy frames and table tops are accompanied by hooks for hanging transfusing fluids.

FT-823 and **FT-825** series are the classic exam tables. Our tables come in several sizes with hydraulic or electric push rod pedestals. These tables can be customized according to hospital space requirements and limitations.

We also have cost effective value tables, **FT-871** and **FT-872** available. We produce tables that will meet your needs and fall within your budget.

Tub Table

If an operation yields large amount of fluid loss or requires for fluids for cleaning, we recommend.

Fully rounded corner design, visible drop height of two sides, provide easy flow of waste liquid.

Scale Table

The weight of pet directly determines the accurate dosage of medication. When administering and prescribing medications to small pets, a slight miscalculation can endanger the animal. It is imperative to accurately weigh the pet to ensure proper dosages are administered. Safety is key.

FT-822 uses precision metallic sensors for accuracy in weighing. Weights ranging from 20g ($\frac{3}{4}$ oz.) to 100kg (220 lbs.) can be weighed without hesitation.

Our floor mounted walk-on scale **WOS-4824** is designed for larger animals.





Complete Dental Operation Station with Delivery System!

Preparation Table Dental Table with Stainless Steel Tub

FT-851


First Class Preparation Table.

Fashionable styling, practical function and detailed finish make the FT-851 preparation tables the best choice for your clinic. SUS304 stainless steel basin with slope top, flat surface grid and drainage on lower side provides very useful preparation space. Drawers and cabinet under top basin afford handy storage for the tools and supplies. Left side, right side and middle table configurations offer a multitude of space and design options.



1. Versatile cabinet/drawer combinations.
2. Right hand, left hand and middle versions.
3. Stainless steel hair strainer.
4. Top resting grid.
5. Round cornered tub with slope.
6. Comfortable legroom underneath round corners.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-851R/L	72.50	133.00	120×60×90	47-1/4×23-5/8 ×35-7/16	125× 65×104H	0.85
FT-851M	69.50	127.70	120×60×90	47-1/4×23-5/8 ×35-7/16	125× 65×104H	0.85

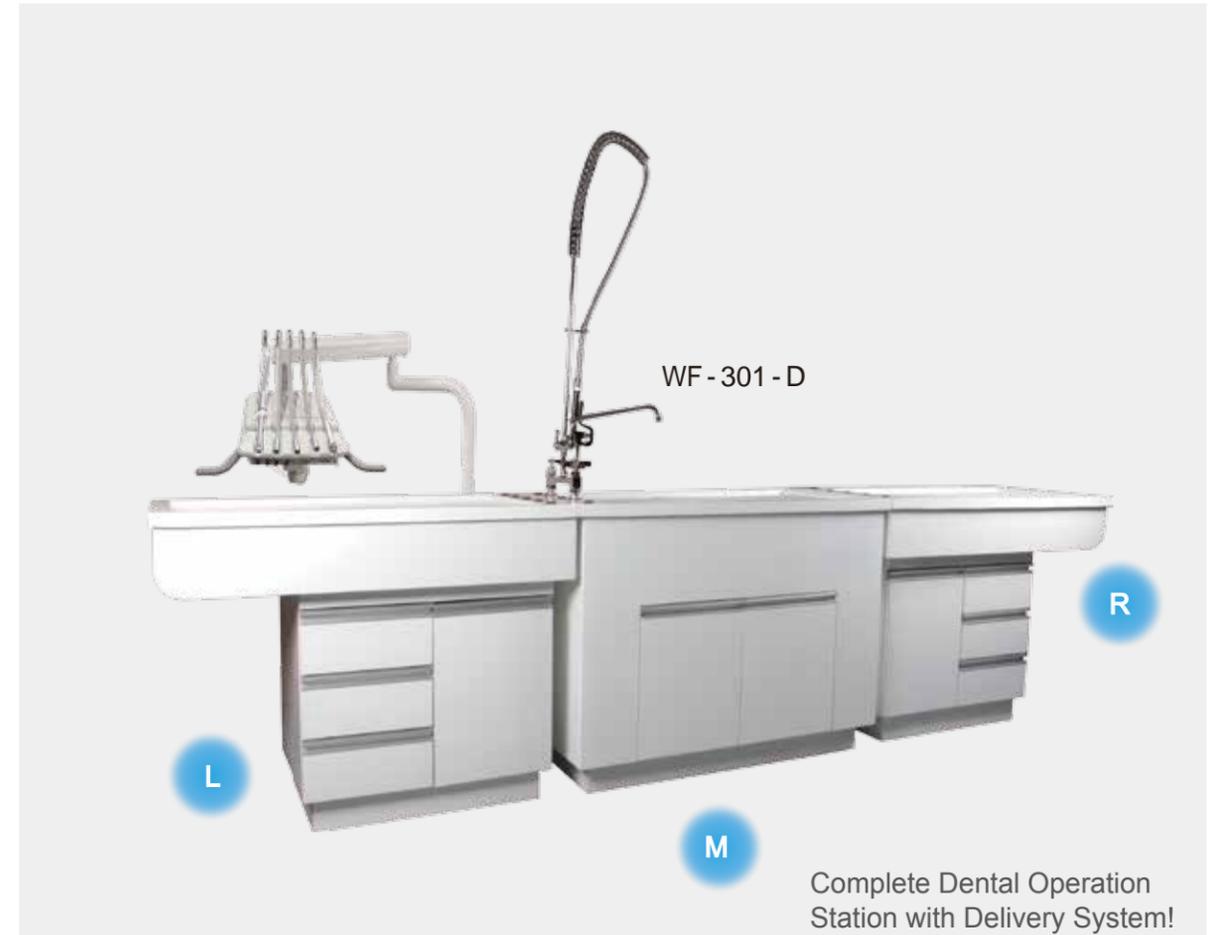

WF-301-D

L

M

R

Complete Dental Operation Station with Delivery System!

Preparation Table Dental Table with Stainless Steel Tub

FT-851-AS

The first choice in Preparation Tables.

This beautiful dental table; made from Artificial Stone Tub and top quality cabinets underneath will be a focal point in your clinic. Less clinical in look and feel than stainless steel, this table offers a softer esthetic. The attractive eggshell white tub and supporting cabinet provide a warm feel to the examination area. The sturdy plastic coated white grid blends perfectly with the white tub and cabinet underneath. Competitively priced, the AS table is a cost effective option when compared to the stainless steel tub version. PVC cabinet is available as an option.



1. Solid One-piece Tub
2. Plastic Dip-in Treated Grid
3. Pre-drilled Holes for Tap Installation
4. Plastic Waste Pipe
5. Smooth Sloped Drainage
6. Stainless Drain with Filter

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-851M-AS	72.50	133.00	120.5×60.5×91	47-7/16×23-13/16 ×35-53/64	125× 65×105H	0.85
FT-851R/L-AS	69.50	127.70	120.5×60.5×91	47-7/16×23-13/16 ×35-53/64	125× 65×105H	0.85



Complete Dental Operation Station with Delivery System!

Preparation Table Dental Table with Stainless Steel Tub | FT-851-SS

Fully constructed out of 304 stainless steel for durability.


We encourage you to offer us your special requirements.

We will customize the tables to fit your specifications.

Can be mounted to working island or wall. Left handed, right handed and middle versions available.

High-pressure tap and faucet available, please refer to our WF-301 series.

Could also be equipped with flat counter top, please refer to 850SS-R/M/L-SS.

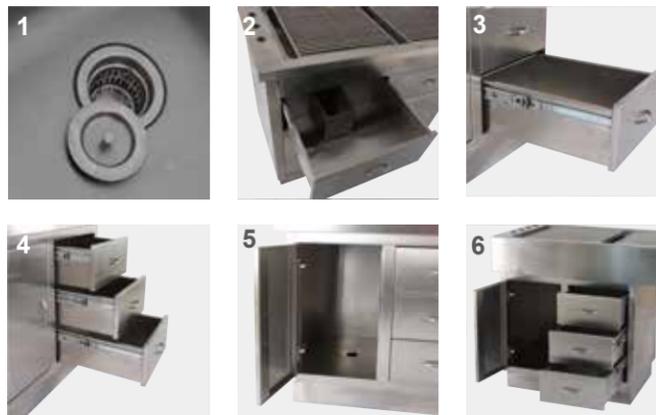


- 1. S/S hair strainer
- 2. Useful storage drawers
- 3. Storage cabinet
- 4. Storage cabinet
- 5. Slope on top tub for fluent draining
- 6. Slope on top tub for fluent draining


ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-851R/L	72.50	133.00	120×60×90	47-1/4×23-5/8 ×35-7/16	125× 65×104H	0.85
FT-851M	69.50	127.70	120×60×90	47-1/4×23-5/8 ×35-7/16	125× 65×104H	0.85

Size


L-48", W-24"

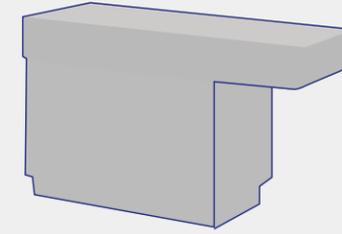


L-60", W-24"

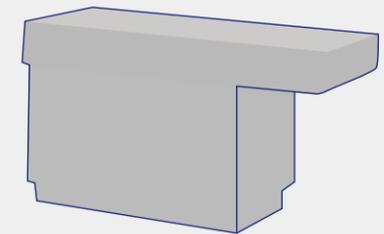

Knee Space

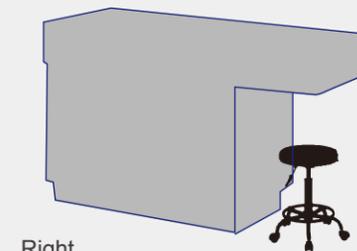

Right Front View

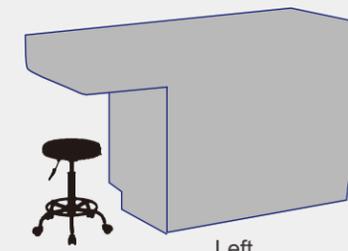

Left Front View

- Right Side
- Left Side

Knee Space determines which side of the table storage will be built.


Cabinet Configuration



- 1 Door + 3 Drawers
- 6 Drawers
- 2 Door

Not available for Tub Tables

Table Top Material

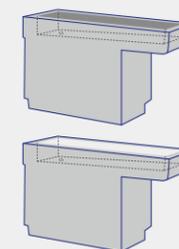

- Stainless Steel
- Artificial stone

Cabinetry Material



- Stainless Steel
- Painted Wooden Board



Consultant Room Island

Combined by columns and several preparation tables.

Versatile, can be configured with varied elements like prep tables, sinks, tubs, and fold-down tables.

The whole combinations, including length and width, can be modified to meet you specs. Electrical wires, plumbing and oxygen lines can all be surrounded if the columns are built from floor to ceiling.

All side panels are functional to accept accessories like anesthesia station, Otoscope/Ophthalmoscope sets, exam lights, glove boxes, tissue trays, etc.

The columns and tables can be tailor made to your requirements, just let our sales representatives know your needs; we can offer prompt professional design assistance.



Custom order welcome
Great flexibility





FT-828H

V-top Operation Table | FT-828/828H

A two panel table top system that adjusts from horizontal to near vertical to hold the pets in the desired position during surgery.

Overall size of the top: 20"W× 55"L.

Powerful and durable electric lifting actuator holds 120kg (240lbs) on the surgery top.

Two sliding, quick release easy controls are provided on each side of the rail to offer secure pet restraint.

10cm (4") stainless steel trough encourages fluid run off into a drain.

Table can be lifted from 28" to 47".

Worktop Swivel Angle: 0°- 45°.



Sliding drain trough



Operation surface angled



Angle adjuster mechanism



Side panel handle



FT-828



Buffered hinges and air spring added to cushion the dropping impact of the whole top and V-top side panels



Improved waste channel and drain hole design allow for fluidly drain flow. Close fitting translucent bucket improves waste collection.



Buffer hinge



New 3" casters installed for smoother movement and easier braking.



Heating system with accurate temperature monitoring and control (Only for FT-828H)

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-828	70	85	140×50×70-120	55-1/8×19-11/16×27-9/16×47-1/4	146×58×85H	0.73
FT-828H	80	95	140×50×70-120	55-1/8×19-11/16×27-9/16×47-1/4	146×58×85H	0.73



FT-886H

Universal Surgical Table | FT-886/886H

V-top Operation Table

Powered by sturdy and reliable Danish made lift column.
 Telescoping pedestal style, slim lifting column leaves extra leg space underneath, and provides room for other equipment that needs to be close to the operating area (i.e. X-Ray equipment).
 V-top table with buffer hinges and air spring secures the top where positioned.
 Heating system with accurate temperature monitoring and control.
 The IV pole mount is welded on.
 Improved waste channel and drain hole design promote drain flow. Translucent bucket provides visibility of waste collection.
 Water resistant foot pedals.



Sliding drain trough



Angle adjuster mechanism



Side panel handle



Operation surface angled



FT-886



Heating system with accurate temperature monitoring and control (Only for FT-828H)



Improved waste channel and drain hole design allow for fluidly drain flow. Close fitting translucent bucket improves waste collection.



Buffer hinge



New 3" casters installed for smoother movement and easier braking.



Buffered hinges and air spring added to cushion the dropping impact of the whole top and V-top side panels

ITEM	Net Weight	Gross Weight	Dimension (cm)		Volume (cbm)
	(kg)	(kg)	Stand	Package	
FT-886	70	85	140×50×61-101	146×58×76H	0.65
FT-886H	80	95	140×50×61-101	146×58×76H	0.65



FT-887H

New Universal V-top Surgery Table | FT-887

Optional heated top disperses heat evenly across the tabletop.
 Accurately controlled temperature is the best way to prevent hypothermia, without any risk of overheating.
 Transparent drainage bucket collects waste fluid, and offers direct vision of waste liquid volume.
 Reliable electric lift column (powered by Linak), provides long service life.
 Removable foot pedal set provides easy and comfortable table height adjustments.
 Tight and secured V-top and table tilting adjustments.
 Unique swivel tray design helps to keep your tools very organized and provides a handy space to collect tissues/organs.



Slim base leaves good legroom



Transparent waste bucket



Easy tilting adjustments



V-top panels setting



Reliable electric lift column



Swivel tools trays



Functional Cutting for laparoscopy operations



Separate temperature control for both panels



Triangle extension board for small animals



Clear waste fluid collection bucket

V-top Deluxe Operation Table | FT-888

More ergonomic working position, allows veterinarians to get even closer to patients and avoid knee or back pains/injuries caused by bad postures.

The special cutout contour provides even more comfortable position for surgeons, with matching inserts available to put it back to regular table top.

Telescoping pedestal column with ultra slim size leaves extra big legroom underneath.

Powerful loading capacity, hold up to 150kg patients.

Heating Tops with accurately control temperature, the best solution to avoid hypothermia on patients.

Omron Thermostats ensures ease of use and long termed reliability.

Electric tilting system makes the top suit Trendelenburg position perfectly.

Tapering extension tray allows easy access for even the very small patients.

Adjustable top panels, forming comfortable V shapes at 0 to 45 degrees, with buffering system for gently drop-off.

Easy access with both slide-out trough and exposable V-top panels for deeply cleaning purpose of whole fluid draining system.

Surrounding Rail System, with three sides of the table boarders fully covered by the neatly finished rails, good for both animal restraint cleats and accessory holders.

Universal accessory holder design, one easy step to attach all the functional accessories, including drape, surgical tray, IV post, infusion pump holder, monitor support post, etc.



Electric driving system for table top tilting



Foot pedal control for height adjustment



Multi-functional Operation Table | FT-826

There are two options for this wonderful table, FT-826 hydraulic, FT-826E electric.

This unique table is very versatile.

The three panel operation top makes it possible for you to enlarge the operation platform size from 19" wide 28 1/2" wide. This unique feature is not available elsewhere.

Two side panels can adjust from horizontal to near vertical to hold the patient in the desired position during surgery.

Overall size of the top: 19"W (to 28 1/2"W) × 58"L.

Worm bar driving system allows you to move and tilt the top, and adjust the angle at 0 - 45 degrees as needed.

Stainless steel operating table top with tilting function, convenient for various breeds and sizes of animals.

Superior quality table is rock-solid at any height, adjustable from 21.5" to 42".

Ergonomically designed, each side panel is 9", based upon average forearm length. Three waste liquid troughs under the top boards. Specially mounted fixture allows for easy restraint.



Sliding Side Panels And Removable Middle Panel Extend Operating Surface



Sliding Side Panels And Removable Middle Panel Extend Operating Surface



Swiveling Side Panels With Adjustable Angle From 0-45 Degrees



Smart Worm Bar Driving System Ensures Variable Angle for Work Face



Smart Worm Bar Driving System Ensures Variable Angle for Work Face



Smart Worm Bar Driving System Ensure Variable Angle for Work Face

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume
	(kg)	(kg)	W×D×H(cm)	Stand	W×D×H (inch)	Package (cbm)
FT-826	100	133	145×57×55-107	57-3/32×22-7/16 ×21-21/32×42-1/8	153×57×82H	0.68



Ophthalmology & Dental Table | FT-853

Dual functional table.

Fully constructed from 1.2mm thick 304 grade stainless steel; the table is very durable and reliable.

Strong base incorporates a high quality electric lifting system.

Smart design V-top side panel system, adjusts easily in 3 stages.

Tilting tabletop ensures a better overview of the patient. It also can utilize gravity to help reposition organs.

Rails run the full length of the table allowing cleats to be positioned anywhere along the way.

Waste liquid runs through a stainless steel trough. It slides out to be cleaned easily.

S/S diversion pipe for both soft pipe connection and waste bucket.

Optional end tray & grid for ophthalmology and dental surgery available.



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH(cm)	Stand	Package	
FT-853M	45	65	90×45×73-123	35-7/16×17-23/32×28-47/64×48-27/64	93×48×75H	0.33
FT-853L	47	67	130×55×73-123	51-3/16×21-21/32×28-47/64×48-27/64	133×58×75H	0.58
FT-853M-ECO	45	65	90×45×73-123	35-7/16×17-23/32×28-47/64×48-27/64	93×48×75H	0.33
FT-853L-ECO	47	67	130×55×73-123	51-3/16×21-21/32×28-47/64×48-27/64	133×58×75H	0.58



Electric Lift Operation Table | FT-854

Sturdy Z style stainless steel base ensures stability.
 Luxurious satin finish tabletop surface is easy on the eyes (even under bright surgery lights).
 One piece, reinforced tabletop, made from 1.2mm 304 stainless steel, has unique V-shape drain trough in the middle to ensure no-spill fluid collection.
 3 section removable panels with tiny holes allow for moisture to drain away from the operation site.
 Add four castors to convert this unit into a mobile operation table or emergency table.



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-854	54	74	120×60×52-101	47-1/4×23-5/8 ×20-15/32 ×39-49/64	125×65×66	0.54



Tilt Top Exam Table | FT-855

This is a versatile, utilitarian table that will serve you well and take a beating. Made from top grade 304 stainless steel.
 Adjustable top tilts to position a patient correctly.
 Full-length rail and easy use cleats make animal restraint simple.
 Water resistant foot pedals and S/S power box case provide safety during daily use.



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-855	62	82	120×60×69-171	47-1/4×23-5/8 ×27-11/64×67-21/64	125×65×85	0.69



Low-Low Operation Table | FT-862

Dropping to within 10" of the floor makes this operation table an exceptional tool. Unbeatable low height provides easy access.

One piece reinforced tabletop, made from 304 grade stainless steel, is designed for long term use.

V-shape drain trough in the middle provides no-spill fluid collection.

Scissor action is supported by noise free nylon bearing in specially milled channel.

Scissor action has nylon washers at every joint that prevents friction, resulting in the infrequent need for maintenance.

The top section is made up of three grid panels. They have tiny holes that allow for constant fluid runoff.

Add 4 casters and this table converts to a mobile station.



Low-Low Dental Table | FT-863

Heavy duty lifting mechanism and very durable linear actuator system, able to lift up to 250lbs.

Low-low profile design, drop down to 10" high at the bottom, making it relatively easy for larger and heavier patients to be handled.

Tapered top board shape, leaves more room for vets to be closer to patients.

Removable top grids, resting on lip inside tub for an easy flush work surface, no hooks hanging over to snag clothing.

Also could be used for other purposes, like transportation, exam or normal treatment.





Rotating Operation Table | FT-827

49" (L) x 24" (W) flattop.
 Table can be lifted from 55cm to 100cm (22"-40").
 Flat top revolves 360 degrees, and tilts to let the water drain into the hanging waste receptacle.
 Multi-Directional spring with a pull-up control panel secures the top in almost any position, from horizontal to near vertical.
 Ball bearing transfer mechanism makes the rotation smooth.
 Fully-automated operation improves the work effort and manages time more efficiently.
 Powerful and durable electric lifting actuator holds up to 120kg (240lbs) on the surgery top.

Dual sliding, quick release control attachments mounted on each side of the rail provide secure pet restraint.
 Slightly raised center encourages fluid run-off into a perimeter trough (to a drain outlet).



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH (cm)	Stand WxDxH (inch)	Package	
FT-827	60	63	125x60x55-100	49-7/32x23-5/8x21-21/32x39-3/8	130x66x71H	0.62



FT-825E

Hydraulic or Electrical Stainless Steel Operation Table | FT-825 Series

Classic lifting mechanism:
 FT-825 hydraulic.
 FT-825E electric.
 FT-824 is medium size at 90cmx60cm (36"x24").
 FT-825 and FT-825E are large size at 120cmx60cm (48"x24").
 Sturdy tables with stainless steel surface and reinforced construction that discourages wobbling and tipping.
 Hydraulic foot pedals control the height from 51.5 - 101cm (20-1/4" - 39-3/4").
 Steel basin at one side of the surface.



Table top grids are optional parts

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH (cm)	Stand WxDxH (inch)	Package	
FT-825	47	67	120x60x51.5-101	47-1/4x23-5/8 x 20-9/32 x 39-49/64	124x64x71H	0.55
FT-825E	46	66	120x60x65.5-106	47-1/4x23-5/8 x 25-25/32 x 41-47/64	124x64x81H	0.65



FT-871E

Multi-functional Economic Electric Operation Table | FT-871/871E

This table provides valuable service at cost effective prices. Perfect for quick procedures, exams, changing dressing, suturing, dental exams, etc. Optional end tray is available.



Optional tool tray for your choice.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-871E	35	38	120×60×52-100	47-7/16×23-13/16 ×20-15/32 ×39-3/8	124×64×52H	0.41



FT-872E

Multi-functional Lifting Operation Table | FT-872/872E

Two options: FT-872 hydraulic and FT-872E electric. Very versatile tables; cost effective tables. Powder coated or stainless, hydraulic or electric. The choice is yours. Economic, sturdy, well built. Optional end tray is available.



Optional tool tray for your choice.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-872	40	43	120×60×51-101	47-7/16×23-13/16 ×20-5/64 ×39-49/64	125×65×52H	0.42
FT-872E	40	43	120×60×65-106	47-7/16×23-13/16 ×25-19/32×41-47/64	125×65×65H	0.53



Transportation Table

| FT-835

A unique one-act-lift system with automatic locking latch, allowing quick release operation too.

Light weighted, compact in transit, good for ambulance and accident rescue.

Stainless steel exam top is 48" by 24", with removable stretcher.



1. Removable stretcher on top
2. Super easy one act lift operation
3. Self-adapting lock latch
4. Medical class casters



Wall Mount Exam Table

| FT-843

Maximize valuable space with this wonderful fold-up wall mount table.

The best solution to balance the working capacity between peak and normal seasons.

Radius corners eliminate the risk of user injuries.

Big loading capacity, up to 200lbs.

Fluid retaining border design, no worry about fluid overflows.

Easy installation and assembly works required.

Table top measures 36" by 24".





Stainless Steel Bi-level Vet Tub | TT-201

Manufactured with seamless welded construction from SUS304 stainless steel. This vet tub is rust resistant. Four sturdy stainless steel legs support this versatile table. The feet adjust to provide level footing and a secure platform on which to work. The shallow portion of the tub is 4.5" deep and drops to a 21.5" depth. The inner dimensions measure 20"W x 57"L; outer dimensions are 24"W x 61"L. The height of the tub top is 39.5", and the bottom of the tub sitting 12" off the floor. The drain opening is designed for a 3.5" diameter standard drain. All holes for the faucet and drains are pre-drilled. Stainless steel flat surface grill included.



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH(cm)	Stand WxDxH (inch)	Package	
TT-201	60	81	155x60x90	61-1/32x23-5/8x35-7/16	164x65x101H	1.08



Tub Table | TT-202

Cost effective and versatile, designed for multiple uses; tub table/ dental table and preparation table. It's 100% stainless steel made and suitable for any area. It easily tolerates high humidity and wet conditions. Flat surface stainless steel grills included. These are sturdy, secure and supportive grills. 4" deep on one end and 6" on the other, 2" incline drop speeds up waste liquid flow. This table is flat packed and shipped with the legs and brackets disassembled. This makes it very cost effective. Assembly is basic, quick and easy. Four fixture holes with universal diameter are pre-drilled, making installation of water faucets a simple task.



ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH(cm)	Stand WxDxH (inch)	Package	
TT-202	30	33	140x65x90	55-1/8x25-19/32x27-9/16	144x69x34H	0.34



Universal Surgical Table

FT-881

Tabletop is 120cmL x 60cmW (47 1/4"L x 23 1/2"W); very practical size.
 Extended adjustable height range from 530mm to 930mm (21" – 36 1/2")
 Telescoping pedestal style, slim lifting column provides extra leg space underneath, and provides room for other equipment that need to be close to the operating area (i.e. X-Ray equipment).
 Wide stance stable base and water resistant proof foot pedals are key features to this marvelous table.
 Powered by Linak brand electric lifting column, offering world class quality and durability.
 Equipped IV pole holder and pet restraint hooks.
 Add 4" casters and easily convert this table to a mobile surgical unit.



Electric Lift Examination Table

FT-861 Series

SUS304 stainless steel table top and electric lifting base provide a great table you'll rely on for years.
 Lower down to 31cm (12 1/4"), enabling older, infirm and larger dogs to easily walk on.
 Drainage hole at corner makes waste liquid flow away easily.
 Incredibly stable X shape lifting base, discourages wobbling at any height.
 Equipped with a height adjustable IV pole.
 Brand new model FT-861S: weighing functional for purchases (optional).

FT-861S



Universal Operation Table

FT-882

Neatly welded and nice polished table top, constructed from SUS304 stainless steel.
 Daily examination and simple operation could be executed on this table perfectly.



Electric Lift Dental Table

FT-852

SUS304 stainless steel table top and electric lifting base will provide you with great service for years to come. Pet friendly, comfortable flat surface grid on top. Sloped tub to allow waste fluid to drain easily.
 Ultra-stable X shape lifting base, discourages wobbling at any height. A very useful table.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH (cm)	Stand WxDxH (inch)	Package	
FT-852	54	87	140x65x41-110	55-118x25-19x32 x16-9x64 x43-5/16	145x70x53H	0.53
FT-861	48	80	130x65x31-98	51-3/16x25-19x32 x12-13x64 x38-37/64	135x70x44H	0.41



Exam Table

FT-870L

Sturdy, versatile, all purpose fixed stand table.
 Stainless steel table top, aluminum alloy legs.
 Diagonal embossed tabletop surface aids in dispersing liquid waste efficiently.
 Terrific value. Flat packed. Some assembly required.



Small/Medium/Large Exam Table

FT-823S/M/L

Flat packed for economical shipping, assembly required. Tabletop constructed from 1.0mm quality SUS304 stainless steel, supported by sturdy S/S frame.
 The frame has anchor hooks welded underneath, providing options for securing and restraining pets.
 Drainage hole at one side of surface to collect waste residue into a stainless steel basin.
 Constructed of quality 304 stainless steel board; sturdy, easy to clean up.
 Includes 170cm tall IV pole and stainless steel waste bucket.



Economical Examination Table

FT-873L

Sturdy pedestal and wide stance base make the most of your space and provide plenty of operating room.
 Stainless steel table top with diagonal embossing structure, disperses liquid waste efficiently.
 Economically priced and provides high-end performance. Flat packed for ease of transportation.
 Some assembly required.



Economical Examination Table

FT-874X

Stainless steel top with wooden support. Industrial strength powder coated legs/frame. Pre-drilled for IV attachment; includes folding IV pole.
 Shipped flat saves transportation cost.
 Some assembly required.

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH (cm)	Stand WxDxH (inch)	Package	
FT-870L	12	14	120×60×82	47-1/4×23-5/8 ×32-9/32	125×65×12H	0.09
FT-873L	20	23	120×60×90	47-1/4×23-5/8 ×35-7/16	125× 65×11H	0.09
					12×12×80H	0.01

ITEM	Net Weight (kg)	Gross Weight (kg)	Dimension (cm)			Volume (cbm)
			WxDxH (cm)	Stand WxDxH (inch)	Package	
FT-874X	26	28.5	120×60×93	55-1/8×23-5/8 ×36-39/64	122×64×20H	0.16
FT-823S	10	14.5	75×50×85	29-17/32 ×19-11/16 ×33-15/32	80×55×18.5H	0.08
FT-823M	12	17.0	90×55×80	35-7/16×21-21/32 ×31-1/2	95×60×18.5H	0.10
FT-823L	14	19.5	120×60×80	47-1/4×23-5/8 ×31-1/2	125×65×18.5H	0.15



Deluxe Mobile Tool Cart w/ Cabinet

KB-521

Stainless steel construction, practical and functional design. Two big drawers and 9 stainless steel cup holders for tools storage. Smooth, easy to clean surface. Stainless steel tray for tools storage; waste bin hook on the side of cart. 4" high quality casters with brake, easy to control.



Mobile Tool Cart

KB-524

Stainless Steel construction. Stable, durable and scratch resistant. Smooth, easy to clean. Waste bin hook on one side, stainless steel tool tray on the other side. 4 high-quality casters, 2 of which are with brakes.



Mobile Tool Cart

KB-520


Quality 304 stainless steel. Smooth, flat, easy to clean surface. Two shelves, strong enough to hold all your operating tools.



Mayo Table

KB-531

Compact yet large enough for most surgical procedures. Stainless steel tray measures 17.5"×13.5"×1.5"d (44cm×34cm×4cm). Easy-slide extension posts allow height adjustments from 31.5" (80cm) to 44" (112cm).

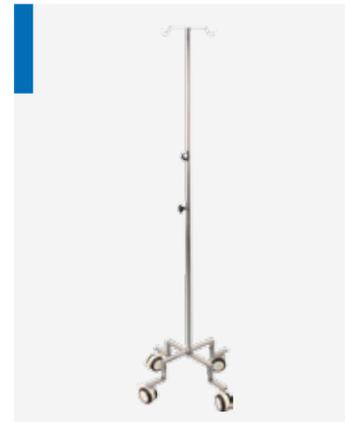

IV POLE

IP-101

1.3m (4 ¼") to 2.2m (7 ¼") adjustable height. Stand size: 120×53×12cm. Package size: 136×54×12cm.



Heavy Duty Infusion Stand

IP-102

Strong and versatile, enabling easy infusions and drips wherever you need it. Die cast heavy duty base lowers the center of gravity and provides great stability. Two hooks on top standard, 4 hooks top version available too. 5 durable ball bearing casters at the bottom for smooth movement. Height adjusts from 128 to 194cm.



WF-301-D


WF-301-W

WF-302-D

WF-302-W

Professional Multi-functional Faucet | WF-301/302

For cleaning wounds and irrigating during dental procedures. Swiveling faucet pipe and cold/hot water control. Made from cooper, with polished chrome plating. Pressure water spray and stainless steel hose. Package size: 65×35×45cm. Gross weight: 6.5kg. A lot more options available, please check our website to know the latest items. D means Deck Mounted and W for Wall Mounted.

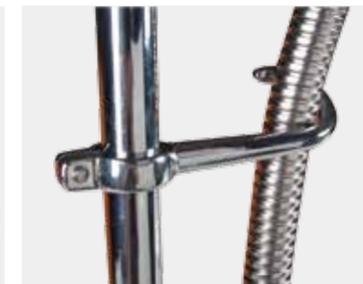



Scale Table | FT-822M

This extremely accurate Scale Table is ideal for your clinic

Platform size 90cm (35 1/2") x60cm (23 1/2"); weight range from 20g (3/4 oz) to 100kg (220lbs).

Pinpoint accuracy is our foremost concern when weighing pets. Treatment and proper administration of medicines are determined by weight.

SUS304 stainless steel top that is fully sealed (water resistant) provides durability and accuracy.

Pre-drilled drainage hole on platform. Transparent skid-proof mat provides a non-slip surface.



- 1. Quick Fitting Plug
- 2. Control Panel
- 3. Hexagon Bolts Installation
- 4. Drainage Hole
- 5. Flat Packing
- 6. Adjustable Leveler

ITEM	Net Weight	Gross Weight	Dimension (cm)			Volume (cbm)
	(kg)	(kg)	W×D×H(cm)	Stand W×D×H (inch)	Package	
FT-822M	48	67	90×60×87	35-7/16×23-5/8 ×34-1/4	92× 65×39H	0.25



Walk-on Scale
WOS-4824

Large size at 120cm by 60cm (48"×24"). 4" thick; easy walk on access for large animal. Capable of handling large patients (up to 150kg). Sturdy scale top constructed from premium 304 stainless steel.



Walk on Scale
WOS-3624

Medium Size: 90-60cm 4" thick; easy walk on access for large animal. Capable of handling large patients (up to 150kg). Sturdy scale top constructed from premium 304 stainless steel.



Walk-on Scale
WOS-FL-2010

Weighing capacity up to 20kg. Powered by 4 standard AA batteries. LCD display, compact, portable. Product size: 60x34x9cm.



Walk-on Scale
WOS-FL-4010

Weighing capacity up to 40kg. Powered by 2 AA batteries. Thin top tray. Size: 60.5×32×7.5mm



Medical Pendant

This medical pendant system is designed to facilitate movement, support loads and supply medical equipment with electrical, gas, or data.

Varied medical gas outlets, including oxygen, vacuum, air and NO2.


Best positioning of medical equipment in the area, also makes the sterilization works much easier.

Finished with antibacterial powder coating.

A super easy access to high-tech medical devices put all the precious devices in well organizing.

Safe to use and feature a long service life with minimum maintenance demands.



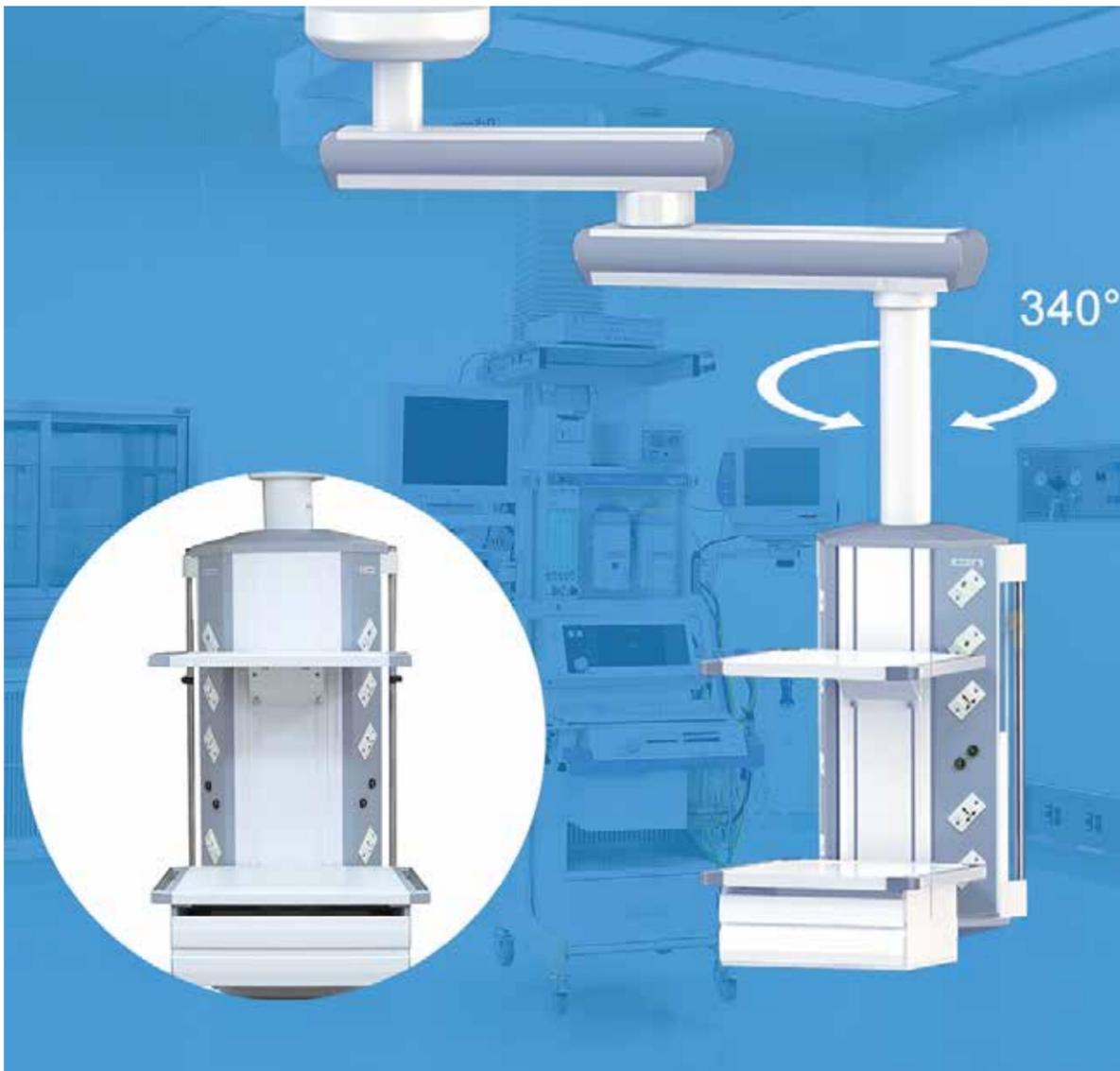

INTERFACE

Different types of interfaces for customers to choose from Telephone / Internet / Video Communication Interface

ELECTRICAL SOCKETS

Different types of adaptors for customers to choose from

French	Chinese	American	Universal
German	British	Australian	Indian

