




Our Solutions


With superior image quality and intuitive imaging software, the ExamVue Duo digital imaging system delivers fast image acquisition, easy workflow management, and lower dose for your patients.

The ExamVue Duo system comes equipped with a strong suite of chiropractic tools, and is available in any configuration for your practice needs.

JPI ExamVue Duo Package

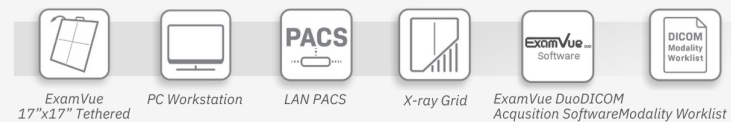

Optional Upgrade

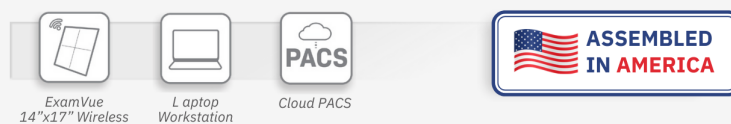

5-5-5 Warranty Program

JPI is not your average imaging solutions provider. Reach out to us! Our Account Managers are happy to discuss the imaging solutions and coverage options that are best for you.

- X **5 years DROP COVERAGE PROGRAM**
- X **5 years REMOTE SOFTWARE SUPPORT**
- X **5 years WARRANTY**




Contact Us Today

Sales: (787) 918-1080 | **Service:** (787) 631-7157



www.wmedsolutions.com

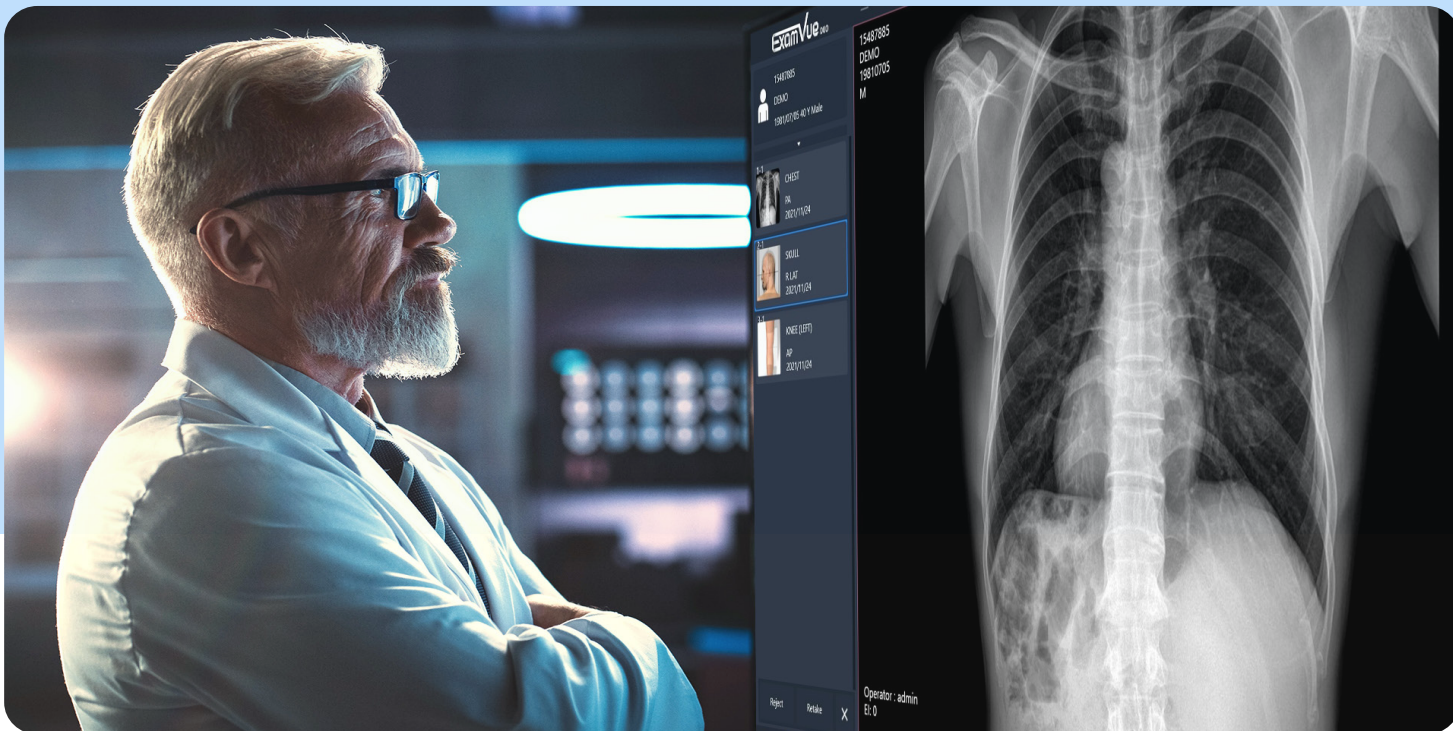

SEE MORE. DIAGNOSE WITH CONFIDENCE.


Less than 100 μm resolution.

Designed for the discerning hospital market, the **ExamVue Apex** flat panel detectors offer unparalleled performance and versatility. Compatible with all x-ray generators, these detectors seamlessly integrate into your existing infrastructure, delivering exceptional image quality. Leveraging Cesium direct-deposit technology and high resolution, the ExamVue Apex ensures exceptional image clarity for all patient sizes and eliminates the need for generator wiring.

The ExamVue Apex's **glassless design** ensures a lightweight and streamlined workflow. Equipped with built-in **drop sensors**, the device operates independently on a dedicated battery with a six-month charge, eliminating the need for constant panel power. These sensors actively monitor for sudden movements, stops, or impacts, recording any incidents for review.


LIGHTWEIGHT. BIG IMPACT.

17"x17" Wireless

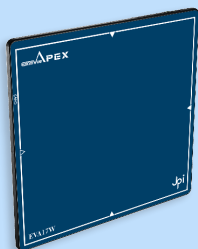


6.6 lbs with battery

14"x17" Wireless

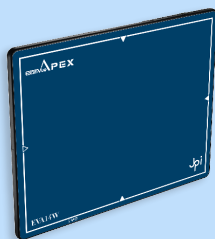

5.7 lbs with battery

10"x12" Wireless

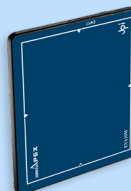


3.5 lbs with battery



ExamVue DUO

Advanced Tool-set

Center Point
Orthogonal Line
Torg's Ratio
Pelvic Obliquil
Baseline Measurement
Jackson Angle
Ishihara Index
George's Line
Cervical Curve
Lumbar Curve
Ilium Analysis
Spine Analysis
Spine Label

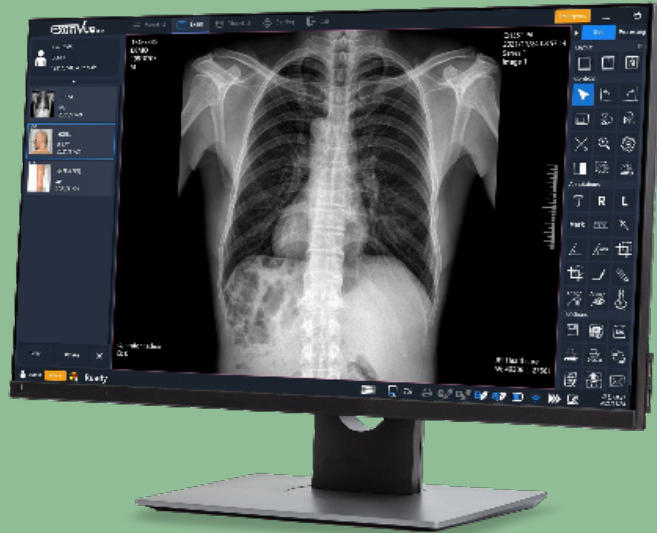

User-friendly
Design

Simple & Fast
Workflow

Easy To Use


Software At A Glance



X-ray Acquisition



Images Preview

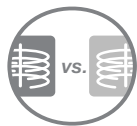

Compare Side-by-side



Stitching Software



Grid Line Remover



Patient Registration



Emergency Mode



Send DICOM to PACS & Print

Additional Features

- ▶ **Seamless Exam Acquisition:** Perform complete exams from start to finish without computer interruptions.
- ▶ **Automated Image Processing:** Enjoy automatic processing for all views with manual fine-tuning options.
- ▶ **Real-time Status Updates:** Stay informed about panel connection, PACS, MWL, and hard drive space.
- ▶ **Integrated Workflow:** Benefit from MWL and Procedure Code Matching for seamless integration with your RIS/HIS or our standalone MWL generator.
- ▶ **Versatile File Formats:** Export and import DICOM and JPG files for easy sharing and archiving.
- ▶ **Comprehensive Support:** Access integrated technique charts, positioning guides, and exposure index for all programmed views.

Optimized Workflow

